


Welcome!

VITLAB has over 100 years of tradition. The company VITRI GmbH & Co. KG was established in 1908 in Mühltal, and the laboratory division was spun off in 1989 as VITLAB. Today, VITLAB is one of the leading manufacturers of liquid handling instruments and performance plastic laboratory products for industrial and scientific applications. We develop and manufacture these laboratory products at our own production facility.

Our extensive range of products provides optimal support in your laboratory work in a wide variety of application ranges. Regardless of whether your work involves volume measurement, sampling or storage: VITLAB products will facilitate it while continuing to ensure that you achieve perfect results.

We hope that this new catalogue will serve as a valuable resource to aid in your lab work. We would be glad to respond to your enquiries, and look forward to receiving suggestions and ideas from our users.

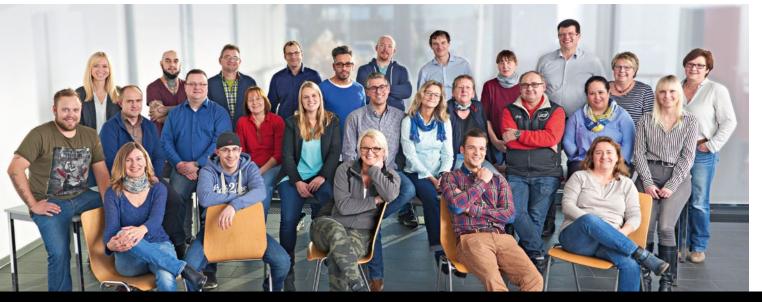



This catalogue describes our products and provides all essential information. You will find detailed data regarding the various plastics clearly presented in the chapter entitled "General and Technical Information". To simplify your search, our product range has been categorised into the following areas of application: dispensing, pipetting, titration, volume measurement, measuring and transferring, sample preparation, saving and storing, lab assistants.

Under "Volume Measurement", for example, you will find a wide variety of classical volumetric instruments such as volumetric flasks, measuring cylinders, and associated accessories.

In addition to the range of products that appears in our catalogue, we also produce plastic products according to customer specification. For example, bottles and beakers required for special tasks can be commissioned with a specified geometry and thickness. Optionally, the products can be printed with an individualised scale or with customised labelling. Labware for promotional use can be designed and printed with a company name and logo. We can even accommodate requests for individualised packaging formats, materials, and designs. Further information is available in the VITLAB® Promotional chapter.

Many possibilities can be realised – don't hesitate to ask what we can do for you!



VITLAB Your reliable partner


Certified quality

Independent inspections and routine internal audits guarantee the effectiveness of VITLAB's quality management system throughout the entire company, from development to shipment. As a result, the phrase 'Made by VITLAB' has become synonymous with quality.

Over 98% of our product line is made in Germany. Supplemental procedures such as tempering and volume testing are conducted in our own facilities, which guarantees the highest possible product quality

and measurement accuracy. Our continuous improvement paradigm supports our goal of 0% failure.

The VITLAB Quality Management System has been continuously certified since January 1994, according to DIN EN ISO 9001. Active stewardship of the environment is an equally strong pillar of our business philosophy. VITLAB has been certified according to DIN EN ISO 14001 since May, 1999.

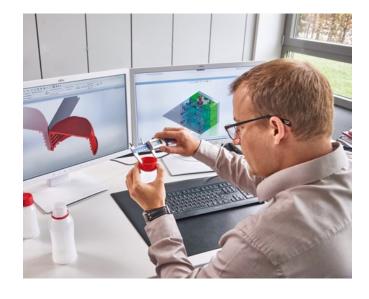



Prompt deliveries Competent customer service

The highly efficient logistics in the Großostheim production facility facilitate the shortest possible delivery times for all products listed in the catalogue. For the standard articles, we strive for an availability of over 94%.

Due to its intensive partnerships with distributors in over 70 countries, VITLAB can offer sound on-site advice, individual support, and quick answers to your questions. Our qualified product training sessions provide comprehensive technical and applicationoriented information on using our products. Should problems arise, our expert repair service keeps downtime as short as possible.

VITLAB products can be ordered from specialist dealers worldwide. Our authorised sales partners can be found on the internet at:

www.vitlab.com

Or contact us directly.





For your information



Your contact Customer Service

Our Customer Service staff is at your service to provide you with competent advice and answers to all your queries and questions concerning offers, orders and deliveries. Our Product Management and Sales Team are at your disposal – also "on site" – with any technical information or assistance that you might require for your application.

VITLAB GmbH, Customer Service

tel: +49 6026 9 77 99-0

fax: +49 6026 9 77 99-30

e-mail: info@vitlab.com web: www.vitlab.com catalogue numbers and designs may change after this catalogue is published. The illustrations used are for representation only, and the details may vary from the description. All measurements, with the exception of exact tolerances, should be understood as approximate values. Please keep in mind that the actual testing and measuring results can be influenced by a variety of factors that are beyond our control. Therefore, you should carefully check the transferability of the data applying it to a particular application.

Please understand that technical specifications,

The packaging units (PU) correspond to the minimum order quantities. All up-to-date information is also available on the internet at www.vitlab.com.

If you need additional details, please call us.

California Residents: For more information concerning California Proposition 65, please refer to www.vitlab.com/calprop65.

VITLAB®, maneus®, pipeo®, VITsafe™, VITgrip™ are brands of VITLAB GmbH.

Table of contents

Dosing Bottle-top dispensers and accessories	Starting on page 9
Pipetting Micropipettes, tips and accessories	Starting on page 21
Titration Bottle-top burettes, burettes and accessories	Starting on page 29
Volume measurement Volumetric flasks, graduated cylinders, pipettes, pipette helpers and accessories	Starting on page 35
Measuring and transferring Graduated beakers, scoops, spatula, forceps, funnels and accessories	Starting on page 53
Sample preparation Wash bottles, dropping bottles, dropping pipettes, spray bottles, Griffin beakers, watch glasses, Erlenmeyer flasks, magnetic stirring bars, mortars, pestles, equipment for material separation, desiccators, sample vials, staining cuvettes and accessories	Starting on page 63
Saving and storing Lab bottles, reagent bottles, storing bottles and transport containers	Starting on page 99
Lab assistants Connectors, valves, trays, drawer organisers and drying racks	Starting on page 113
VITLAB® Promotional	Starting on page 117
General and technical information	Starting on page 121

Clear product statements

Our aim is to provide you with a clear and comprehensive presentation of all relevant product information. For quick reference, we use the following symbols:



DIN ISO-compliant Class A volumetric instruments



DE-M marking for conformity certified products according to the German Measurement and Calibration Regulation



Food-safe products according to EU Directive No. 10/2011



Products with high protection for light-sensitive substances



Products individually packaged in PE bags, labelled with the article number, description and EAN code



Products that can be autoclaved at 121 °C (2 bar) according to DIN EN 285. Note restrictions!



CE mark according to EU Guideline 2004/108/EC, 93/68/EEC; 73/23/EEC, and 93/68/EEC



CE-IVD mark according to EU Guideline 98/79/EC

Perfection in Liquid Handling HIGHEST RELIABILITY IN DISPENSING



VITLAB® Dispenser line: genius², simplex², and TA²

VITLAB® genius² and simplex² bottle-top dispensers are a family of instruments with proven precision that offer many advantages in routine liquid-handling operations. VITLAB® genius² and simplex² instruments can be used for practically any task and are suitable for organic and inorganic solutions, while VITLAB® TA² dispensers have been specially developed for use in trace analysis and with highly concentrated media. As they are produced from materials with extremely high chemical resistance (e.g. PTFE, PFA, FEP, borosilicate glass and platinum-iridium), VITLAB® bottle-top dispensers are very robust and reliable and resistant against most acids, bases and organic solvents.



	VITLAB® genius²/simplex²/simplex ² _{fix}	VITLAB® TA²
Applications	Salt solutions, acids, bases, and many organic solvents	Specially for use in trace analysis for dispensing high-purity and highly concentrated acids and alkalis, as well as hydrogen peroxide, bromine and HF
Components in contact with media	Borosilicate glass, Al ₂ O ₃ -ceramic, FEP, ETFE, PFA, PTFE, platinum-iridium, PP (screw cap)	Various fluoroplastics (e.g., ETFE, FEP, PFA, PCTFE, PTFE), Al ₂ O ₃ -sapphire, platinum-iridium or tantalum (depending on the model)
Operating limits	Temperature: +15 °C to +40 °C Steam pressure: max. 600 mbar Viscosity: max. 500 mm²/s Density: max. 2.2 g/cm³	Temperature: +15 °C to +40 °C Steam pressure: max. 600 mbar Viscosity: max. 500 mm²/s Density: max. 3.8 g/cm³

^{*} Dynamic viscosity [mPas] = kinematic viscosity [mm²/s] x density [g/cm³]

General guide for dispenser selection (for the classification of dispenser media, see page 11).

Salt solutions	Acids and bases	Solvents	High-purity and highly concentrated acids and bases	Hydrofluoric acid (HF), bromine, hydrogen peroxide
VITLAB [®] ger	ius²/simplex²	VITLAB® genius²/simplex²		
			VITLA	B® TA²

Dosing

Recommended usage ranges for VITLAB® genius² / simplex² / simplex²_{fix}:

	Medium		Medium		Medium
0	Acetaldehyde	1	Copper sulphate	0	Methyl butyl ether
	Acetic acid, ≤ 96%		Cresol		Methyl ethyl ketone
0	Acetone	0	Cumene (Isopropylbenzene)		Methyl formate
0	Acetonitrile		Cyclohexanone		Methyl propyl ketone
0	Acetylacetone	0	Decane		Mineral oil (Motor oil)
0	Acrylic acid	0	1-Decanol	0	Monochloroacetic acid, ≤ 50%
0	Acrylonitrile	0	Diethylene glycol		Nitric acid, ≤ 60%*/**
0	Adipic acid	0	Dibenzyl ether		Nitrobenzene
0	Allyl alcohol	0	Dichlorobenzene	0	Octane
1	Aluminium chloride	0	Dichloroethane	0	Oleic acid
0	Amino acid	0	Dichloromethane	0	Oxalic acid
1	Ammonia solution, ≤ 20%	0	Diethanolamine	1	Perchloric acid
-1	Ammonium chloride	0	Diethyl ether	0	Petroleum
1	Ammonium fluoride	0	Diethylamine	0	Phenol
-1	Ammonium hydroxide, ≤ 20%	0	1,2 Diethylbenzene	0	Phenylethanol
1	Ammonium sulphate	0	Dimethyl sulphoxide (DMSO)	0	Phenylhydrazine
0	Amyl acetate	0	Dimethylaniline	1	Phosphoric acid, ≤ 85%
0	Amyl alcohol (Pentanol)	0	Dimethylformamide (DMF)	-1	Phosphoric acid, 85% + sulphuric acid, 98%,1:1
0	Amyl chloride (Chloropentane)	0	1,4 Dioxane	0	Piperidine
0	Aniline	0	Diphenyl ether	-1	Potassium chloride
-1	Barium chloride	0	Ethanol	-1	Potassium dichromate
0	Benzaldehyde	0	Ethanolamine	-1	Potassium hydroxide
0	Benzene	0	Ethyl acetate	1	Potassium permanganate
0	Benzoyl chloride	0	Formaldehyde, ≤ 40%		Propanol
0	Benzyl alcohol	0	Formamide	0	Propionic acid
0	Benzyl chloride	0	Formic acid, ≤ 100%	0	Propylene glycol (Propanediol)
0	Benzylamine	0	Gasoline	0	Propylene oxide
-1	Boric acid, ≤ 10%	0	Glacial acetic acid	0	Pyridine
0	Bromobenzene	0	Glycerine	0	Pyruvic acid
0	Bromonaphthalene	0	Glycol (Ethylene glycol)	0	Salicylaldehyde
0	Butanediol	0	Glycolic acid, ≤ 50%	0	Salicylic acid
0	1-Butanol	0	Heating oil (Diesel oil)	0	Silver acetate
0	n-Butyl acetate	0	Hexane	1	Silver nitrate
0	Butyl methyl ether	0	Hexanoic acid	0	Sodium acetate
0	Butylamine	0	Hexanol	1	Sodium chloride
0	Butyric acid	-1	Hydrochloric acid, ≤ 37%**	1	Sodium dichromate
-1	Calcium carbonate	-	Hydroiodic acid, ≤ 57%**	-	Sodium fluoride
-1	Calcium chloride	1	lodine / potassium iodide solution	-	Sodium hydroxide, ≤ 30%
-1	Calcium hydroxide	0	Isoamyl alcohol	-	Sodium hypochlorite
1	Calcium hypochlorite	0	Isobutanol	1	Sulphuric acid, ≤ 98%
0	Chloroacetaldehyde, ≤ 45%	0	Isopropanol (2-propanol)	0	Tartaric acid
0	Chloroacetic acid	0	Isopropyl ether	0	Tetramethylammonium hydroxide
0	Chloroacetone		Lactic acid	0	Toluene
0	Chlorobenzene	-	Magnesium chloride	0	Turpentine
0	Chlorobutane		Mercury chloride	0	Urea
0	Chloronaphthalene	0	Methanol	0	Xylene
1	Chromic acid, ≤ 50%	0	Methoxybenzene	-1	Zinc chloride, ≤ 10%
1	Chromic-sulphuric acid	0	Methyl benzoate	1	Zinc sulphate, ≤ 10%

The above data have been carefully checked and reflect the current state of knowledge. Always follow the instructions for use that accompany the instrument as well as the reagent manufacturer's instruction manual. In addition to the chemicals listed above, solutions of a wide variety of organic or inorganic salts (e.g., biological buffers), biological detergents, and cell culture media can be dispensed. Should you require information on chemicals not listed, please do not hesitate to contact us. Last updated: 09/17.

^{**} Use drying tube

1	Inorganic media
0	Organic media

^{*} Use ETFE/PTFE bottle adapter

VITLAB® genius² / simplex² / simplex_{fix}

VITLAB® dispensers are practically universally applicable and can be used with many **organic and inorganic solutions**. The materials that come into contact with media (borosilicate glass, Al_2O_3 ceramics, FEP, ETFE, PFA, PTFE, platinum iridium and PP) are resistant to most acids, solvents and bases.

The devices are equipped with a positive displacement piston and a fluoroplastic (PFA) sealing lip on the cylinder wall. This latter acts like a windscreen wiper to **prevent crystal build-up on the cylinder wall** from readily crystallisable media. The glass cylinder is also coated with a plastic material that reduces the risk of splashes should breakage occur. The telescopic filling tube can be adjusted smoothly to different bottle heights.

The practical screwing mechanism and inner toothed bar enable a fast and precise volume adjustment (simplex² and genius²). The simple-to-use calibration function helps meet all the requirements for test equipment monitoring with minimal downtime. Reagent loss while bleeding is avoided with the innovative recirculation valve (only genius²). The screwable discharge valve is equipped with an additional safety bulb and closes the dosing channel, if dispensing tube is not installed, so that no medium can escape.

VITLAB® genius², simplex² and simplex $_{fix}^2$ are completely autoclavable at 121 °C (2 bar) according to DIN EN 285.

Also available with DAkkS calibration certificate or individual certificate (at additional cost).



VITLAB® genius²



Bottle-top dispenser with variable volume and recirculation system.

Included in delivery: VITLAB® genius² (GL 45 thread), 3 or 5 thread adapters* made of PP, telescopic filling tube, recirculation tube, mounting tool, quality certificate and operating manual.

Volume ml	Graduation ml	A** ≤ ± %	A** ≤ ± µl	CV** ≤ %	CV** ≤ μI	PU	Cat. No.
0.2 - 2.0	0.05	0.5	10	0.1	2	1	1625503
0.5 - 5.0	0.10	0.5	25	0.1	5	1	1625504
1.0 - 10.0	0.20	0.5	50	0.1	10	1	1625505
2.5 - 25.0	0.50	0.5	125	0.1	25	1	1625506
5.0 - 50.0	1.00	0.5	250	0.1	50	1	1625507
10.0 - 100	.0 1.00	0.5	500	0.1	100	1	1625508

^{*} Nominal volume 2 - 10 ml: with adapters GL 25, GL 28, GL 32, GL 38, S 40 and telescopic intake tube (length 125 - 240 mm). Nominal vollume 25 - 100 ml: with adapters GL 32, GL 38, S 40 and telescopic intake tube (length 170 - 330 mm).

^{**} Error tolerance conforming to DIN EN ISO 8655-5, related to the nominal (maximum) volume marked on the device where the device, environment and distilled H₂O are at the same temperature (20 °C). Checks are done in accordance with DIN EN ISO 8655-6 with the device filled to capacity and with uniform and jolt-free dosing.


VITLAB® simplex²



Bottle-top dispenser with variable volume.

Included in delivery: VITLAB® simplex² (GL 45 thread), 3 respectively 5 thread adapters* made of PP, telescopic filling tube, mounting tool, quality certificate and operating manual.

Volume	Graduation	A**	A**	CV**	CV**	PU	Cat. No.
ml	ml	≤ ± %	≤ ± µl	≤ %	≤ µl		
0.2 - 2.0	0.05	0.5	10	0.1	2	1	1621503
0.5 - 5.0	0.10	0.5	25	0.1	5	1	1621504
1.0 - 10.0	0.20	0.5	50	0.1	10	1	1621505
2.5 - 25.0	0.50	0.5	125	0.1	25	1	1621506
5.0 - 50.0	1.00	0.5	250	0.1	50	1	1621507
10.0 - 100	0.0 1.00	0.5	500	0.1	100	1	1621508

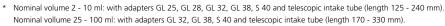

^{**} Error tolerance conforming to DIN EN ISO 8655-5, related to the nominal (maximum) volume marked on the device where the device, environment and distilled H₂O are at the same temperature (20 °C). Checks are done in accordance with DIN EN ISO 8655-6 with the device filled to capacity and with uniform and jolt-free dosing.



VITLAB® simplex_{fix}



Bottle-top dispenser with fixed volume.

Included in delivery: VITLAB® simplex² (GL 45 thread), 5 thread adapters* made of PP, telescopic filling tube, mounting tool, quality certificate and operating manual.

Volume ml	Graduation ml	A** ≤ ± %	A** ≤ ± µl	CV** ≤ %	CV** ≤ µl	PU	Cat. No.
1.0	-	1.0	10	0.2	2	1	1622502
5.0	-	0.5	25	0.1	5	1	1622504
10.0	-	0.5	50	0.1	10	1	1622505

^{*} Nominal volume 1 - 10 ml: with adapters GL 25, GL 28, GL 32, GL 38, S 40 and telescopic intake tube (length 125 - 240 mm).

^{**} Error tolerance conforming to DIN EN ISO 8655-5, related to the nominal (maximum) volume marked on the device where the device, environment and distilled H₂O are at the same temperature (20 °C). Checks are done in accordance with DIN EN ISO 8655-6 with the device filled to capacity and with uniform and jolt-free dosing.





VITLAB® TA²

The VITLAB® TA² dispenser is the dosing device of choice to meet the demanding purity standards required in trace analysis. The high quality parts that come exclusively in contact with the medium and the specially developed and proven cleaning process to be done before use results in a reduced release of trace metal ions to the low ppb range, or, depending on the application, even the ppt range. The parts that are in contact with media are made of various fluoroplastics (e.g. ETFE, FEP, PFA, PTFE, PCTFE), Al₂O₃-sapphire, platinum-iridium or tantalum (depending on model).

Thanks to the excellent chemical resistance of the materials used, the new dispenser can also be deployed with **highly concentrated acids and bases**, such as perchloric, sulphuric and nitric acid. Depending on the application, there is a choice of two different valve spring systems: the VITLAB® ${\rm TA}^2$ with tantalum spring is recommended for dosing of hydrogen peroxide (${\rm H_2O_2}$). For applications using sodium hydroxide (up to a max. concentration of 30%) or hydrogen fluoride (HF) the platinum-iridium spring is recommended. In order to minimize the loss of valuable reagents or sample solutions, VITLAB offers the dispenser with the optional recirculation valve. The practical screwing mechanism and inner toothed bar enable a fast, easy and precise volume adjustment. Also available with DAkkS calibration certificate (at additional cost).

Included in delivery:

VITLAB® TA² dispenser (thread GL 45), 3 bottle adapters (GL 28/S 28 (ETFE), GL 32 (ETFE) and S 40 (PTFE), telescopic filling tube, recirculation valve (optionally), mounting tool, quality certificate and operating manual.

Volume ml	Valve spring	Recircu- lation	Graduation ml	A* ≤ ± %	CV* ≤ %	PU	Cat. No.
1.0 - 10.0	Pt-lr	no	0.2	0.5	0.1	1	1627515
1.0 - 10.0	Pt-lr	yes	0.2	0.5	0.1	1	1627525
1.0 - 10.0	Та	no	0.2	0.5	0.1	1	1627535
1.0 - 10.0	Та	yes	0.2	0.5	0.1	1	1627545

Error tolerance conforming to DIN EN ISO 8655-5, related to the nominal (maximum) volume marked on the device where the device, environment and distilled $\rm H_2O$ are at the same temperature (20 °C). Checks are done in accordance with DIN EN ISO 8655-6 with the device filled to capacity and with uniform and jolt-free dosing.

Recommended dispensing media for VITLAB® TA²

Dispensing medium	Valve spring: Pt-Ir	Valve spring: Ta
Acetic acid	+	+
Ammonia solution	+	+
Bromine	+	+
Hydrochloric acid	+	+
Hydrofluoric acid*	+	-
Hydrogen peroxide	-	+
Nitric acid	+	+
Perchloric acid	+	+
Phosphoric acid	+	+
Sodium hydroxide, 30%	+	-
Sulphuric acid	+	+
Water	+	+

⁺ suitable / - unsuitable

^{*} Note: Hydrofluoric acid reacts slightly with sapphire resulting in slightly increased aluminium levels. To reduce these values we recommend discarding 3-5 dosings of 2 ml each before performing analysis.

Brown glass bottles for VITLAB® genius² and simplex²

Threaded brown glass (soda lime glass) bottles with an ethylene acrylate coating for increased safety, and a screw cap. The plastic coating significantly reduces hazardous glass splintering during breakage. The maximum working temperature for coated bottles is 80 °C. To preserve the coating, do not clean at temperatures exceeding 60 °C.

Volume	Form	Thread	PU	Cat. No.
ml		GL		
250	square	32	1	1671515
500	square	32	1	1671520
1000	square	45	1	1671500
2500	round	45	1	1671510



Suitable plastic bottles can be found starting on page 101.

Plastic stand for VITLAB® dispensers

For secure anchoring, made entirely of polypropylene for contamination-free operation (no metal). Suitable for VITLAB® dispensers with screw coupling GL 45.

Stand rod: 300 mm; base: 220 x 160 mm; weight: 1,130 g.


Description	PU	Cat. No.
Plastic stand	1	1671116



Drying tube for VITLAB® dispensers

PP, transparent, unfilled, with sealing ring (PTFE). Can directly be connected to every dispenser.

Description	PU	Cat. No.
Drying tube, PP, unfilled	1	1671090

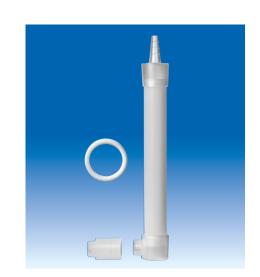



Flexible discharge tubes for VITLAB® dispensers

Coiled, made of FEP, approx. 80 cm length, with handle and recirculation valve made of PTFE. Includes handle and assembly instructions. Not suitable for hydrofluoric acid (HF)!

Description	PU	Cat. No.
Flexible discharge tube for genius ² / simplex ² 2, 5 and 10 ml	1	1678132
Flexible discharge tube for genius ² / simplex ² 25, 50 and 100 ml	1	1678134
Flexible discharge tube for VITLAB® TA² (with grey lever)	1	1678136



Adapter for VITLAB® dispensers

For securely screwing the dispenser onto the reagent bottles with an NS neck, GL screw threading or an S buttress thread. For VITLAB® TA² please use adapter ETFE/PTFE. These are recommended for VITLAB® genius² and simplex² when an increased resistance to chemicals is required (see media table pg. 11).



Description	External thread	Bottle neck threads	PU	Cat. No.
NS-adapter, PP	GL 32	NS 19/26	1	1670066
NS-adapter, PP	GL 32	NS 24/29	1	1670067
NS-adapter, PP	GL 32	NS 29/32	1	1670068
Thread adapter, PP	GL 32	GL 25	1	1670150
Thread adapter, PP	GL 32	GL 28	1	1670155
Thread adapter, PP	GL 45	GL 32	1	1670180
Thread adapter, PP	GL 45	GL 38	1	1670110
Thread adapter, PP	GL 45	S 40	1	1670120
Thread adapter, ETFE	GL 32	GL 25	1	1670072
Thread adapter, ETFE	GL 32	GL 28	1	1670080
Thread adapter, ETFE	GL 32	GL 45	1	1670105
Thread adapter, ETFE	GL 45	GL 32	1	1670100
Thread adapter, ETFE	GL 45	GL 38	1	1670115
Thread adapter, PTFE	GL 45	S 40	1	1670125

Telescopic filling tubes for VITLAB® dispensers

Telescopic filling tube made of FEP, ETFE and PTFE.

Suitable for dispenser with nominal volume	Outer-Ø mm	Length mm	PU	Cat. No.
2/5/10 ml	6.0	70-140	1	1678210
2/5/10 ml	6.0	125-240	1	1678212
2/5/10 ml	6.0	195-350	1	1678214
2/5/10 ml	6.0	250-480	1	1678216
25/50/100 ml	7.6	170-330	1	1678218
25/50/100 ml	7.6	250-480	1	1678220



Sealing ring for the valve block for VITLAB® dispensers

Sealing ring for valve block made of PTFE for dosing highly volatile media.

Description	PU	Cat. No.
Sealing ring for valve block	1	1671683



Ventilation plug for micro filter for VITLAB® genius² and simplex²

PP, with Luer-cone and sealing ring (PTFE).

Description	PU	Cat. No.
Ventilation plug for micro filter	1	1671682



Dispensing cartridge for VITLAB® TA²

Calibrated, including safety ring, with quality certificate. Nominal volume 10 ml.

Description	PU	Cat. No.
Dispensing cartridge	1	1670702





VITLAB® piccolo

For dispensing tiny quantities of liquids in all areas of biochemical and medical research.

Even the **smallest quantities can be dispensed directly from the bottle** with the VITLAB® piccolo - a big help, particularly for serial dispensing operations. Special advantage: Pipette tips are unnecessary. This reduces costs.

The ergonomic design makes dispensing effortless and stress-free. The VITLAB® piccolo **can be operated with only one hand**. Use the thumb to depress the volume dispensing button, just as with a pipette, and a reset mechanism refills the volume automatically. The discharge tube can be rotated over 360° so that it is always optimally situated with

respect to the bottle label.

The VITLAB® bottletop dispensers piccolo 1 and piccolo 2 are used mainly in connection with aqueous and highly diluted agents.

VITLAB® piccolo 1 with a fixed volume VITLAB® piccolo 2 with two fixed volumes

Included in delivery:

VITLAB® piccolo 1 or 2 (GL 28 threads), mounting tool, quality certificate and operating manual.

Туре	Volume	A*	CV*	PU	Cat. No.
	μl	≤ ± %	≤ %		
piccolo 1	100	3.0	0.4	1	1610501
piccolo 1	200	2.5	0.4	1	1610502
piccolo 1	250	2.0	0.4	1	1610503
piccolo 1	500	1.5	0.3	1	1610504
piccolo 1	1000	1.0	0.2	1	1610506
piccolo 2	100 / 250	2.0	0.4	1	1611503
piccolo 2	500 / 1000	1.0	0.2	1	1611506
piccolo 2	1000 / 2000	1.0	0.2	1	1611508

^{*} Error tolerance conforming to DIN EN ISO 8655-5, related to the nominal (maximum) volume marked on the device where the device, environment and distilled ${\rm H_2O}$ are at the same temperature (20 °C). Checks are done in accordance with DIN EN ISO 8655-6 with the device filled to capacity and with uniform and jolt-free dosing.

Adapter for VITLAB® piccolo

For securely screwing the dispenser onto reagent bottles with GL screw threading.

Description	External thread	Bottle neck threads	PU	Cat. No.
Thread adapter, PP	GL 28	GL 32	1	1670145



Bottles for VITLAB® piccolo, PE-HD

Transparent. With screw cap made of PP. Space-saving due to the square cross-section and the high shoulders.

Volume ml	Thread	Height mm	Dimension mm	PU	Cat. No.
100	GL 32	78	46 x 46	24	92489
250	GL 28	80	80 x 80	24	91989
500	GL 32	106	90 x 90	12	92089
1000	GL 32	187	80 x 80	12	92189



Brown glass bottles for VITLAB® piccolo

Threaded brown glass (soda lime glass) bottles with an ethylene acrylate coating for increased safety, and a screw cap. The plastic coating significantly reduces hazardous glass splintering during breakage. The maximum working temperature for coated bottles is 80 °C. To preserve the coating, do not clean at temperatures exceeding 60 °C.

Volume ml	Form	Thread GL	PU	Cat. No.
100	round	28	1	1671505
100	square	32	1	1671506



Perfection in Liquid Handling PRECISE AND CONVENIENT PIPETTING





VITLAB® micropipette



The VITLAB® piston-operated pipettes are the ideal manual pipettes for demanding laboratory applications, and have all the features required by users: robust, with ergonomic shape and simple operation, completely autoclavable, highly accurate with simple calibration for long-lasting reliability.

The large, central pipetting button provides a uniform and smooth movement of the piston. For rapid replacement of the tips, the ergonomic eject button is placed easily accessible to the thumb on the front side. The VITLAB® micropipette is easy to use for both right-and left-handers. The 4-digit volume display with integrated zoom function and vertical arrangement of the numbers (top to bottom reading direction) ensures an **optimal readability of the volume** at all times. The desired volume can be set by rotating the volume-setting wheel with ease and precision. The clearly visible colour-coded frame of the volume display allows easy selection of the right pipette tip.

If necessary, e.g. for applications with non-aqueous solutions, the **integrated calibration function allows an adjustment without tools directly in the laboratory**. The corrosion-resistant piston and ejector ensure a long product life.

The micropipette is completely autoclavable at $121 \, ^{\circ}\text{C}$ (2 bar) according to DIN EN 285. Also available with DAkkS calibration certificate (at additional cost).

Included in delivery: VITLAB® micropipette, silicone grease, sample bag with pipette tips, quality certificate, and operating manual.

Volume	A*	A*	CV*	CV*	Tip	PU	Cat. No.
μl	≤ ± %	≤ ± µl	≤ %	≤ µl	μl		
0.5 - 10	1.0	0.1	0.5	0.05	20	1	1641000
2 - 20	8.0	0.16	0.4	0.08	200	1	1641002
10 - 100	0.6	0.6	0.2	0.2	200/300	1	1641004
20 - 200	0.6	1.2	0.2	0.4	200/300	1	1641006
100 - 1000	0.6	6	0.2	2	1000	1	1641008
500 - 5000	0.6	30	0.2	10	5000	1	1641010
1000 - 10000	0.6	60	0.2	20	10000	1	1641012

^{*} Calibrated to deliver ,Ex'. * Accuracy and coefficient of variation based on the nominal volume (= maximum volume) printed on the instrument, if instrument, environment and distilled water are at the same temperature (20 °C), as well as uniform, jerk-free handling. The margins of error are under those specified in DIN EN ISO 8655-2.

Pipetting

VITLAB® micropipette Starter-Sets

Each VITLAB® Starter Set includes 3 variable VITLAB® micropipettes with different volumes and associated, color-coded tip boxes, as well as 3 rack mounts for appropriate storage of your new VITLAB® micropipettes.

Our micropipettes are completely autoclavable at 121 °C (2 bar) according to DIN EN 285.



Starter Set "Mini"

Scope of delivery:

- VITLAB® micropipette 0.5 10 μ l
- VITLAB® micropipette 10 100 μl
- VITLAB® micropipette 100 1000 μl
- Tip-Box 0.5 20 μl
- Tip-Box 2 200 μl
- Tip-Box 50 1000 μl
- Rack mount (3x)

Cat. No.: 33331

Starter Set "Classic"

Scope of delivery:

- VITLAB® micropipette 2 20 μl
- VITLAB[®] micropipette 20 200 μl
- VITLAB® micropipette 100 1000 μl
- Tip-Box 2 200 µl (2x)
- Tip-Box 50 1000 μl
- Rack mount (3x)

Cat. No.: 33332

Starter Set "Maxi"

Scope of delivery:

- VITLAB® micropipette 100 1000 μl
- VITLAB® micropipette 500 5000 μl
- VITLAB® micropipette 1000 10000 μl
- Tip-Box 50 1000 μl
- Tip-Box 0.5 5 ml
- Tip-Box 1 10 ml
- Rack mount (3x)

Cat. No.: 33333




Accessories for VITLAB® micropipettes

With the practical rack mount and freely rotatable bench-top stand, VITLAB® micropipettes can be stored safely and ready to use.

Description	PU	Cat. No.
Wall mount for 1 pipette	1	1672000
Bench-top stand for 6 pipettes	1	1672002
Filter for pipette, 5 ml	25	1672010
Filter for pipette, 10 ml	25	1672012
Silicone grease for pipettes, up to 1000 µl	1	1672015
Silicone grease for pipettes 5 ml / 10 ml	1	1672016

Pipette tip selection guide Which pipette tip will fit my VITLAB® micropipette?

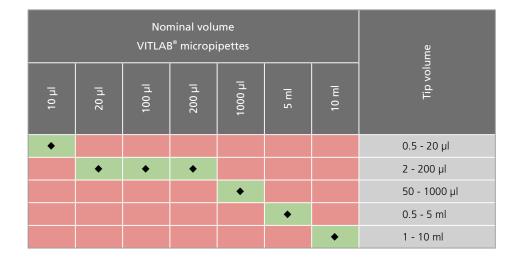



VITLAB® pipette tips

VITLAB® pipette tips are made from high-quality polypropylene and are autoclavable at 121 °C (2 bar) according to DIN EN 285. The raw material used is free from additives such as DiHEMDA (di(2-hydroxyethyl)methyldodecyl-ammonium) and oleamide (9-octadecenamide) that often cause interference, particularly in biological labs. All palletized pipette tips up to 1000 μ l are free of DNA (< 40 fg), RNase (< 8.6 fg), endotoxins (< 1 pg) and ATP (< 1 fg). The pipette tips are CE marked according to the IVD guideline 98/79 EC and optimally suited for VITLAB® micropipettes.

Furthermore, the tips are compatible with most pipette models from BRAND, GILSON®, Thermo Fisher Scientific FINNPIPETTE®, Eppendorf® and sartorius® Biohit®. The 5 ml tip is only tested for VITLAB, BRAND and Thermo Fisher Scientific FINNPIPETTE®. The 10 ml tip is only tested for VITLAB, BRAND and Eppendorf®. Note: Pipette shafts are subject to modification and should be checked before use. The fit depends on the manufacturer, pipette type, serial number, and date of manufacture, among other things.

Packaging variations

VITLAB® pipette tips are available in the variations palletized in the Tip-Box and packed in bags. Additionally, empty Tip-Boxes for self-filling are available. Tip-Boxes up to 1000 μ l are stackable and their format conforms to the common 96 (8x12) unit format.



Resealable bags

All tips up to 1000 µl are produced under cleanroom conditions, automatically shrink-wrapped in reclosable bags and packaged in cartons. The article number, volume range and lot number of the tips are printed on every bag.



Tip-Box (up to 1000 µl)

PP box with functional hinged and snap-on lid. For all volume ranges up to 1000 µl in practical 8x12 format. Stackable and autoclavable at 121 °C according to DIN EN 285.



Tip-Box 5/10 ml

PP box with fitted lid. Filled with 5 ml (28 pcs.) or 10 ml (18 pcs.) tips. The box is autoclavable at 121 °C according to DIN EN 285.



Pipette tips, $0.5 - 20 \mu l$



PP, non-sterile with graduation at 2 and 10 μ l. Length: 46 mm. Slim tip for contact-free pipetting into microtiter plates. Tip-Box with gray mounting plate for easy identification; palletized tips are colourless.

Variation	Packaging	PU	Cat. No.
Bag, Standard	2 bags with 1000 tips	2000	148894
Bag, Maxi	10 bags with 1000 tips	10000	155494
Tip-Box, filled	Box with 96 tips on gray mounting plate	5	149794
Tip-Box, empty	Box with gray mounting plate, without tips	1	155400



Pipette tips, $2 - 200 \mu l$



PP, non-sterile with graduation at 20 and 100 μ l. Length: 50 mm. Tip-Box with yellow mounting plate for easy identification; palletized tips are colourless. Tips in bags are coloured yellow.

Variation	Packaging	PU	Cat. No.
Bag, Standard	1 bag with 1000 tips	1000	148994
Bag, Maxi	10 bags with 1000 tips	10000	155694
Tip-Box, filled	Box with 96 tips on yellow mounting plate	5	149994
Tip-Box, empty	Box with yellow mounting plate, without tips	1	155600

Pipette tips, 50 - 1000 μl



PP, non-sterile with graduation at 250, 500 and 1000 μ l. Length: 70 mm. Tip-Box with blue mounting plate; palletized tips are colourless. Tips in bags are coloured blue.

Variation	Packaging	PU	Cat. No.
Bag, Standard	2 bags with 500 tips	1000	149194
Bag, Maxi	10 bags with 500 tips	5000	155994
Tip-Box, filled	Box with 96 tips on blue mounting plate	5	150194
Tip-Box, empty	Box with blue mounting plate, without tips	1	155900



Pipette tips, 0.5 - 5 ml



PP, non-sterile. Length: 160 mm. Diameter: approx. 9.6 mm. Slim shape for pipetting into narrow vessels such as measuring flasks with NS 12/21.

Variation	Packaging	PU	Cat. No.
Bag, Standard	1 bag with 200 tips	200	146294
Tip-Box, filled	Box with 28 tips	1	150294

Pipette tips, 1 -10 ml



PP, non-sterile. Length: 156.5 mm. Diameter: approx. 15 mm.

Variation	Packaging	PU	Cat. No.
Bag, Standard	2 bags with 100 tips	200	146494
Tip-Box, filled	Box with 18 tips	1	150394

Perfection in Liquid Handling

RAPID AND ACCURATE TITRATION











VITLAB® continuous E/RS

CE

The VITLAB® continuous bottle-top burette (Figure 1) enables continuous titration, which leads to fast, convenient and accurate results. The angled display shows 4-position titration volume in large, easily read numbers (Figure 2), which simplifies operation. Turning the two hand wheels supplies the titration medium in a **continuous and pulse-free** manner via the specially developed double-piston pump (Figure 3). Filling procedures are not necessary. This innovative technology increases safety; its compact design and low centre of gravity reduce risk of falling over, especially with smaller bottles. The height and length of the discharge tube can be adjusted, making it possible to work safely with both short and tall bottles. The innovative recirculation system (Figure 4) **prevents the loss of valuable reagent** and reduces the risk of splashes. With its simple-to-use calibration function, VITLAB® continuous fulfils the corresponding requirements for test equipment monitoring without instrument downtime. Margins of error are under those specified in the DIN EN ISO 8655-3 standard, even for partial volumes. Also available with DAkkS calibration certificate (at additional cost).

Included in delivery:

VITLAB® continuous E/RS (GL 45 thread), 3 PP thread adapters (GL 32, GL 38 and S 40), telescopic filling tube (200 - 350 mm), telescopic discharge tube (140 - 220 mm), two microbatteries 1.5 V (LR 03/AAA), quality certificate and operating manual.

Туре	Nominal volume ml	A* ≤ ± %	CV* ≤ %	PU	Cat. No.
Е	25	0.2 at 25 ml	0.1 at 25 ml	1	1620506
RS	50	0.2 at 50 ml	0.1 at 50 ml	1	1620507

^{*} Error tolerance conforming to DIN EN ISO 8655-3, related to the nominal (maximum) volume marked on the device where the device, environment and distilled H₂O are at the same temperature (20 °C). Checks are done in accordance with DIN EN ISO 8655-6 with the device filled to capacity and with uniform and jolt-free dosing.

The VITLAB® continuous E/RS bottle-top burette can be used for the following titrants up to a concentration of **1 mol/L**:

Acetic acid	Potassium hydroxide
Ammonium iron (II) sulphate solution	Potassium iodate solution
Ammonium thiocyanate solution	Potassium permanganate solution
Barium chloride solution	Potassium thiocyanate solution
Bromide bromate solution	Silver nitrate solution
Cerium (IV) sulphate solution	Sodium arsenite solution
EDTA solution	Sodium carbonate solution
Hydrochloric acid	Sodium chloride solution
Iron (II) sulphate solution	Sodium hydroxide
Nitric acid	Sodium nitrite solution
Oxalic acid solution	Sodium thiosulphate solution
Perchloric acid	Sulphuric acid
Potassium bromate solution	Tetra-n-butylammonium hydroxide solution
Potassium bromide / bromate solution	Zinc sulphate solution
Potassium dichromate solution	

The recommendations in this table have been carefully tested and reflect the most current information available. Always follow the instruction manual for the instrument as well as the reagent manufacturer's specifications. Should you require information on chemicals not listed, please do not hesitate to contact us. As at 03/17.

Adapter for VITLAB® continuous E/RS

For secure screwing of the burettes onto reagent bottles with an NS neck, GL screw threading or an S buttress thread.

Description	External thread	Bottle neck threads	PU	Cat. No.
NS-adapter, PP	GL 32	NS 19/26	1	1670066
NS-adapter, PP	GL 32	NS 24/29	1	1670067
NS-adapter, PP	GL 32	NS 29/32	1	1670068
Thread adapter, PP	GL 32	GL 28	1	1670155
Thread adapter, PP	GL 45	GL 32	1	1670180
Thread adapter, PP	GL 45	GL 38	1	1670110
Thread adapter, PP	GL 45	S 40	1	1670120
Thread adapter, ETFE	GL 32	GL 28	1	1670080
Thread adapter, ETFE	GL 45	GL 32	1	1670100
Thread adapter, ETFE	GL 45	GL 38	1	1670115
Thread adapter, PTFE	GL 45	S 40	1	1670125





Drying tube for VITLAB® continuous E/RS

PP, transparent, unfilled. Can be connected directly to the burette.

Description	PU	Cat. No.
Drying tube, PP, unfilled	1	1671095



Telescopic filling tube for VITLAB® continuous E/RS

For the filling of titration medium from bottles of different heights.

Description	Length mm	PU	Cat. No.
Telescopic filling tube, FEP, ETFE, PTFE	200 - 350	1	1671085



Bottles for VITLAB® continuous E/RS

Threaded brown glass (soda lime glass) bottles with an ethylene acrylate coating.

Volume ml	Thread GL	Shape	PU	Cat. No.
1000	45	square	1	1671500
2500	45	round	1	1671510





VITLAB® Dr. Schilling burettes

Burettes made of borosilicate glass 3.3, tolerances according to DIN ISO 384 Class B, with high-contrast black markings. Calibrated to deliver 'Ex'. Automatic zero setting. The burette stopcock turns easily and enables fine titration. The holding device for the riser pipe serves as additional shock-proofing.

The patented VITLAB® symbiotic (DE 10 2005 034 963) offers in addition to the impact protection a thermally stable plastic coating of the glass burette tube. This provides additional breaking resistance and acts as a splinter protection.

Materials: burette of borosilicate glass 3.3, PP filling tube, PMP/PTFE burette stopcock, and PE-LD reservoir bottle.

Volume ml	Tolerance ± ml	Graduation ml	Height mm	Bottle ml	PU	Cat. No.
VITLAB® symbiotic, with Schellbach stripes (blue/white)						
25	0.05	0.05	900	1000	1	106599
50	0.10	0.10	900	1000	1	106699
Burette with Schellbach stripes (blue/white), without plastic-coating						
25	0.05	0.05	900	1000	1	106399
50	0.10	0.10	900	1000	1	106499
Burette made of brown glass, without Schellbach stripes, without plastic-coating						
25	0.05	0.05	900	1000	1	106799
50	0.10	0.10	900	1000	1	106899



Burettes, borosilicate glass 3.3

Plastic-coated burette made of borosilicate glass 3.3, tolerances corresponding to DIN ISO 384, Class B. With Schellbach stripes (blue/white) and high-contrast black marking. Calibrated to deliver 'Ex'. The burette stopcock turns easily and enables fine titration. The temperature-stable plastic coating on the tube provides splinter protection.

Materials: Burette of borosilicate glass 3.3, with PMP/PTFE burette stopcock.

Volume ml	Tolerance ± ml	Graduation ml	Length mm	PU	Cat. No.
25	0.05	0.05	800	2	105599
50	0.10	0.10	800	2	105699

Titration

Burette stopcocks, PMP/PTFE

Stopcocks made of PMP. Plugs of PTFE with polished surfaces, turn easily. Insert with two seals.

Cat. No. 105799: For 25 ml tubes with an inner diameter of 7.75 \pm 0.1 mm.

Cat. No. 105899: For 50 ml tubes with an inner diameter of 11.5 ± 0.1 mm.

For burettes	Tip inner Ø	PU	Cat. No.
ml	mm		
25	1.25	5	105799
50	1.25	5	105899




Burette clamps, PP

Practical holder for anchoring burettes vertically on the support stand. The volume scale remains clearly readable.

Burette clamps with slip-resistant rubber-coated tips and springs made of stainless steel. With stand clamps for securing to rods of \emptyset 8-14 mm.

Type	PU	Cat. No.
For 1 burette	5	80139
For 2 burettes	5	80140

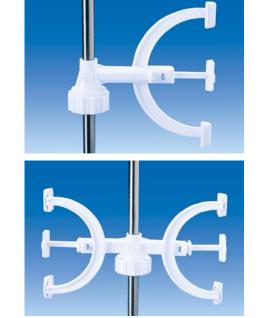

Calibration certificates

For all volumetric instruments that are subject to test equipment monitoring, a written documentation about the regularly calibration resp. volume control is necessary. The documentation should contain the values for accuracy and coefficient of variation as well as the testing procedure and test frequency. A distinction is made between:

- Quality certificates (factory calibration report)
- Official calibration certificates (Bureau of standards, DAkkS)

Quality certificates

The VITLAB quality certificate is a factory calibration report on the basis of the quality assurance system according to DIN EN ISO 9001. Quality certificates are available as a batch or individual certificate. Devices from one production batch have the same lot number as the quality certificate. The certificate records for the specific batch the mean value, standard deviation and day of issue. In the case of an individual certificate*, the volumetric instrument and the certificate bear an individual serial number in addition to the lot number. The certificate records besides the day of issue also the measured volume and the measurement uncertainty.

DAkkS calibration certificate

The DAkkS calibration certificate* documents officially the traceability of measuring results to national and international standards as required by the standards DIN EN ISO 9001 and DIN EN ISO / IEC 17025 for the monitoring of measuring instruments. A major difference between factory calibration services and DAkkS laboratories is the accurate determination of the respective uncertainty of measurement guaranteed by the accredited laboratory and supervised by the DAkkS. DAkkS calibration certificates are appropriate in uses in which calibrations of an accredited laboratory are requested, where high level calibrations are demanded and for calibration of reference standards and instruments for comparative measurements.

Calibration service

VITLAB offers a repair, maintenance and calibration service (incl. DAkkS calibration) for all Liquid Handling devices made by VITLAB. The calibration laboratory accredited by the "Deutsche Akkreditierungsstelle GmbH" (DAkkS) is authorized to issue DAkkS Calibration certificates for the following instruments: Liquid Handling products like VITLAB piston-operated pipettes and burettes, VITLAB dispensers and VITLAB volumetric plastic labware.

* Available at additional cost

Competence in Volume Measurement





The very highes volumetric instr

Volume measurement is a routine laboratory operation. Therefore, volumetric instruments such as volumetric flasks, measuring cylinders and pipettes are standard equipment in any analytical laboratory.

The importance of the standard of measurement accuracy in your routine laboratory operations cannot be overstated. VITLAB has decades of experience in the development and production of laboratory products which are used to measure volumes. VITLAB is the first manufacturer to produce Class A measuring cylinders from PMP that are certified compliant according to DIN 12681.

All Class A PMP volumetric flasks are optionally available in transparent or UV-absorbing variations for light-sensitive substances.

t precision in uments

Calibration

Type "Ex": The delivered quantity of liquid corresponds to the volume printed on the instrument (pipettes and burettes).

Type "In": The contained quantity of liquid corresponds to the volume printed on the instrument (volumetric flasks and measuring cylinders).

VITLAB calibrates each individual volumetric flask "to contain" (In) at a reference temperature of 20 °C. The hydrophobic characteristics of the materials in plastic volumetric instruments lead to the measured volume being the same as the delivered quantity ("In" = "Ex") for aqueous solutions.

Accuracy classes

Class A: The volume tolerances lie within the limits specified by DIN and ISO.

Class B: The volume tolerances are twice the error limits for Class A specified by DIN and ISO. Detailed explanations on "accuracy in volume measurement" are available in the chapter on "General and Technical Information".

Certificate of conformity

The DE-M marking is VITLAB's guarantee that the respective products comply with the German Measurement and Calibration Regulation. The special manufacturing process developed by VITLAB, and the proven VITLAB quality management system, ensure compliance with the volume tolerances specified by the standards.





Volumetric flasks, PFA, Class A, with screw cap, PFA



Highly transparent.

Ring mark individually calibrated to 'In'.

Class A tolerances according to DIN EN ISO 1042.

The PFA screw cap guards against contamination.

Outstanding chemical resistance, can be used with strong oxidants, highly concentrated acids and alkalis, hydrocarbons, and ketones.

With laser-engraved lot number and batch certificate. Thermal stress up to 121 °C (autoclaving) does not permanently exceed the tolerance limit.

To preserve the ring mark, do not clean at temperatures exceeding 60 °C.

Also available with DAkkS calibration certificate or individual quality certificate (at additional cost).

The advantages of PFA:

- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no cross-contamination due to the extremely hydrophobic, anti-adhesive and smooth surfaces
- High thermal stability, from -200 °C to +260 °C
- Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability, suitable for volumetric instruments
- Easy to clean
- Use of high purity raw materials

Further information on PFA can be found starting on page 122.

Volume ml	Tolerance ± ml	Height* mm	Thread GL	PU	Cat. No.
10	0.04	90	18	1	107097
25	0.04	115	18	1	107197
50	0.06	150	18	1	107297
100	0.10	180	18	1	107397
250	0.15	235	25	1	107497
500	0.25	270	25	1	107597
* Height wit	hout screw cap				

Compare: VITLAB® volumetric flasks ...

- ... have a circular, precisely calibrated ring mark with which the meniscus can be read accurately from any position
- ... have a straight neck for precise volume measurement
- ... have a specially formed bottom for superior stability

... are MADE IN GERMANY

Replacement screw caps can be found on page 103.

VITLAB[®] UV-protect volumetric flasks, PMP, Class A with NS stoppers, PP



UV-absorbing, highly transparent. For storage of light-sensitive substances.

With ring mark individually calibrated to 'In'.

Class A tolerances according to DIN EN ISO 1042.

With printed lot number and batch certificate.

Thermal stress up to 121 °C (autoclaving) does not permanently exceed the tolerance limit. To preserve markings, cleaning at no higher than 60 °C is recommended.

Also available with DAkkS calibration certificate or individual quality certificate (at additional cost). More information on VITLAB® UV-protect can be found on page 124.

Volume ml	Tolerance ± ml	Height* mm	Neck NS	PU	Cat. No.
10	0.04	90	10/19	2	670950
25	0.04	115	10/19	2	671950
50	0.06	150	12/21	2	672950
100	0.10	180	14/23	2	673950
250	0.15	235	19/26	2	674950
500	0.25	270	19/26	2	675950
1000	0.40	310	24/29	1	676950
* Height wit	hout stopper				

Volumetric flasks, PMP, Class A with NS stoppers, PP





Highly transparent.

With ring mark individually calibrated to 'In'.

Class A tolerances according to DIN EN ISO 1042.

With printed lot number and batch certificate.

Thermal stress up to 121 °C (autoclaving) does not permanently exceed the tolerance limit. To preserve markings, cleaning at no higher than 60 °C is recommended.

Also available with DAkkS calibration certificate or individual quality certificate (at additional cost).

Volume	Tolerance	Height*	Neck	PU	Cat. No.
ml	± ml	mm	NS		
10	0.04	90	10/19	6	67704
25	0.04	115	10/19	6	67104
50	0.06	150	12/21	6	67204
100	0.10	180	14/23	6	67304
250	0.15	235	19/26	5	67404
500	0.25	270	19/26	4	67504
1000	0.40	310	24/29	3	67604
* Height wit	hout stopper				



Replacement stoppers can be found on page 109.

VITLAB® UV-protect replaces brown glass and is...

- ... substantially lighter in weight ... practically unbreakable
- ... practically impermeable in the UV region
- ... comparable to a light protection factor of 20



Replacement stoppers can be found on page 109.



Volumetric flasks, PMP, Class B with NS stoppers, PP



Highly transparent.

With ring mark individually calibrated to 'In'.

Class B tolerances according to DIN EN ISO 1042.

Thermal stress up to 121 °C (autoclaving) does not permanently exceed the tolerance limit. To preserve markings, cleaning at no higher than 60 °C is recommended.

Volume	Tolerance	Height*	Neck	PU	Cat. No.
ml	± ml	mm	NS		
10	0.08	90	10/19	6	67795
25	0.08	115	10/19	6	67195
50	0.12	150	12/21	6	67295
100	0.20	180	14/23	6	67395
250	0.30	235	19/26	5	67495
500	0.50	270	19/26	4	67595
1000	0.80	310	24/29	3	67695
* Height wit	hout stopper				

Replacement stoppers can be found on page 109.

Volumetric flasks, PMP, Class B with screw caps, PP



Highly transparent.

With ring mark individually calibrated to 'In'.

Class B tolerances according to DIN EN ISO 1042.

Thermal stress up to 121 $^{\circ}$ C (autoclaving) does not permanently exceed the tolerance limit. To preserve markings, cleaning at no higher than 60 $^{\circ}$ C is recommended.



Replacement screw caps can be found on page 105.

Volume	Tolerance	Height*	Thread	PU	Cat. No.
ml	± ml	mm	GL		
10	0.08	90	18	6	677895
25	0.08	115	18	6	671895
50	0.12	150	18	6	672895
100	0.20	180	18	6	673895
250	0.30	235	25	5	674895
500	0.50	270	25	4	675895
1000	0.80	310	32	3	676895
* Height witl	hout screw cap			<u>-</u>	

Volume measurement

Volumetric flasks, PP, Class B with NS stoppers, PP



Highly transparent.

With ring mark individually calibrated to 'In'.

Class B tolerances according to DIN EN ISO 1042.

Thermal stress up to 60 °C does not permanently exceed the tolerance limits.

To preserve markings, cleaning at no higher than 60 °C is recommended.

Volume	Tolerance	Height*	Neck	PU	Cat. No.
ml	± ml	mm	NS		
10	0.08	90	10/19	6	677941
25	0.08	115	10/19	6	671941
50	0.12	150	12/21	6	672941
100	0.20	180	14/23	6	673941
250	0.30	235	19/26	5	674941
500	0.50	270	19/26	4	675941
1000	0.80	310	24/29	3	676941
* Height wit	hout stopper				



Replacement stoppers can be found on page 109.

Volumetric flasks, PP, Class B, with screw cap, PP



Highly transparent.

With ring mark individually calibrated to 'In'.

Class B tolerances according to DIN EN ISO 1042.

Thermal stress up to 60 °C does not permanently exceed the tolerance limits.

To preserve markings, cleaning at no higher than 60 °C is recommended.

Volume	Tolerance	Height*	Thread	PU	Cat. No.
ml	± ml	mm	GL		
10	0.08	90	18	6	677891
25	0.08	115	18	6	671891
50	0.12	150	18	6	672891
100	0.20	180	18	6	673891
250	0.30	235	25	5	674891
500	0.50	270	25	4	675891
1000	0.80	310	32	3	676891
* Height wit	hout screw cap				



Replacement screw caps can be found on page 105.



Graduated cylinders, PMP, Class A, tall form, red printed graduations







Highly transparent. DE-M marked.

With red printed graduations and ring marks at the primary scale points, calibrated 'In'. The lot certificate supplied bears the batch number and the actual nominal value ascertained under the test conditions. The resulting deviations from the nominal value fall well under the allowed tolerances of Class A according to DIN 12681 and ISO 6706. With printed batch number and year of production. Also available with DAkkS calibration certificate or individual quality certificate (at additional cost).

Hexagonal base with bottom studs provides high stability. To preserve markings, do not clean at temperatures exceeding 60 °C. Thus, conditionally autoclavable at 121 °C (2 bar) according to DIN EN 285. For autoclaving we recommend the design with molded graduations (Cat.-No. 64604 - 65304).

Volume	Tolerance	Graduation	Height	Ø	PU	Cat. No.
ml	± ml	ml	mm	mm		
10	0.10	0.20	145	15	2	64614
25	0.25	0.50	170	22	2	64714
50	0.50	1.00	200	27	2	64814
100	0.50	1.00	250	33	2	64914
250	1.00	2.00	315	44	2	65014
500	2.50	5.00	360	58	1	65114
1000	5.00	10.00	440	69	1	65214
2000	10.00	20.00	535	97	1	65414



Graduated cylinders, PMP, Class A, tall form, molded graduations









Highly transparent. DE-M marked.

With molded graduations and ring marks at the primary scale points, calibrated 'In'. The lot certificate supplied bears the batch number and the actual nominal value ascertained under the test conditions. The resulting deviations from the nominal value fall well under the allowed tolerances of Class A according to DIN 12681 and ISO 6706. With the laser engraved batch number and the year of manufacture. Also available with DAkkS calibration certificate or individual quality certificate (at additional cost).

Hexagonal base with bottom studs provides high stability. Thermal stress up to 121 °C (autoclaving) does not cause tolerance limits to be permanently exceeded.

Volume	Tolerance	Graduation	Height	Ø	PU	Cat. No.
ml	± ml	ml	mm	mm		
10	0.10	0.20	145	15	2	64604
25	0.25	0.50	170	22	2	64704
50	0.50	1.00	200	27	2	64804
100	0.50	1.00	250	33	2	64904
250	1.00	2.00	315	44	2	65004
500	2.50	5.00	360	58	1	65104
1000	5.00	10.00	440	69	1	65204
2000	10.00	20.00	482	97	1	65304

Graduated cylinders, PP, Class B, tall form, with molded blue graduations



Highly transparent.

With easily readable, molded, embossed blue graduations and ring marks at the primary scale points. Calibrated 'In'. Class B tolerances according to DIN 12681 / ISO 6706.

Hexagonal base with bottom studs provides high stability. Thermal stress up to 80 $^{\circ}$ C does not cause tolerance limits to be permanently exceeded. To preserve markings, do not clean at temperatures exceeding 60 $^{\circ}$ C.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Tolerance	Graduation	Height	Ø	PU	Cat. No.
ml	± ml	ml	mm	mm		
10	0.20	0.20	145	15	12	646081
25	0.50	0.50	170	22	12	647081
50	1.00	1.00	200	27	12	648081
100	1.00	1.00	250	33	12	649081
250	2.00	2.00	315	44	6	650081
500	5.00	5.00	360	58	6	651081
1000	10.00	10.00	440	69	6	652081
2000	20.00	20.00	482	97	3	653081



Graduated cylinders, PP, Class B tall form, with molded graduations



Highly transparent.

With molded graduations and ring marks at the primary scale points, calibrated 'ln'.

Class B tolerances according to DIN 12681 / ISO 6706.

Hexagonal base with bottom studs provides high stability. Thermal stress up to 80 °C does not cause tolerance limits to be permanently exceeded.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Tolerance	Graduation	Height	Ø	PU	Cat. No.
ml	± ml	ml	mm	mm		
10	0.20	0.20	145	15	12	646941
25	0.50	0.50	170	22	12	647941
50	1.00	1.00	200	27	12	648941
100	1.00	1.00	250	33	12	649941
250	2.00	2.00	315	44	6	650941
500	5.00	5.00	360	58	6	651941
1000	10.00	10.00	440	69	6	652941
2000	20.00	20.00	482	97	3	653941





Graduated cylinders, SAN, Class B tall form, with molded graduations



Crystal clear.

With molded graduations and ring marks at the primary scale points, calibrated 'In'. Class B tolerances according to DIN 12681 / ISO 6706.

Hexagonal base with bottom studs provides high stability. Thermal stress up to 60 °C does not cause tolerance limits to be permanently exceeded.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Tolerance ± ml	Graduation ml	Height mm	Ø mm	PU	Cat. No.
50	1.00	1.00	199	28	12	64891
100	1.00	1.00	260	34	12	64991
250	2.00	2.00	315	47	6	65091
500	5.00	5.00	350	61	6	65191
1000	10.00	10.00	415	76	6	65291



Graduated cylinders, PP, Class B short form, with molded graduations



Highly transparent.

With molded graduations, calibrated 'In'.

Thermal stress up to 80 °C does not cause tolerance limits to be permanently exceeded. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Tolerance ± ml	Graduation ml	Height mm	Ø mm	PU	Cat. No.
25	0.50	0.50	122	22	12	640941
50	1.00	1.00	142	27	12	641941
100	2.00	2.00	163	37	12	642941
250	5.00	5.00	192	51	6	643941
500	10.00	10.00	218	67	6	644941
1000	20.00	20.00	285	78	6	645941

Graduated cylinders, SAN, Class B, short form, with molded graduations



Crystal clear.

With molded graduations, calibrated 'In'.

Thermal stress up to 60 °C does not cause tolerance limits to be permanently exceeded. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Tolerance ± ml	Graduation ml	Height mm	Ø mm	PU	Cat. No.
25	0.50	0.50	122	22	12	64091
50	1.00	1.00	142	27	12	64191
100	2.00	2.00	163	37	12	64291
250	5.00	5.00	192	51	6	64391
500	10.00	10.00	218	67	6	64491
1000	20.00	20.00	285	78	6	64591



Compare: VITLAB® graduated cylinders...

... have guaranteed seamless interiors, which mean the analysis is unaffected by residues and carryover ... have precise calibration ring marks at the primary scale points, with which the meniscus can be read accurately ... a sturdy, even stand for precise volume measurement

... are MADE IN GERMANY

Hydrometer cylinder, PP

Highly transparent, with spout and overflow vessel. For density measurements using a hydrometer. Hydrometer can be read through the overflow vessel with a completely filled cylinder.

With molded graduations and ring marks at the primary scale points, calibrated 'In'. Class B tolerances according to DIN 12681 / ISO 6706.

Hexagonal base with bottom studs provides high stability. Thermal stress up to 80 $^{\circ}$ C does not cause tolerance limits to be permanently exceeded.

Volume ml	Graduation ml	Height mm	Ø mm	PU	Cat. No.
500	5.00	351	73	1	760941





Bulb pipettes, PP, Class B

Calibrated to deliver 'Ex'.

Highly transparent. With high-contrast, blue markings.

Break-resistant.

High chemical resistance.

Exposure to temperatures above 60 °C can lead to volume changes.

Recommended cleaning with mild alkaline detergents up to 60 °C.

Volume	Tolerance	Length	PU	Cat. No.
ml	± ml	mm		
1	0.02	300	12	164094
2	0.02	300	12	164194
5	0.03	300	6	164294
10	0.04	440	6	164394
25	0.05	450	6	164494
50	0.10	460	6	164594



Graduated pipettes, PP, Class B

Calibrated to deliver 'Ex'.

Highly transparent. With high-contrast, blue markings.

Break-resistant.

High chemical resistance.

Outer diameter of suction tube: max. 8 mm.

Exposure to temperatures above 60 °C can lead to volume changes.

Recommended cleaning with mild alkaline detergents up to 60 °C.

Volume	Tolerance	Graduation	Length	PU	Cat. No.
ml	± ml	ml	mm		
1	0.02	0.1	300	12	163094
2	0.02	0.1	300	12	163194
5	0.05	0.1	330	12	163294
10*	0.10	0.1	330	12	163394
10	0.10	0.1	320	12	163594
* Suction tub	ne outer diame	ter: 10 mm			_


Volume measurement

Disposable pipettes, PS, sterile



Crystal clear, graduated, individual sterile packaging, pyrogen-free. Identified by bar-code. With cotton-wool filter.

Volume ml	Graduation ml	Length mm	PU	Cat. No.
1	0.01	272	1000	160110
2	0.01	272	1000	160210
5	0.10	320	500	160510
10	0.10	320	500	161010
25	0.20	345	250	162510

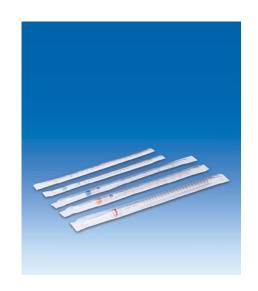

Disposable pipettes, PS, non-sterile

Crystal clear, graduated, non-sterile. Identified by bar-code. With cotton-wool filter.

Volume ml	Graduation ml	Length mm	PU	Cat. No.
1	0.01	272	200	160119
2	0.01	272	200	160219
5	0.10	320	100	160519
10	0.10	320	100	161019



Pipette helpers can be found on pages 48 - 50.



VITLAB pipeo®

CE

For all pipettes from 0.1 to 200 ml.

With the VITLAB pipeo® pipette controller, pipette handling is simple and comfortable. The ergonomic handle - **very light weight** at about 190 grams - and excellent balance all contribute to ease of operation. The speed can be adjusted easily, continuously and exactly with one hand using two buttons. A 50 ml pipette can be filled comfortably in less than ten seconds. The liquid release can be done either by gravity delivery when calibrated 'Ex' (to deliver), or in blow out mode using the battery-operated motor.

Pipettes are held securely and tightly in the exchangeable adapter. Liquid vapours are purged directly to protect the instrument.

One full charge of the nickel-metal hydride battery allows 8 hours of non-stop pipetting. The charge level of the recyclable battery is shown by the LED indicator. Defective batteries are easily replaced. To avoid surprises, the LED light changes from green to red two hours before the battery must be recharged. The VITLAB pipeo® can still be operated while the battery is being recharged.

Included in delivery:

VITLAB pipeo $^{\circ}$, battery charger (100 - 240 V, 50/60 Hz), four plug adapters (EU, UK, US/J, AUS), battery, battery compartment cover, two replacement 0.2 μ m membrane filters and operating manual.

Description	PU	Cat. No.
pipeo®	1	1631500

VITLAB maneus®



The VITLAB maneus® Pipette Helper enables both left- and right-handers to operate all current volumetric and graduated pipettes from 0.1 to 200 ml easily and fatigue-free. Its safe and easy handling allows even inexperienced users **to adjust the meniscus precisely**.

With the design, unscrewing the adapter enables easy replacement of the hydrophobic membrane filter, which **protects the instrument against fluid penetration**.

The valve system is optimised so that liquids can be drawn up simply, without exerting pressure. The highly sensitive filling and discharge of liquids are controlled gently by the pipetting knob. Thus, the suction element provides rapid filling of the pipette (capacity: 50 ml in less than 10 seconds). The discharge bellows are used for the emptying (blow-out) of the pipette. The specially moulded intake cone ensures secure seating for all normal bulb and graduated pipettes (0.1 to 200 ml).

The VITLAB maneus® is simple to dismantle, easy to clean, and completely autoclavable at 121 °C (2 bar) according to DIN EN 285.

Included in delivery: VITLAB maneus $^{\circ}$, replacement 3 μm membrane filter und operating manual.

Description	PU	Cat. No.
maneus [®]	1	1630500



Accessories for VITLAB pipeo® & maneus®

Description	PU	Cat. No.
Membrane filter, 0.2 μm, sterile, VITLAB pipeo®	11	1670647
Membrane filter, 0.2 μm, non-sterile, VITLAB pipeo®	10	1670648
Membrane filter, 3 μm, non-sterile, VITLAB pipeo®, VITLAB maneus®	10	1670650
Wall rack, VITLAB pipeo®	1	1670660





Pipette fillers, NR

Classic accessory for pipetting with volumetric or measuring pipettes. With 3 valves. Valve A: Air release, Valve S: Liquid filling, Valve E: Liquid dispensing.

Туре	PU	Cat. No.
Universal model, for pipettes up to 10 ml	1	104099
Universal model, for pipettes up to 100 ml	1	104199



Pipette fillers

For pipetting liquids, fit all glass and plastic pipettes. Slow rotation of the actuator-wheels draws liquid into the pipette. Pressing the air bleed valve automatically empties the pipette without returning the piston.

For pipettes ml	Colour	PU	Cat. No.
2	Blue	10	324594
10	Green	10	324694
25	Red	10	324794



Pipette stand, PP

Upper portion with 94 bore holes of different diameters for secure placement of volumetric and measuring pipettes of any size.

The stable base has a rotatable, ribbed base plate in which the pipette tips can be gently seated.

The racks are supplied unassembled, and can easily be assembled according to the accompanying assembly instructions.

Ø	Height	PU	Cat. No.
mm	mm		
230	470	2	79194

Volume measurement

Pipette washer, PE-HD

For simple and basic cleaning of pipettes. With discharge siphon for an automatic water exchange.

The complete washing system includes the pipette washer, pipette jar (for pre-cleaning) and pipette basket (for dipping pipettes into the pipette washer or pipette jar). Pipette jars and pipette baskets need to be ordered separately.

Suitable for the use with pipette baskets (cat. nos. 80219 and 80222).

Ø	Height	Effective length	PU	Cat. No.
mm	mm	mm		
170	734	600	1	80217
170	990	840	1	80215



Pipette jars, PE-HD

For pre-cleaning pipettes in detergent solutions.

Suitable for the use with pipette baskets (cat. nos. 80219 and 80222).

Ø mm	Height mm	PU	Cat. No.
162	503	1	80221
162	650	1	80218



Pipette baskets, PE-HD

For dipping pipettes into the pipette jar or pipette-washer and for transferring pipettes. Basket height 300 mm.

With the extension piece, the total height of the pipette basket (cat. no. 80219) increases from 650 to 870 mm.

Description	Ø	Overall height	PU	Cat. No.
	mm	mm		
Pipette basket	145	648	1	80219
Pipette basket	145	497	1	80222
Extension piece for the handle (pipette basket 80219) 2				81219



Volume measurement

Competence in Plastic Labware

MEASURING AND TRANSFERRING





Graduated beakers, PP, molded blue graduations



Highly transparent. With easily readable, molded, embossed blue graduations, and stable, easy-grip handle. To preserve markings, do not clean at temperatures exceeding 60 °C. Conditionally autoclavable at 121 °C (2 bar) according to DIN EN 285. For autoclaving, we recommend the design with molded graduations (cat. nos. 440941 - 447941). Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Divisions	Height	Ø	PU	Cat. No.
ml	ml	mm	mm		
50	2	70	40	24	446081
100	2	80	50	24	447081
250	5	120	74	12	440081
500	10	140	92	12	441081
1000	10	181	117	6	442081
2000	20	213	152	6	443081
3000	50	242	172	6	444081
5000	50	270	204	6	445081



Graduated beakers, PP, molded graduations





Highly transparent. With molded graduations and stable, easy-grip handle. Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Divisions	Height	Ø	PU	Cat. No.
ml	ml	mm	mm		
50	2	70	40	24	446941
100	2	80	50	24	447941
250	5	120	74	12	440941
500	10	140	92	12	441941
1000	10	181	117	6	442941
2000	20	213	152	6	443941
3000	50	242	172	6	444941
5000	50	270	204	6	445941

Measuring and transferring

Graduated beakers, PP, nesting



Highly transparent. With stable handle and easily readable, printed black graduations on both sides. Therefore, the volume is equally visible for left- and right handers. With recess in the handle for better water drainage in the dishwasher. To preserve markings, cleaning at no higher than $60\,^{\circ}\text{C}$ is recommended. For autoclaving we recommend the design with molded graduations (cat. nos. 440941 - 447941).

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Divisions ml	Height mm	Ø mm	PU	Cat. No.
250	5	115	75	12	480941
500	10	140	100	12	481941
1000	10	167	125	12	482941
2000	20	212	148	12	483941
3000	50	242	170	12	484941



Graduated beakers, PP, nesting, coloured



Graduated beakers in four different colours. Transparent, with stable handle and easily readable, printed graduations on both sides. Therefore, the volume is equally visible for left- and right handers. With recess in the handle for better water drainage in the dishwasher. To preserve markings, cleaning at no higher than 60 °C is recommended. For autoclaving we recommend the design with molded graduations (cat. nos. 440941 - 447941). Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Colour [Divisions	Height	Ø	PU	Cat. No.
ml		ml	mm	mm		
500	blue	10	140	100	12	481942
500	yellow	10	140	100	12	481943
500	red	10	140	100	12	481944
500	green	10	140	100	12	481945
500	Set: blue, yellow, red, green (1 item e	each) 10	140	100	1	4811111
1000	blue	10	167	125	12	482942
1000	yellow	10	167	125	12	482943
1000	red	10	167	125	12	482944
1000	green	10	167	125	12	482945
1000	Set: blue, yellow, red, green (1 item e	each) 10	167	125	1	4821111





Graduated beakers, SAN



Crystal clear.

With molded graduations and stable, easy-grip handle.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Divisions ml	Height mm	Ø mm	PU	Cat. No.
250	5	120	70	12	44091
500	10	133	91	12	44191
1000	10	170	116	6	44291
2000	20	215	150	6	44391
3000	50	242	170	6	44491



Collectors, PP or SAN



With molded graduations. Volume: 2000 ml. With stable, easy-grip handle and white PC lid.

Diameter: 150 mm; height: 220 mm.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Description	Divisions ml	PU	Cat. No.
SAN, molded graduations (Picture 1)	20	6	97891
PP, molded graduations	20	6	978941
PP, molded, blue embossed graduations (Pic	ture 2) 20	6	978081
Accessories for collectors			
Lid, PC		6	97791



Measuring and transferring

Buckets, PE-HD



White. Without spout. With division into 1 liter segments.

Stable handle with reinforcement in the middle for comfortable carrying.

Tightly closing, transparent snap-on lid made of PE-LD - please order separately.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Description	Volume I	Divisions I	Height mm	Ø mm	PU	Cat. No.
Bucket	5	1	240	250	1	96093
Bucket	10	1	300	290	1	96393
Lid	for 5 L				1	96293
Lid	for 10 L				1	96593



Buckets with spout, PP



Transparent. With division into 1 liter segments.

With stable handle and spout for easy emptying.

Highly resistant to chemicals.

Without lid.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Divisions	Height	Ø	PU	Cat. No.
1	1	mm	mm		
12	1	330	310	1	96694
15	1	370	310	1	96794





Measuring scoops, PP



White. Also suitable as weighing scoops. With precision formed filling edge and comfortable, stable handle. Easily readable volume quantities on the upper side of the handle.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Length	PU	Cat. No.
ml	mm		
2	60	12	39194
5	82	12	39294
10	100	12	39394
25	135	12	39494
50	160	12	39594
100	200	12	39694
250	260	6	39794
500	315	6	39894
1000	385	6	39994



Measuring scoops, PP, coloured



Measuring scoops in different colours. Also suitable as weighing scoops. With precision formed filling edge and comfortable, stable handle. Easily readable volume quantities on the upper side of the handle.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Cat. No.
395940
395950
396940
396943
396944
396946
396950
396952
396955
396956
397940
397950
3961111

Measuring and transferring

Scoops, PE-HD



Conical in shape with tapered filling edge.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Length	Colour	PU	Cat. No.
ml	mm			
15	115	natural	12	40093
25	135	natural	12	40193
65	185	natural	12	40293
110	215	natural	12	40393
150	250	natural	12	40493
350	310	natural	6	40593
750	350	natural	6	40693
750	350	ultramarine	6	406950
750	350	black	6	406944
1250	400	natural	6	40793
1250	400	ultramarine	6	407950
1250	400	black	6	407944






Spatula, PA

Glass-fibre reinforced. Double spatula or spatula spoon, with stable, easy to hold handle in the middle.

Description	Length mm	PU	Cat. No.
Double spatula	150	10	80594
Double spatula	180	10	80595
Spatula spoon	180	10	80596
Spatula spoon	210	10	80593

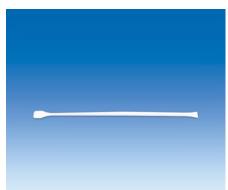

Stirring rod, PP

Spatula-shaped extension for effective manual stirring of small volumes.

Length	PU	Cat. No.
mm		
245	10	80828



Forceps, POM

Yellow, blunt, elastic, very good resilience. With grooves on the outside for optimum handling and grip.

Length	PU	Cat. No.
mm		
Forceps without grooves on the inside of the tip		
115	50	68099
145	50	68199
Forceps with grooves on the inside of the tip		
180	25	68299
250	25	68399



Forceps, PMP

White, pointed, elastic, very good resilience.

Length mm	PU	Cat. No.
115	10	67895
145	10	67995

Measuring and transferring

Funnels, PP



Transparent. Rapid flow due to a steep 60° angle.

Practical handle with loop for hanging.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Ø	Length	Inner stem Ø	Stem length	PU	Cat. No.
approx. ml	mm	mm	mm	mm		
5	30	45	1.5	25	24	40894
6	30	47	4	25	24	41094
14	40	65	4	35	24	41194
32	50	85	7	43	24	41294
88	75	108	7.2	55	12	41394
222	100	155	8	77	12	41494
342	120	180	11	90	12	41594
817	150	220	15	95	12	41694



Powder funnels, PP



Transparent. With short, wide stem and practical tab for hanging. For transfer of powdered and granular substances. Rapid flow due to a steep 60° angle.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Ø mm	Length mm	Inner stem Ø mm	Stem length mm	PU	Cat. No.
65	70	15.5	26	10	70794
80	75	21	26	10	70894
100	92	24	23	10	70994
120	105	27.5	22	10	71094
150*	138	28	22	5	71194
* Without	tab				





Large funnels, PP



Transparent. Rapid flow due to a steep 60° angle. Practical handle for hanging. (Size 12500 ml without handle.) Suitable for filling large amounts of liquids. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011. Optional accessories available: Stainless steel and aluminium sieve insert; however, not permissible for use with foodstuffs.

Volume	Ø	Length	Inner stem Ø	PU	Cat. No.
approx. ml	mm	mm	mm		
1300	203	252	22	6	41794
3200	260	305	23	6	41894
12500	350	440	35	1	41994
Sieve insert Ø: 5	0 mm, for fur	nels no. 41794,	41894	1	42099



Large funnels, PE-HD

Transparent. Rapid flow due to a steep 60° angle. Practical handle for hanging. Suitable for filling large amounts of liquids.

Volume approx. ml	Ø mm	Length mm	Inner stem Ø mm	PU	Cat. No.
12500	400	365	42	1	42294
17500	430	420	37	1	42393



Standard joint funnels, PP



Transparent. For multi-neck flasks, laterally flattened, suitable for standard joint necks of various sizes. Suitable for the filling of liquid or powdered reagents into a reaction flask, especially for loading of multi-neck flasks during a reaction.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

NS	Length mm	Wide opening mm	Stem length mm	PU	Cat. No.
14/23	75	40	17	10	70494
19/26	95	50	23	10	70594
29/32	135	75	30	5	70694


Competence in Plastic Labware



VITsafe[™] - the Safety Wash Bottle


Working with chemical substances, which are sometimes dangerous, requires a high level of responsibility and concentration. With the VITsafe™ safety wash bottles VITLAB offers laboratory equipment that provides safety in the highest degree.

VENT-CAP virtually prevents leakage

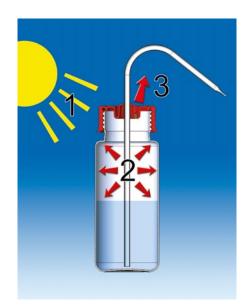

Temperature changes (1) in a laboratory can often cause conventional wash bottles to leak or drip due to the resulting increase in internal pressure (2). The patent-registered, metal-free VENT-CAP screw closure of the VITsafe™ safety wash bottle virtually prevents this. The expansion of the gas due to the rise in pressure can escape through an integrated capillary (3), thus dissipating the static over pressure. In addition, the lack of a spray insert allows turbulences to be almost completely avoided. The smooth and finely drawn tip of the spray bottle allows a precise spray jet and optimizes the fluid backflow. Dripping is subsequently almost entirely prevented.

Clear identification due to safety imprint

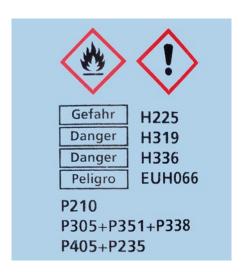

The permanent imprint according to directive (EC) No. 1272/2008 **(GHS)** offers even more safety.

It contains all essential information:

- Substance name in German, English, French and Spanish
- Chemical formula and CAS number
- Hazard pictogram with signal word
- Hazard (H) and precautionary (P) statements
- as well as the U.S.-based NFPA diamond

The VITsafeTM safety wash bottles are available as narrow-mouth or wide-mouth type. The particularly large opening of the wide-mouth bottles allows filling even without a funnel. Select the safety wash bottles to fit your applications from among **17 different substance names** and three volumes (250/500/1000 ml).

Sample preparation

VITsafe™ safety wash bottles, narrow-mouth



Bottle made of PE-LD (or PP for acetone and MEK), PP spray tube. More safety due to the durable safety imprint in accordance to Directive (EC) No. 1272/2008 (GHS), as well as with all important information:

- Material name in German, English, French and Spanish
- Chemical formula, CAS No., hazard pictogram, signal word
- Risk phrases (H phrases), safety phrases (P phrases), as well as NFPA Code

Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube. Practically no leakage or dripping with the VENT-CAP screw cap, the design of which prevents almost all static overpressure.

Imprint	Colour	Volume	Thread	Height*	Ø	PU	Cat. No.
	VENT-CAP	ml	GL	mm	mm		
						_	
Acetic acid	Red	500	25	180	74	6	1332979
Acetone	Red	250	25	135	65	12	1431829
Acetone	Red	500	25	180	74	12	1432829
Acetone	Red	1000	32	215	92	12	1433829
Acetonitrile	Red	500	25	180	74	6	1332969
Distilled Water	Blue	250	25	135	65	12	1331819
Distilled Water	Blue	500	25	180	74	12	1332819
Distilled Water	Red	1000	32	221	92	12	1333819
Ethanol	Orange	250	25	135	65	12	1331869
Ethanol	Orange	500	25	180	74	12	1332869
Ethanol	Red	1000	32	221	92	12	1333869
Ethyl acetate	Red	250	25	135	65	12	1331859
Ethyl acetate	Red	500	25	180	74	12	1332859
Ethyl acetate	Red	1000	32	221	92	12	1333859
Heptane	Red	500	25	180	74	6	1332899
i-Hexane	Red	500	25	180	74	6	1332909
Isopropanol	Yellow	250	25	135	65	12	1331849
Isopropanol	Yellow	500	25	180	74	12	1332849
Isopropanol	Red	1000	32	221	92	12	1333849
Methanol	Green	250	25	135	65	12	1331839
Methanol	Green	500	25	180	74	12	1332839
Methanol	Red	1000	32	221	92	12	1333839
Methyl ethyl ketone (MEK) Red	500	25	180	74	6	1432989
Methylene chloride	Transparent	500	25	180	74	6	1332879
N,N-Dimethylformam	ide Red	500	25	180	74	6	1332889
Pentane	Red	500	25	180	74	6	1433959
Tetrahydrofuran (THF)	Red	500	25	180	74	6	1332939
Toluene	Red	500	25	180	74	6	1332949
Xylene	Red	500	25	180	74	6	1332959
* Height without spra	ay tube						



VENT-CAP wash bottle caps can be found on page 67.



VENT-CAP wash bottle caps can be found on page 67.

VITsafe™ safety wash bottles, wide-mouth



Bottle made of PE-LD (or PP for acetone and MEK), PP spray tube. More safety due to the durable safety imprint in accordance to Directive (EC) No. 1272/2008 (GHS), as well as with all important information:

- Material name in German, English, French and Spanish
- chemical formula, CAS No., hazard pictogram, signal word
- Risk phrases (H phrases), safety phrases (P phrases), as well as NFPA Code

Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube. Practically no leakage or dripping with the VENT-CAP screw cap, the design of which prevents almost all static overpressure.

Imprint	Colour	Volume	Thread	Height*	Ø	PU	Cat. No.
	VENT-CAP	ml	GL	mm	mm		
	5 1	500	45	4.5.5	7.0	-	4252270
Acetic acid	Red	500	45	166	76	6	1352979
Acetone	Red	250	45	146	58	12	1451829
Acetone	Red	500	45	166	76	12	1452829
Acetone	Red	1000	63	226	91	12	1453829
Acetonitrile	Red	500	45	166	76	6	1352969
Distilled Water	Blue	250	45	146	58	12	1351819
Distilled Water	Blue	500	45	166	76	12	1352819
Distilled Water	Red	1000	63	226	91	12	1353819
Ethanol	Orange	250	45	146	58	12	1351869
Ethanol	Orange	500	45	166	76	12	1352869
Ethanol	Red	1000	63	226	91	12	1353869
Ethyl acetate	Red	250	45	146	58	12	1351859
Ethyl acetate	Red	500	45	166	76	12	1352859
Ethyl acetate	Red	1000	63	226	91	12	1353859
Heptane	Red	500	45	166	76	6	1352899
i-Hexane	Red	500	45	166	76	6	1352909
Isopropanol	Yellow	250	45	146	58	12	1351849
Isopropanol	Yellow	500	45	166	76	12	1352849
Isopropanol	Red	1000	63	226	91	12	1353849
Methanol	Green	250	45	146	58	12	1351839
Methanol	Green	500	45	166	76	12	1352839
Methanol	Red	1000	63	226	91	12	1353839
Methyl ethyl ketone (MEK) Red	500	45	166	76	6	1452989
Methylene chloride	Transparent	500	45	166	76	6	1352879
N,N-Dimethylformam	ide Red	500	45	166	76	6	1352889
Pentane	Red	500	45	166	76	6	1453959
Tetrahydrofuran (THF)) Red	500	45	166	76	6	1352939
Toluene	Red	500	45	166	76	6	1352949
Xylene	Red	500	45	166	76	6	1352959
* Height without spr	ay tube						

Other variations available upon request.

Sample preparation

VENT-CAP wash bottle caps, PP


Screw cap and spray tube from PP.

Practically no leakage or dripping with the VENT-CAP screw cap,

the design of which prevents almost all static overpressure.

Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube.

For bottle size ml	Thread GL	Colour	PU	Cat. No.
250	25	Red	12	834111
250	25	Blue	12	834121
250	25	Green	12	834131
250	25	Orange	12	834141
250	25	Yellow	12	834151
250	45	Red	12	834311
250	45	Blue	12	834321
250	45	Green	12	834331
250	45	Orange	12	834341
250	45	Yellow	12	834351
500	25	Transparent	12	834102
500	25	Red	12	834112
500	25	Blue	12	834122
500	25	Green	12	834132
500	25	Orange	12	834142
500	25	Yellow	12	834152
500	45	Transparent	12	834302
500	45	Red	12	834312
500	45	Blue	12	834322
500	45	Green	12	834332
500	45	Orange	12	834342
500	45	Yellow	12	834352
1000	32	Red	12	834213
1000	63	Red	12	834413

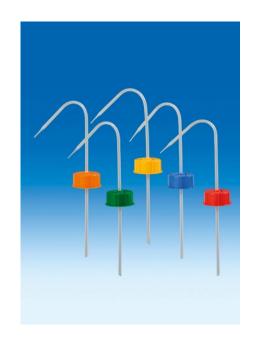



Wash bottles with imprint, PE-LD/PP



Narrow-/wide-mouth bottles made of PE-LD, transparent.

Screw cap and spray tube from PP.

Imprinted with "Distilled Water" in German, English, French and Spanish.

Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube.

Upon request, also with other imprints for non-hazardous media according to the REACh Directive.

Volume ml	Thread GL	Height* mm	Ø mm	PU	Cat. No.
250	25	135	65	12	133181
250	45	146	58	12	135181
500	25	180	74	12	133281
500	45	166	76	12	135281
1000	32	221	92	12	133381
1000	63	226	91	12	135381
* Height with	nout spray tube	<u>.</u>			

Wash bottles made from PE-LD/PP with no imprint, transparent and coloured, can be found on page 69.

Wash-bottles, PP



Narrow-/wide-mouth bottles made of PP, transparent. Screw cap and spray tube from PP. Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube.



Volume	Thread	Height*	Ø	PU	Cat. No.
ml	GL	mm	mm		
250	25	135	65	12	94993
250	45	146	58	12	93793
500	25	180	74	12	95093
500	45	166	76	12	93993
1000	32	215	92	12	95193
1000	63	226	91	12	94193
* Height with	out spray tube				

Height without spray tube

Sample preparation

Wash-bottles, PE-LD/PP



Narrow- / wide-mouth bottles made of PE-LD, transparent. Screw cap and spray tube from PP.

Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube.

Volume ml	Thread GL	Height* mm	Ø mm	PU	Cat. No.
50	18	85	37	24	94588
100	18	114	43	24	94688
250	25	135	65	12	94988
250	45	146	58	12	93788
500	25	180	74	12	95088
500	45	166	76	12	93988
1000	32	221	92	12	95188
1000	63	226	91	12	94188
* Height with	out spray tube				



Wash bottles with imprint can be found on page 68.

Wash-bottles, coloured, PE-LD/PP

Narrow-mouth bottles, made from PE-LD. Available in four different colours to facilitate ready identification. Screw cap and spray tube from PP.

Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube.

Colour	Volume	Thread	Height*	Ø	PU	Cat. No.		
	ml	GL	mm	mm				
red	250	25	135	65	5	132603		
red	500	25	180	74	5	132703		
red	1000	32	221	92	5	132803		
green	250	25	135	65	5	132605		
green	500	25	180	74	5	132705		
green	1000	32	221	92	5	132805		
yellow	250	25	135	65	5	132606		
yellow	500	25	180	74	5	132706		
yellow	1000	32	221	92	5	132806		
blue	250	25	135	65	5	132608		
blue	500	25	180	74	5	132708		
blue	1000	32	221	92	5	132808		
Set: red, green,	500	25	180	74	1	1327111		
yellow, blue (1 ite	m each)							
Set: red, green,	1000	32	221	92	1	1328111		
yellow, blue (1 ite	yellow, blue (1 item each)							
* Height without spray tube								






Wash-bottles, PE-LD

Narrow-mouth bottles, transparent, with screw cap. Spray tube and spray-tube insert made from PE-LD. The classic model, made from soft material with good restoring force.

Volume	Thread	Height*	Ø	PU	Cat. No.
ml	GL	mm	mm		
100	18	106	45	50	134293
250	25	140	59	50	134393
500	25	180	75	50	134493
1000	28	212	94	25	134593
* Height with	out spray tube				

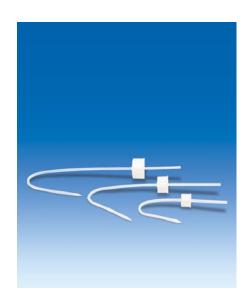

Wash bottle caps, PP

Screw cap and spray tube with drawn-out tip, made of PP. Precise spray jet and optimised medium backflow through the continuous and finely drawn tip of the spray tube.

Thread	PU	Cat. No.
GL		
_18	24	83300
25	12	83301
32	12	83302
45 63	12	83303
63	12	83304

Dropping bottles, PE-LD/PE-HD



Narrow-mouth bottles made of PE-LD, transparent, with dropper insert and screw cap made from PE-HD.

Extra long, fine dropping tip for accurate dispensing.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Thread GL	Height mm	Ø mm	PU	Cat. No.
50	18	129	37	24	94587
100	18	155	43	24	94687
250	25	183	65	12	94987
500	25	228	74	12	95087
1000	32	269	92	12	95187



Caps with dropper inserts, PE-HD



For bottles with GL threads. Cap with dropper insert, complete with screw cap made from PE-HD.

Extra long, fine dropping tip for accurate dispensing.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Thread GL	PU	Cat. No.
18	24	83306
25	12	83307
32	12	83308




Dropping bottles, PE-LD

Narrow-mouth bottles, transparent, with dropper insert and red screw cap made of PE-LD. The classic model, made from soft material with good restoring force.

Volume ml	Thread GL	Height mm	Ø mm	PU	Cat. No.
_20	14	88	31	100	132193
30	14	96	34	100	132293
50	18	115	39	100	132393
100	18	136	45	50	132493
250	25	170	59	50	132593
500	25	209	75	50	132693
1000	28	240	94	25	132793



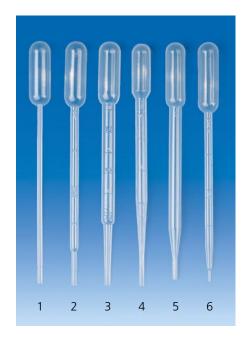

Pasteur pipettes, PE-LD

Disposable. Very good reproducibility of the number of drops per milliliter, thus ideal for distributing aliquots of liquid portions. Pasteur pipettes can be deep-frozen when filled, or if needed, be converted into sealed vessels through heat-sealing. The integrated suction bulb can readily be compressed. Thus, finger fatigue from frequent pipetting is avoided. Can be sterilised with gas or gamma radiation.

Figure No.	Nominal volume ml	Withdraw volume with ball ml	Graduated to ml	Outer Ø tip mm	Length mm	Number of drops per ml	PU	Cat. No.
1		4.4	without	23	148	27-29	5000	148893
2	1	5	1.0	26	150	23-25	5000	148993
3	3	7	3.0	33	150	19-21	5000	149093
4	2	5.9	2	27	153	25	5000	149193
5	4	7	without	32	147	25	5000	149293
6	1	5.5	0.25	12	150	70	5000	149393




Dropping pipettes, PE-LD

With integrated bellows.

For sampling and decanting of infectious or toxic liquids. Graduated.

Volume	Length	PU	Cat. No.
ml	mm		
1.5	134	100	149893
5	195	100	149993

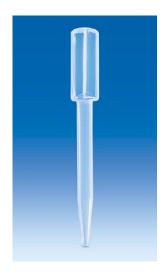

Dropping pipettes, PE-LD

With integrated pipetting bulb.

For sampling and decanting of infectious or toxic liquids.

Without graduations.

Volume ml	Length mm	PU	Cat. No.
1.8	98	250	149693

Sample preparation

Spray bottles

White or transparent bottles made from PP and PE-LD.

Sprayer insert with stable, smoothly operated pump trigger and adjustable spray nozzle, which can be regulated from the finest mist (nebulising) to a precise liquid jet.

Range: approx. 3-4 meters.

For spraying detergents or disinfectants, especially into difficultly accessible areas, as well as applications in thin layer chromatography.

Volume ml	Colour	Material	PU	Cat. No.	
400	white	PP	5	53510	
850	white	PP	5	53610	
1000	transparent	PP	5	95286	
1000	transparent with imprint "Ethanol"*	PE-LD	5	952861	
* More information on imprint on page 64					









Griffin beakers, PFA





Transparent. With molded graduations. Excellent chemical resistance and very high thermal stability from -200 to +260 °C.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

The advantages of PFA:

- Especially suitable for use in trace analysis
- Ideal for sensitive and valuable samples
- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no contamination due to the extremely hydrophobic, anti-adhesive and smooth surfaces
- High thermal stability, from -200 °C to +260 °C
- Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability
- Easy to clean
- Use of high purity raw materials

Further information on PFA can be found starting on page 122.

Volume ml	Divisions ml	Height mm	Ø mm	PU	Cat. No.
25	5	50	32	1	110205
50	10	59	39	1	110305
100	20	72	50	1	110405
250	50	96	67	1	110605
500	100	122	88	1	110905
1000	100	141	109	1	111005

Sample preparation

Griffin beakers, ETFE



Transparent. With easily readable, printed black graduations. According to ISO 7056. Very good chemical resistance and very high thermal stability from -100 to +150 $^{\circ}$ C. To preserve markings, cleaning at no higher than 60 $^{\circ}$ C is recommended.

Volume ml	Divisions ml	Height mm	Ø mm	PU	Cat. No.
25	5	50	32	1	110204
50	10	59	39	1	110304
100	20	72	50	1	110404
250	50	96	67	1	110604
400	50	109	77	1	110704
500	100	122	88	1	110904
600	100	125	91	1	110804
1000	100	143	105	1	111004



Griffin beakers, PTFE



White, opaque, thick-walled. Without scale.

Excellent chemical resistance and high thermal stability from -200 to +260 °C.

Also suitable for microwave ovens.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Volume ml	Wall thickness mm	Height mm	Ø mm	PU	Cat. No.
5	2	26	20	1	112197
10	2	33	24	1	112297
25	2	47	32	1	112397
50	2	60	43	1	112497
100	3	68	54	1	112597
250	3	97	66	1	112697
500	4	125	80	1	112797
1000	4	155	100	1	112897



Watch glasses can be found on page 78.



Griffin beakers, PMP, printed red graduations

Crystal clear. With easily readable, printed red graduations.

According to ISO 7056.

To preserve markings, cleaning at no higher than 60 °C is recommended.

Volume	Divisions	Height	Ø	PU	Cat. No.		
ml	ml	mm	mm				
10*	2	36	30	12	60503		
25	5	50	38	12	60603		
50	10	60	47	12	60703		
100	10	70	55	12	60803		
150*	20	80	66	12	60903		
250	25	95	77	6	61003		
400*	50	112	87	6	61103		
500	50	118	94	6	61803		
600*	50	127	100	6	61203		
1000	100	147	120	6	61403		
2000	200	187	149	6	61503		
3000	250	212	170	4	61603		
5000	500	247	203	4	61703		
* Variant in addition to ISO 7056							



Griffin beakers, PMP, molded graduations



Crystal clear. With molded graduations.

According to ISO 7056.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Volume ml	Divisions ml	Height mm	Ø mm	PU	Cat. No.			
25	5	50	38	12	60695			
50	10	60	47	12	60795			
100	10	70	55	12	60895			
150*	20	80	66	12	60995			
250	25	95	77	6	61095			
400*	50	112	87	6	61195			
500	50	118	94	6	61895			
600*	50	127	100	6	61295			
1000	100	147	120	6	61495			
2000	200	187	149	6	61595			
3000	250	212	170	4	61695			
5000	500	247	203	4	61795			
* Variant in addition to ISO 7056								

* Variant in addition to ISO 7056

Sample preparation

Griffin beakers, PP, molded blue graduations



Highly transparent. With easily readable molded, embossed blue graduations. According to ISO 7056.

To preserve markings, do not clean at temperatures exceeding 60 °C.

Conditionally autoclavable at 121 °C (2 bar) according to DIN EN 285. For autoclaving we recommend the design with molded graduations (Cat. No. 606941 – 617941).

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Divisions	Height	Ø	PU	Cat. No.		
ml	ml	mm	mm				
10*	2	36	30	12	605081**		
25	5	50	38	12	606081**		
50	10	60	47	12	607081**		
100	10	70	55	12	608081		
150*	20	80	66	12	609081		
250	25	95	77	6	610081		
400*	50	112	87	6	611081		
500	50	118	94	6	618081		
600*	50	127	100	6	612081		
1000	100	147	120	6	614081		
2000	200	187	149	6	615081		
3000	250	212	170	4	616081		
5000	500	247	203	4	617081		
* Variant in addition to ISO 7056							



Griffin beakers, PP, molded graduations



Highly transparent. With molded graduations.

** Blue printed graduations, not molded

According to ISO 7056.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

50 10 60 47 12 60794 100 10 70 55 12 60894 150* 20 80 66 12 60994 250 25 95 77 6 61094 400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294	Volume	e Divisions	Height	Ø	PU	Cat. No.
50 10 60 47 12 60794 100 10 70 55 12 60894 150* 20 80 66 12 60994 250 25 95 77 6 61094 400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294	ml	ml	mm	mm		
50 10 60 47 12 60794 100 10 70 55 12 60894 150* 20 80 66 12 60994 250 25 95 77 6 61094 400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294						
100 10 70 55 12 60894 150* 20 80 66 12 60994 250 25 95 77 6 61094 400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294	25	5	50	38	12	606941
150* 20 80 66 12 60994 250 25 95 77 6 61094 400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294	50	10	60	47	12	607941
250 25 95 77 6 61094 400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294	100	10	70	55	12	608941
400* 50 112 87 6 61194 500 50 118 94 6 61894 600* 50 127 100 6 61294	150*	20	80	66	12	609941
500 50 118 94 6 61894 600* 50 127 100 6 61294	250	25	95	77	6	610941
600* 50 127 100 6 61294	400*	50	112	87	6	611941
	500	50	118	94	6	618941
1000 100 147 120 6 61494	600*	50	127	100	6	612941
	1000	100	147	120	6	614941
2000 200 187 149 6 61594	2000	200	187	149	6	615941
3000 250 212 170 4 61694	3000	250	212	170	4	616941
5000 500 247 203 4 61794	5000	500	247	203	4	617941

^{*} Variant in addition to ISO 7056



Watch glasses can be found on page 78.




Measuring cup, PP

Transparent. Very durable. With molded graduations. Can be used as dispensing or lab beaker. Not a medical device.

Fitting lid made of PE - please order separately.

Description	Volume	Divisions	Height	Ø	PU	Cat. No.
	ml	ml	mm	mm		
Massuring sup	20	1	42	37	100	60204
Measuring cup	30	I	42	3/	100	69394
Snap-on lid, white, PE					100	69493

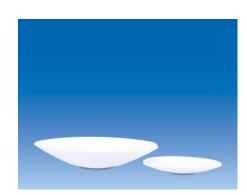

Watch glasses, PTFE



White. Without base.

High thermal stability and chemical resistance.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Ideal for use to cover beakers.

Ø	PU	Cat. No.
mm		
50	1	113197
_75	1	113297
100	1	113397
125	1	113497



Watch glasses, PP



Transparent. With base.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Ideal for use to cover beakers.

Ø	PU	Cat. No.
mm		
60	10	80452
80	10	80454
100	10	80455
118.5	10	80456

Sample preparation

Erlenmeyer flasks, GL 45, PMP with PP screw cap



Transparent.

Ideal for use as a receiving vessel in titrations.

Well suited for storage and cultivation of cell cultures. Far safer than glass flasks for use in incubator shakers due to the break resistance of plastic. Suitable for microwaves.

To preserve markings, cleaning at no higher than 60 °C is recommended.

Volume ml	Divisions ml	PU	Cat. No.
75	10	6	56695
125	20	6	56795
250	50	6	56895
500	100	6	56995
1000	200	4	57095



Erlenmeyer flasks, GL 45, PP with PP screw cap



Transparent.

Well suited for storage and cultivation of cell cultures. Far safer than glass flasks for use in incubator shakers due to the break resistance of plastic. Suitable for microwaves.

To preserve markings, cleaning at no higher than 60 °C is recommended.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Divisions ml	PU	Cat. No.
75	4.0		566044
75	10	6	566941
125	20	6	567941
250	50	6	568941
500	100	6	569941
1000	200	4	570941




Screw caps, GL 45, PP

Screw caps suitable for Erlenmeyer flasks and VITgrip™ lab bottles.

Colour	PU	Cat. No.
red	6	83340
blue	6	83350



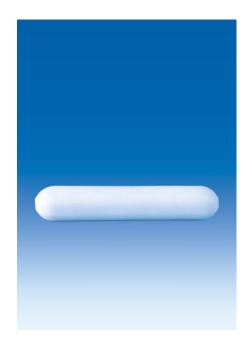

Magnetic stirring bars, polygonal, PTFE



With permanent magnet AlNiCo V core. The angled shape gives rise to significant turbulence, and thus achieves effective mixing, even at low spin rates.

Ø	Length	PU	Cat. No.
mm	mm		
2	-	40	200407
2	5	10	300497
2	7	10	300597
3	8	10	300897
3	10	10	301097
3	13	10	301197
4.5	12	10	301597
6	10	10	301697
6	15	10	301797
6	25	10	301997
6	30	10	302097
7	20	10	301897
7	50	10	302297
7	60	10	302397
8	40	10	302197
10	70	10	302497
10	80	10	302597
27	57	10	303097
27	108	10	303197
27	159	10	303297



Magnetic stirring bars, oval, PTFE



With permanent magnet AlNiCo V core. Highly suitable for vessels with a round bottom, such as round-bottom flasks. The angled side surfaces act to produce high turbulence, and thus achieve effective mixing.

Ø	Length	PU	Cat. No.
mm	mm		
5	10	10	311097
6	15	10	311197
10	20	10	311297
12	25	10	311397
16	30	10	311497
16	35	10	311597
20	40	10	311697
20	50	10	311797

Sample preparation

Magnetic stirring bars, octagonal, PTFE



With rings and permanent magnet AlNiCo V core. The eight-sided shape gives rise to significant turbulence, and thus achieves effective mixing, even at low spin rates. The middle ring also promotes stable centering with convex or uneven bottoms.

Ø mm	Length mm	PU	Cat. No.
8	13	10	307697
8	15	10	307797
8	22	10	307897
8	25	10	307997
8	28	10	308097
8	38	10	308197
8	41	10	308297
8	51	10	308397
8	64	10	308497
10	13	10	308597
10	25	10	308897
10	35	10	308997
10	38	10	309097
10	51	10	309297
10	64	10	309397



Magnetic stirring bars, cross shape, PTFE



With permanent magnet AlNiCo V core. The angled shape gives rise to significant turbulence, and thus achieves effective mixing, even at low spin rates. The cross shape promotes very stable centering.

Size	Height	PU	Cat. No.
mm	mm		
10 x 10	5	10	316097
20 x 20	8	10	316197
25 x 25	9	10	316297
30 x 30	10	10	316397
38 x 38	11	10	316497





Magnetic stirring bars, triangular, PTFE



With permanent magnet AlNiCo V core. The angled shape gives rise to significant turbulence, and thus achieves effective mixing, even at low spin rates.

Edges	Length	PU	Cat. No.
mm	mm		
6	12	10	310197
8	25	10	310297
14	40	10	310397
12	50	10	310497

Magnetic stirring bars, barbell, PTFE



With permanent magnet AlNiCo V core. Excellent centering due to the small contact surface, and effective mixing. Disk diameter: 20 mm; shaft diameter: 8 mm.

Length	PU	Cat. No.
35	10	3125970
55	10	3126970




Magnetic stirring bar retrievers, flexible, PTFE





Flexible magnetic stir bar retriever with a total length of 330 mm. Magnet encapsulated. \emptyset x L: 12.5 x 51 mm. Due to the high flexibility, the magnet stir bar can be retrieved from inaccessible locations, e.g., from the water trap in a laboratory sink. High chemical resistance, simple to clean.

Length	PU	Cat. No.
mm		
330	1	318597

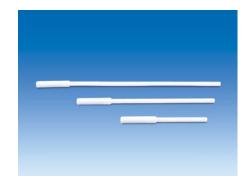

Magnetic stirring bar retrievers, PTFE





With PTFE encapsulated magnetic core. Straight shape. High chemical resistance, simple to clean.

Length	PU	Cat. No.
mm		
150	1	122097
250	1	122197
350	1	122297

Sample preparation

Magnetic stirring bar retrievers, PE



With a permanent magnet on one end and holding ring on the other one. Magnet is completely encapsulated in a PE mantle.



Mortars, MF

White, with spout. Stable circumferential edge. Very stable.


Volume ml	Height mm	Ø mm	PU	Cat. No.
300	75	125	5	72898
500	90	150	5	72998



Pestles, MF

White, heavy design. with ergonomically shaped grip.

Length mm	Head Ø mm	PU	Cat. No.
125	30	5	73498
145	35	5	73598
160	40	5	73698
215	42	1	73898

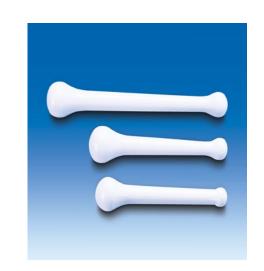



Urbanti funnels, PMP

Crystal clear. The spiral-shaped ribs increase the rate of filtration and prevent the trapping of air between the filter paper and the funnel. With long stem.

Volume	Ø	Length	Stem Ø	Stem length	PU	Cat. No.
approx. ml	mm	mm	mm	mm		
30	51	195	3	150	6	325095
80	70	210	3	150	6	325195
250	100	198	7	108	4	325295
630	140	247	10	132	3	325395
1800	196	315	20	155	2	325495



Analytical funnels, PP

Transparent. With long stem and grooves. Rigidified by a thickened edge. Rapid flow due to a steep 60° angle.

Volume	Ø	Length	Stem Ø	Stem length	PU	Cat. No.
approx. ml	mm	mm	mm	mm		
50	50	194	5	150	10	80162
100	72	208	5	143	10	80164
225	91	227	5	145	10	80165



Büchner funnels, PP

Two parts. Upper and lower parts are detachable to facilitate cleaning.

Volume approx. ml	Filter Ø mm	Length mm	Hole Ø mm	PU	Cat. No.
40	42.5	95	1.2	1	80437
70	55	113	1.1	1	80438
180	70	145	2.0	1	80439
280	80	165	2.0	1	80440
390	90	180	2.5	1	80441
810	110	210	2.5	1	80442
2100	160	280	2.75	1	80443
6000	240	350	3.0	1	80445

Water-jet pump, PP

For generation of a vacuum and to siphon off liquids and vapours (if necessary, with a suction main or condensation trap connected upstream).

Pump fluid: Water

Length of the unit: approx. 210 mm (R 3/4" connector fitted)

Weight: approx. 33 g (R 3/4" connector fitted)

- High chemical resistance, since the pumped media only come into contact with polypropylene, FKM and PTFE.
- Operating temperature up to a maximum of 80 C°.
- Integrated non-return valve increases operating safety.
- Simple operation, and easy to clean.
- Detachable vacuum connection.
- Variety of supplied adapters simplify connections to most water sources, and additional reducing adapters are available.

Very low water consumption:

The flow configuration has been optimised, resulting in a 33% reduction of water consumption (approx. 220 liters/h at 3.5 bar water supply pressure).

Constant discharge pressure:

The discharge pressure of 16 mbar (water temperature: 12 °C) is reachable across a wide range of water supply pressures of from 3 to 6 bar.

High suction capacity:


The suction capacity is approx. 450 L/hour of air (\pm 50 l/h, vs. atmospheric pressure, 12 °C water temperature, 3.5 bar water supply pressure).


Included in delivery:

Water jet pump, including: water supply connection (lock nut R $\frac{3}{4}$ ", reducing adapter R $\frac{1}{2}$ " and hose connection (olive) with outer Ø 10-12 mm), vacuum connection (olive with outer Ø 6-9 mm, detachable with GL 14 screw cap).


Description	PU	Cat. No.
Water-jet pump Accessories:	1	77094
Reducing adapter R3/4 inch to R3/8 inch	1	159665
Reducing adapter R3/4 inch to M 22x1, Perlator thread	1	159670



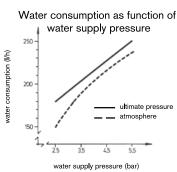

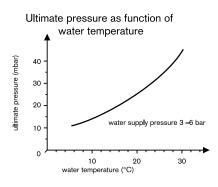

Pumping capacity as function of back pressure

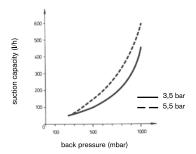



Filtering racks

Funnel holder with base and adjustable height, made from PP, support stand made from aluminium, diameter: 12.7 mm; and, length: 595 mm. To hold from two to four funnels with an upper outer diameter of 50-120 mm.

Positions	Base plate	PU	Cat. No.
	mm		
2	250 x 140	1	78394
4	450 x 140	1	78294

Imhoff Sedimentation Cone, SAN

According to DIN 12 672. Crystal clear, with molded graduations for precise reading of volumes. For simple, basic cleaning and rinsing, the screw coupling on the tip can be removed. Lower breakage risk than for PC or glass containers.

For determination of suspended matter in liquids (e.g., for industrial and municipal wastewater).

Graduation: Divisions:		Error limits:	
0 -	2 ml	0.1 ml	+/- 0.1 ml
2 -	10 ml	0.5 ml	+/- 0.5 ml
10 -	40 ml	1 ml	+/- 1 ml
40 -	100 ml	2 ml	+/- 2 ml
100 - 1	1000 ml	50 ml	+/- 10 ml

Volume ml	PU	Cat. No.
1000	3	75991



Sedimentation rack, PMMA

For 2 Imhoff Sedimentation Cones. Base plate with depression for exact vertical positioning of the sedimentation cone.

L x W x H mm	PU	Cat. No.
150 x 300 x 290	1	81056

Evaporating dishes, PFA



With snap-on lid, PE. For contamination free sample preparation and efficient transportation. Due to a conical depression in the middle of the base very small amounts of a solvent are adequate to absorb the evaporated samples.

Volume ml	Height mm	Ø mm	PU	Cat. No.
25	25	50	1	103297
50	54	50	1	103397







Round-bottom flasks, PFA





Transparent, neck with NS 29/32. Suitable as safety flask for use with rotary evaporators (operation at room temperature) to collect the distilled off liquid. High thermal stability and chemical resistance.

The advantages of PFA:

- Especially suitable for use in trace analysis
- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no cross-contamination due to the extremely hydrophobic, anti-adhesive and smooth surfaces
- High thermal stability, from -200 °C to +260 °C, autoclavable.
- Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability
- Easy to clean
- Use of high purity raw materials

Further information on PFA can be found starting on page 122.

Volume	Height	Ø	PU	Cat. No.
ml	mm	mm		
100	117	65	1	107797
250	147	88	1	107897
500	177	107	1	107997



Round-bottom flask stand, PP



White, for flasks with a round bottom. Excellent chemical resistance. Autoclavable at $121~^{\circ}$ C (2 bar) according to DIN EN 285.

Ø	PU	Cat. No.
mm		
160	5	80271

Gas wash bottles, PFA



Cap with S 40 buttress threads, and screwable frit made of PTFE. A pore size of approx. 3 µm for optimal optimal pearling of the gas into the liquid. A wide field of application is possible due to the use of high-quality fluoroplastic. Suitable only for non-pressurised operation.

The advantages of PFA:

- Especially suitable for use in trace analysis
- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no cross-contamination due to the extremely hydrophobic, anti-adhesive and smooth surfaces
- High thermal stability, from -200 °C to +260 °C, autoclavable
- · Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability
- Easy to clean
- Use of high purity raw materials

Further information on PFA can be found starting on page 122.

Height	Ø	Hose connection	PU	Cat. No.
mm	mm	inner / outer Ø mm		
160	61	4/6	1	159497
190	76	4/6	1	159597
240	96	5/8	1	159697
	mm 160 190	mm mm 160 61 190 76	mm mm inner / outer Ø mm 160 61 4/6 190 76 4/6	mm mm inner / outer Ø mm 160 61 4 / 6 1 190 76 4 / 6 1



Desiccators with stopcock, PC

Crystal clear, with stopcock for evacuation. Lower sections can be filled with desiccants. The materials to be dried are placed on a perforated disc made of PP. Lid is sealed by a neoprene gasket. Ideal for use in educational laboratories.

Ø	Disc Ø	Height	PU	Cat. No.
mm	mm	mm		
171	140	206	1	326496
230	190	260	1	326596
273	230	311	1	326696





Desiccators, PP/PC

Lower portion made from PP can be filled with desiccants. The materials to be dried are placed on a perforated disc made of PP. The lid made from PC is sealed with a neoprene gasket. Ideal for use in educational laboratories.

Ø	Disc Ø	Height	PU	Cat. No.
mm	mm	mm		
171	140	206	1	326094
230	190	260	1	326194
273	230	311	1	326294



Desiccators with stopcock, PP/PC

With a bleed valve and an O-ring seal between the dome and the lower portion. The desiccators are suitable for vacuum and are provided with a non-return valve. Hot crucibles should only be placed on a porcelain plate, and should not come into too close contact with the rim of the desiccator. The insert made of PP serves to hold the drying agent. Desiccator plates should be purchased separately.

Disc Ø	Height	PU	Cat. No.
mm	mm		
140.5	190	1	80550
189	230	1	80230
238	300	1	80554
	mm 140.5 189	mm mm 140.5 190 189 230	mm mm 140.5 190 1 189 230 1

Desiccator plates, PP and porcelain

The PP plates are ideally suited for storage of crucibles and other containers at room temperature. The porcelain plates can also be used for hot crucibles.


Material	For desiccator Ø	Ø	PU	Cat. No.
	mm	mm		
PP	150	140.5	1	80551
PP	200	189	1	80231
PP	250	238	1	80553
Porcelain	150	140	1	65965
Porcelain	200	190	1	65975
Porcelain	250	235	1	65980



Replacement parts for desiccators

Replacement O-rings and valves for desiccators (cat. nos. 80550, 80230, 80554).

Description	PU	Cat. No.
O-ring for desiccator no. 80550	1	80555
O-ring for desiccator no. 80230	1	80556
O-ring for desiccator no. 80554	1	80557
Valve, PC, for desiccator nos. 80550, 80230 and 80554	1	80229

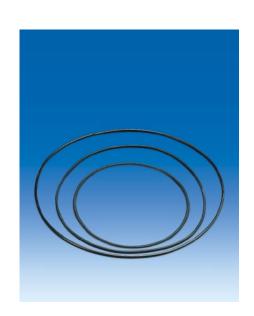



Sample containers, PFA





With screw cap made of PFA. Cylindrical, tall shape. Ideal for sample collection, transport and storage of samples.

The advantages of PFA:

- Especially suitable for use in trace analysis
- Ideal for sensitive and valuable samples
- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no cross-contamination due to the extremely hydrophobic, anti-adhesive, smooth surfaces
- High thermal stability from -200 °C to +260 °C, autoclavable at 121 °C (2 bar) according to DIN EN 285
- Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability
- Easy to clean
- Use of high purity raw material

Further information on PFA can be found starting on page 122.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	GL	mm	mm		
30	40	54	38	1	130297
60	40	90	38	1	130397
90	56	62	54	1	130497
180	56	112	54	1	130597

Sample containers, PP



Transparent. With screw cap made of PP. Cylindrical, tall shape. Ideal for sample collection, transport and storage of samples. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	GL	mm	mm		
30	40	54	38	10	130294
60	40	90	38	10	130394
90	56	62	54	10	130494
180	56	112	54	10	130594



Sample containers, PP



Transparent. With snap-on lid made of PE-LD. Conical shape. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Height mm	Ø mm	PU	Cat. No.
5	25	20	25	68594
18	57	22	25	68894
50	97	30	10	69194
160	110	50	10	69294



Sample containers, PE-LD

Transparent. With attached plug-seal cap made of PE-LD.

Volume ml	Height mm	Ø mm	PU	Cat. No.
1	32	8	500	80730
2.5	32	14	100	80731
5	50	15	100	80737
8	57	17	100	80732
7	32	23	100	80733
20	74	25	100	80734
25	52	31	50	80736
35	74	31	50	80735
-				





Weighing jars, PP

Transparent. With knobbed lid. Cylindrical shape.

Volume	Height	Ø	PU	Cat. No.
ml	mm	mm		
20	30	40	10	80342
23	50	30	10	80340
30	30	50	10	80345
50	35	60	10	80346
60	69	40	10	80343
190	90	60	10	80347
360	120	70	10	80348



Sample tubes, PFA



Sample tubes made from PFA for sample preparation and for use in autosampler racks. With or without individually calibrated ring mark at 10 ml with GL 25 screw cap made from PFA or PE stopper (see Table).

The advantages of PFA:

- Especially suitable for use in trace analysis
- No memory effects
- Practically no carryover due to the extremely hydrophobic, anti-adhesive, smooth surfaces
- High thermal stability, from -200 °C to +260 °C, autoclavable
- Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability
- Easy to clean

Further information on PFA can be found starting on page 122.

Figure Type No.		Volume ml	Height mm	Ø mm	PU	Cat. No.
1	With ring mark and screw cap	15	110	22	1	103897
-	Without ring mark	15	110	22	1	1038971
2	With ring mark and stopper	12	110	16	1	1037979
3	Without ring mark	12	110	16	1	103797

Sample vials, PFA



Sample vials made of PFA with conical interior and molded graduation (5 ml subdivisions). Available in two different types, depending on application:

- Cored outside bottom
- Flat surface on bottom of vial for improved heat transfer (Recommended for use with hot plates)

Both 50 ml sizes fit in common Autosampler racks.

Scope of delivery is without screw cap. Please order the screw cap (Cat. No. 104997) separately.

Volume ml	Type outside bottom	Ø mm	Height* mm	PU	Cat. No.
15	Flat	29	39	1	104197
15	Cored	29	42	1	104097
25	Flat	29	69	1	104397
25	Cored	29	72	1	104297
50	Flat	29	117	1	104597
50	Cored	29	120	1	104497
Screw cap, 33 mm, PFA (suitable for sample vials (104097 – 104597)					104997
* Height w	vith thread				



Autosampler-vials, PFA



Molded graduation with 1 ml subdivisions.

Translucent material for optimum visibility of the liquid contained in the vial.

Conical interior for use with autosamplers. The outer design allows easy handling of the autosampler-vials with forceps. Optionally available with plug-seal cap for long term storage or lid with knob for fast opening and closing (dust protection) of the vial.

Scope of delivery is without lid or cap. Please order the fitting lid (Cat. No. 105597 resp. 105697) separately.

Volume	Ø	Height	PU	Cat. No.
ml	mm	mm		
1.5	13.5	24	1	105097
2.5	13.5	36	1	105197
4	14	52	1	105297
Plug-seal cap, PFA	18	5	1	105697
Lid with knob, PFA	16	9	1	105597





Reagent tube racks, PP

White. For reagent tubes with a diameter of 20 mm.

With three levels for precise, vertical positioning of the reagent tubes.

Working temperatures of -20 to +90 $^{\circ}$ C.

For \emptyset up to	Positions	LxWxH	PU	Cat. No.
mm		mm		
20	2 x 6	190 x 60 x 80	5	80560
20	2 x 12	375 x 65 x 95	5	80562



Reagent tube racks, PP

White. The special shape makes it possible to check the amounts present in the reagent tubes.

With two side-mounted handle straps.

For Ø up to	Positions	LxWxH	PU	Cat. No.
mm		mm		
16	10	220 x 55 x 65	4	80130
18	9	220 x 55 x 65	4	80131
Base plate for 2 reager	nt tube racks	202 x 156 x 13.5	4	80134
(Cat. No. 80130, 8013	1)			

Test tube racks, PP, coloured

Stackable, simple, and small footprint. Alphanumerically identified positions. Suitable for tempering in a water bath as well as storage of samples in the refrigerator and incubation in a climate chamber. The racks are supplied as folded out flat, and can be firmly and inseparably joined together in just a few steps. Working temperatures of -20 to +90 $^{\circ}$ C. Base area: 265 x 126 mm.


For Ø up to mm	Positions	Height mm	Colour	PU	Cat. No.
13	6 x 14	75	White	5	3190940
16	5 x 11	75	White	5	3191940
18	5 x 11	75	White	5	3192940
20	4 x 10	75	White	5	3193940
25	4 x 8	88	White	5	3194940
30	3 x 7	88	White	5	3195940
13	6 x 14	75	Blue	5	3190948
16	5 x 11	75	Blue	5	3191948
18	5 x 11	75	Blue	5	3192948
20	4 x 10	75	Blue	5	3193948
25	4 x 8	88	Blue	5	3194948
30	3 x 7	88	Blue	5	3195948
13	6 x 14	75	Red	5	3190943
16	5 x 11	75	Red	5	3191943
18	5 x 11	75	Red	5	3192943
20	4 x 10	75	Red	5	3193943
25	4 x 8	88	Red	5	3194943
30	3 x 7	88	Red	5	3195943



Microtube racks, PP, coloured

Stackable racks for micro- or cryotubes. Alphanumerically identified positions. Suitable for tempering in a water bath. The racks are supplied as folded out flat, and can be firmly and inseparably joined together in just a few steps. Working temperatures of -20 to +90 $^{\circ}$ C. Base area: 265 x 126 mm.

For Ø up to mm	Positions	Colour	Height mm	PU	Cat. No.
11	8 x 16	White	38	5	3197940
13	6 x 14	White	38	5	3198940
11	8 x 16	Blue	38	5	3197948
13	6 x 14	Blue	38	5	3198948
11	8 x 16	Red	38	5	3197943
13	6 x 14	Red	38	5	3198943

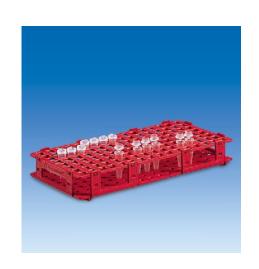



Staining chamber for slides, POM

Consists of staining jar and a staining rack for 25 slides 76 x 26 mm.

LxWxH	PU	Cat. No.
mm		
100 x 87 x 51	5	99099



Staining rack for slides, POM

For serial staining of 25 slides 76 x 26 mm. Usable in the staining jar, cat. no. 99199 (see below).

LxWxH	PU	Cat. No.
mm		
91 x 79 x 38	10	99299



Staining jar for slides, POM

For combination with staining rack, cat. no. 99299 (see above). Can also be used as a storage box for 25 slides.

LxWxH	PU	Cat. No.
mm		
100 x 87 x 51	5	99199



Coplin staining chamber, PP

With screw cap made of PP, GL 45. For serial staining of 10 slides 76 x 26 mm.

Height mm	PU	Cat. No.
94	10	136693

Competence in Plastic Labware SAVING AND STORAGE WITH CONFIDENCE



VITgrip™ – The Allround Lab Bottle

The VITgrip™ lab bottles and screw caps of high-quality polypropylene (PP) are "Made in Germany". VITgrip™ lab bottles are in many aspects a safe alternative to glass, as for example, the higher break resistance reduces the risk of injury and together with outstanding chemical resistance provides a long period of use.

Leakproof* and break resistant



Sometimes it happens very quickly: one inattentive moment and, by accident, the lab bottle is knocked over. Breakage of glass can be dangerous because of possible injury due to glass splinters and/or spilled liquid. VITgripTM lab bottles made of plastic provide a higher level of safety in the lab because the VITgripTM has a significantly higher break resistance and is leakproof*. The bottle thread and the associated screw cap are an ideally matched pair. Together, they form a reliable sealing system without the need of an additional seal that can wear, corrode or cause contamination. Both components are subject to a detailed quality inspection prior to delivery.

Safe storage



The VITgrip™ lab bottles are supplied with a tamper-evident closure; i. e. a ring, which is attached at the lower end of the screw cap will tear off upon the first opening of the closed bottle. It reliably signals, if the bottle is still sealed before opening. Thus, an intact tamper-evident closure can ensure safe storage of e. g. reference samples or safe transfer of samples between sampling site and lab. After the ring is torn off, the closure can be used as a regular screw cap. All VITgrip™ lab bottles have a GL 45 thread and an evenly formed neck area that allow controlled, smooth pouring of liquid.

VITgrip™ lab bottles, PP, GL 45 with tamper-evident cap, PP



Everyday use bottles made of plastic for sampling and storing liquids in the lab.

Due to the innovative design and the ergonomic shape the bottles have an outstanding easy-grip feature. The slim and tapered shape greatly improves handling in comparison to conventional lab bottles. Furthermore, the molded volume graduations provide texture for an enhanced grip, especially when working with gloves.

Through the optimised sealing system of bottle threads and screw cap the bottle is leak-proof* and offers an optimum pouring behaviour and easy cleaning due to the hydrophobic material and round shape.


With double-sided molded graduations (accuracy \pm 5%), the volume inside the bottle is easy to read, even during use.

The included tamper-evident cap reliably signals if the bottle is still sealed before opening.

The bottle has very good chemical resistance against most acids, bases and alcoholic solutions.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Graduation	Height**	Bottom-ø	PU	CatNo.
ml	ml	mm	mm		
125	12.5	103	54	6	110194
250	25	149	64	6	110294
500	25	192	77	6	110394
1000	50	234	97	6	110494
2000	100	278	126	1	110594
Replaceme	ent tamper-evident	cap, PP, GL 45		6	83330
Starter-Set (3 x VITgrip™ (250 / 500 / 1000 ml) 111					
+ 3 x tamper-evident cap					

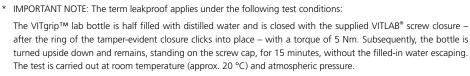

^{**} Height without screw cap.







Narrow-mouth bottles, PFA





Transparent.

With screw cap with buttress threads made of PFA. Ideal for long-term storage of high-purity oxidants, acids, alkalis, as well as hydrocarbons, trace analysis solvents and standards.

The advantages of PFA:

- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no cross-contamination due to the extremely hydrophobic, anti-adhesive, smooth surfaces
- High thermal stability, from -200 °C to +260 °C, autoclavable
- Chemical inertness against nearly all chemicals
- · Good transparency and dimensional stability
- Easy to clean
- Use of high purity raw materials

Further information on PFA can be found starting on page 122.

Volume	Thread	Height	Ø	PU	Cat. No.
ml		mm	mm		
50	S 28	86	37	1	109297
100	S 28	120	45	1	109397
250	S 28	160	61	1	108297
500	S 28	190	76	1	108397
1000	S 28	240	96	1	108497



Wide-mouth bottles, PFA





Transparent.

With screw cap made of PFA with buttress threads. Ideal for long-term storage of high-purity oxidants, acids, alkalis, as well as hydrocarbons, trace analysis solvents and standards.

Volume ml	Thread	Height mm	Ø mm	PU	Cat. No.
250	S 40	150	61	1	109497
500	S 40	179	76	1	109597
1000	S 40	217	96	1	109697
2000	S 40	245	130	1	109797

Saving and storing

Screw caps, PFA



Transparent. For sealing all PFA containers with GL threads or buttress threads. Autoclavable at 121 $^{\circ}$ C (2 bar) according to DIN EN 285.

The advantages of PFA:

- Especially suitable for use in trace analysis
- Ideal for sensitive and valuable samples
- Long-term maintenance of low-concentration reference materials in PFA containers
- No memory effects
- Practically no carryover, no cross-contamination due to the extremely hydrophobic, anti-adhesive and smooth surfaces
- High thermal stability, from -200 °C to +260 °C, autoclavable
- Chemical inertness against nearly all chemicals
- Good transparency and dimensional stability
- Easy to clean
- Highly pure starting materials used

Further information on PFA can be found starting on page 122.

Thread	PU	Cat. No.
S 28	1	102697
S 40	1	102897





Wide-mouth bottles, PTFE





White. Opaque. Thick walled.

With screw cap made of PTFE.

Very high thermal stability and chemical resistance.

With very wide mouth, ideal for filling with powders and paste-like materials.

Volume ml	Thread mm	Height mm	Ø mm	PU	Cat. No.
10	12	50	26	1	122597
25	19	61	33	1	122697
50	25	76	43	1	122797
100	35	88	52	1	122897

Narrow-mouth bottles, PP





Transparent. With high shoulders.

With screw cap made of PP.

Good chemical resistance, ideal for long-term storage of liquids.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Thread GL	Height mm	Ø mm	PU	Cat. No.
250	25	135	65	12	94994
500	25	180	74	12	95094
1000	32	215	92	12	95194



Wide-mouth bottles, PP





Transparent.

With screw cap made of PP.

Good chemical resistance, ideal for long-term storage of liquids.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Wider mouth for simple filling, also ideal for powders and paste-like materials.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	GL	mm	mm		
250	45	146	58	12	93794
500	45	166	76	12	93994
1000	63	226	91	12	94194

Narrow-mouth bottles, PE-LD



Transparent. With high shoulders.

With screw cap made of PP.

Flexible material with good resilience.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Thread GL	Height mm	Ø mm	PU	Cat. No.
50	18	85	37	24	94589
100	18	114	43	24	94689
250	25	135	65	12	94989
500	25	180	74	12	95089
1000	32	221	92	12	95189



Wide-mouth bottles, PE-LD



Transparent.

With screw cap made of PP.

Flexible material with good resilience.

Wider mouth for simple filling, also ideal for powders and paste-like materials.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume ml	Thread GL	Height mm	Ø mm	PU	Cat. No.
50	32	87	39	24	93389
100	32	94	47	24	93489
250	45	146	58	12	93789
500	45	166	76	12	93989
1000	63	226	91	12	94189



Screw caps, PP



Transparent. Autoclavable at 121 $^{\circ}\text{C}$ (2 bar) according to DIN EN 285.

Thread GL	PU	Cat. No.
18	24	83310
25	12	83311
32	12	83312
40	12	83315
45	12	83313
52	12	83316
56	12	83317
63	12	83314





Narrow-mouth bottles, PE-LD

Transparent. With flat shoulders.

With screw cap made of PE-LD.

Flexible material with good resilience.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	GL	mm	mm		
10	14	50	26	100	138093
20	14	58	31	100	138193
30	14	66	34	100	138293
50	18	85	39	100	138393
100	18	106	45	50	138493
250	25	140	59	50	138593
500	25	180	75	50	138693
1000	28	212	94	25	138793
2000	28	264	117	25	138893



Wide-mouth bottles, PE-LD

Transparent.

With screw cap made of PE-LD.

Flexible material with good resilience.

Wider mouth for simple filling, also ideal for powders and paste-like materials.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	GL	mm	mm		
50	32	80	38	100	139393
100	32	94	48	50	139493
250	40	126	62	50	139593
500	50	155	76	50	139693
1000	65	208	93	25	139793
2000	65	246	120	25	139893

Saving and storing

Narrow-mouth bottles, PE-HD

Transparent.

With screw cap made of PE-LD.

Small footprint due to the square cross-section and the high shoulders.

Volume ml	Thread GL	Height mm	Size mm	PU	Cat. No.
100	25	76	43 x 43	24	91789
250	28	80	80 x 80	24	91989
500	32	106	90 x 90	12	92089
1000	32	187	80 x 80	12	92189



Wide-mouth bottles, PE-HD

Transparent.

With screw cap made of PE-LD.

Small footprint due to the square cross-section and the high shoulders.

Wider mouth for simple filling, also ideal for powders and paste-like materials.

Volume ml	Thread GL	Height mm	Size mm	PU	Cat. No.
100	32	78	46 x 46	24	92489
250	50	83	80 x 80	24	92689
500	65	120	90 x 90	12	92789
1000	65	168	90 x 90	12	92889



Wide-mouth bottles, PE-LD, with eye closure

Transparent.

With eyes on the bottle and the screw cap for sealing.

With sealing plug and screw cap made of PP.

Wider mouth for simple filling, also ideal for powders and paste-like materials.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	mm	mm	mm		
5 0	2.4	7.5	40	25	00.400
50	24	75	40	25	80408
100	24	90	50	25	80409
250	36	130	60	25	80410
500	36	160	75	10	80411
1000	50	200	95	10	80412
2000	50	250	115	10	80413





Reagent bottles, PP



Transparent.

With screw cap made of PP.

Good chemical resistance, ideal for storage of liquids.

Autoclavable at 121 °C (2 bar) in compliance with DIN EN 285, except for sizes 5000 and 10000 ml.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Thread	Height	Ø	PU	Cat. No.
ml	GL	mm	mm		
100	18	100	52	20	100389
100	32	96	55	20	101589
250	25	132	70	20	100489
250	45	132	73	20	101689
500	25	165	87	10	100589
500	45	172	87	10	101789
1000	32	202	108	10	100689
1000	45	197	105	10	102089
1000	63	204	108	10	101889
2000	32	245	131	6	100789
2000	45	241	131	6	102189
2000	63	243	131	6	101989
5000*	45	315	178	1	100889
10000**	63	394	222	1	100989
* with hand	dle, PE-HD				
** with PF fo	nam seal and tw	o handles PF-H	D.		

with PE foam seal and two handles, PE-HD

Saving and storing

Reagent bottles, PP



Transparent.

With NS stopper made of PP.

Stopper type A: With square-knob cap and red core.

Stopper type B: With octagonal-knob cap and red core.

Good chemical resistance, ideal for long-term storage of liquids.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

Volume	Neck	Height	Ø	Stopper	PU	Cat. No.
ml	NS	mm	mm			
100	14/23	106	52	А	20	100394
100	29/32	111	55	В	20	101594
250	19/26	138	70	А	20	100494
250	34/35	144	73	В	20	101694
500	24/29	172	87	А	10	100594
500	45/40	183	87	В	10	101794
1000	29/32	213	108	А	10	100694
1000	60/46	214	108	В	10	101894



Standard joint stoppers, PP



Stopper type A: With square-knob cap and red core.

Stopper type B: With octagonal-knob cap and red core.

Autoclavable at 121 °C (2 bar) according to DIN EN 285.

NS	Model	PU	Cat. No.
10/19	А	10	90694
12/21	A	10	90794
14/23	А	10	90894
19/26	А	10	90994
24/29	А	10	91094
29/32	А	10	91194
29/32	В	10	92194
34/35	В	10	91294
45/40	В	10	91394
60/46	В	10	91494





Storage bottles, PE-HD, without tap

Transparent.

With stable carrying handle and screw cap.

Available in wide- and narrow-mouth models.

Volume	Neck inner Ø	Height	Ø	PU	Cat. No.
I	mm	mm	mm		
5	90	318	163	1	81640
5	45	335	163	1	81644
10	120	390	205	1	81642
10	55	415	205	1	81646



Storage bottles, PE-HD, with tap

Transparent. Narrow-mouth model.

With stable carrying handle and screw cap. The 25 and 50 l sizes come equipped with two carrying handles.

Complete with exchangeable, easily operated tap made from PP with a ¾" pipe fitting.

Volume	Neck inner Ø	Height	Ø	PU	Cat. No.
1	mm	mm	mm		
5	45	335	163	1	81660
10	55	415	205	1	81662
25	79.5	565	280	1	81664
50	79.5	700	350	1	81666



Tap for storage bottles, PP

Replacement tap for storage bottles made from PP (Cat. No. 81660 to 81666). Complete with $\frac{3}{4}$ " pipe fitting and rubber ring.

Description	PU	Cat. No.
Tap for storage bottles	1	80375

Saving and storing

Chemical waste disposal system, PE/PP

For collection of liquid chemicals in the laboratory. The inlet hopper made from PE-HD contains a self-closing float, overfill protection, and a splash guard. Additionally, a screw cap (GL 63) with sealing ring is included.

Volume	Height	Ø	PU	Cat. No.
I	mm	mm		
10	560	222	1	151594



Bowl, PP, with lid



White. Rectangular shape.

Broad, stable, easy to grip edge.

Especially easy to clean due to the rounded corners and edges and the smooth surfaces. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume I	L x W x H mm	PU	Cat. No.
17	430 x 331 x 195	1	43610



Transport containers, PE-HD

Transparent.

Easy stackable.

With reinforcing ribs and integrated carrying handles.

Broad, stable edge.

Volume	LxWxH	PU	Cat. No.
I	mm		
20	420 x 310 x 205	1	80602
46	600 x 365 x 260	1	80603
72	700 x 420 x 310	1	80604





Multi-purpose container, SAN

Crystal clear, with fitted lid. Planar bottom inside, reinforced edge outside for stable placement on the lab bench.

Ideal for dust-proof storage of small components, instruments and utensils.

Volume ml	L x W x H mm	PU	Cat. No.
4000	340 x 230 x 94	1	36491



Laboratory trays / catchment trays, PP



White. All-purpose. Robust design. Very good chemical resistance Rounded corners and edges. Smooth surfaces, easy to clean. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Inner bottom	Edge dimensions	Height	PU	Cat. No.
dimensions mm	mm	mm		
130 x 180	180 x 230	42	1	165094
180 x 240	250 x 310	65	1	165194
240 x 300	310 x 370	75	1	165294
300 x 400	420 x 520	120	1	165394
400 x 500	534 x 634	140	1	165494
500 x 700	648 x 846	160	1	165594



Bowls, PP




White. Round. With broad, stable edge and circumferential standing ring on the bottom. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

Volume	Height	Ø	PU	Cat. No.
1	mm	mm		
0.9	70	160	5	42594
1.7	80	200	5	42694
2.9	100	240	5	42794
4.3	120	280	5	42894
6.6	130	320	3	42994
9.2	150	360	3	43094
13.4	180	400	3	43194

Competence in Plastic Labware



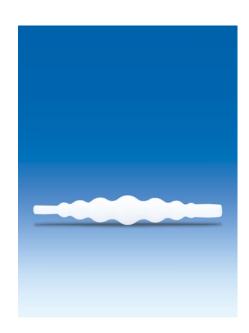

Connectors, universal, PP

Due to its special shape, it can be used for a wide variety of different tubing inner diameters.

For tubing with an inner \emptyset of	PU	Cat. No.
mm		
5 - 15	10	78794



Non-return valve, PE-HD

With valve disc made from FKM. Not suitable for pressure applications.

For tubing with an inner \emptyset of	PU	Cat. No.
mm		
6 - 9	10	78593



PTFE-tape

For sealing and wrapping of threads and other connections. High chemical resistance. Working temperatures up to 250 $^{\circ}$ C.

Length	Width	PU	Cat. No.
m	mm		
4.2	42	10	424007
12	12	10	131097

Lab assistants

Trays, MF



White. Flat shape. Rounded corners. Smooth surfaces, easy to clean. Practical tray for instruments, tools, and sensitive utensils. Stable and self-supporting. Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

LxWxH	PU	Cat. No.
mm		
190 x 150 x 17	5	71598
240 x 180 x 17	5	71698
268 x 208 x 17	5	71798
355 x 240 x 17	5	71898
428 x 288 x 17	5	71998



Trays, MF



White. Tall shape. Rounded corners. Smooth surfaces, easy to clean. Practical tray for instruments, tools, and sensitive utensils. Stable and self-supporting. Fitting lid made of PS, please order separately.

Suitable for contact with foodstuffs according to regulation (EU) No. 10/2011.

LxWxH	PU	Cat. No.
mm		
190 x 150 x 40	5	72098
290 x 160 x 35	5	72198
290 x 160 x 60	5	72398
340 x 245 x 100	5	72498
350 x 250 x 40	5	72298



Lids for trays, PS

Crystal clear. With handle. Rounded corners. Smooth surfaces, easy to clean. Protects the contents of the instrument trays from dust and contamination. The contents remain readily visible.

Size mm	For tray, MF No.	PU	Cat. No.
190 x 150	72098	5	79790
290 x 160	72198, 72398	5	79890
340 x 245	72498	5	79990*
* without handle			





Drying rack

Back plate and trough made from PVC with drainage nozzle.

With 75 metal pegs (length: 10 cm) with PE coating for hanging various sizes of apparatus. With two bore holes for simple wall mounting.

Delivered without installation hardware.

Size	PU	Cat. No.
mm		
450 x 630	1	76299



Drying rack, PS

With wide draining trough and drainage nozzles.

Rack with 72 pegs 95 x 15 mm. For drying larger apparatus, some of the pegs can be removed and the bore holes closed at the rear.

Delivered complete with drainage tube and accessories for the wall installation.

In addition, 11 pegs (95 x 6 mm) are included for objects having a smaller diameter, such as reagent tubes.

Description	Size	PU	Cat. No.
	mm		
Drying rack	450 x 630	1	80213
Pegs	95 x 6	11	81213

VITLAB® Promotional

ADVERTISE WITH YOUR GOOD NAME



Your good name in daily use

Precision is usually of great significance when it comes to ensuring the effective use of granulates, powders or liquids. The transportation, storage and decanting of small volumes often require special containers. VITLAB is one of the leading manufacturers of high-grade plastic labware and specializes in high-precision printing on plastic products with superior chemical and break resistance.

This offers a great advantage: by having your name and logo printed on these products, you will work "hand in hand" with your clients and always maintain a high presence. The products can be used wherever people work with granulates, powders or liquids; for example, in agriculture, laboratories, the medical sector, the food industry and the cleaning



business, as well as when using colours and chemicals.

Plastic labware by VITLAB guarantees that you have the best manufacturing quality and optimal functionality associated with your good name, thus ensuring a sustained positive echo.



A positive echo through individuality

VITLAB develops and manufactures its products at its own production facilities. This allows us to produce and print plastic labware according to your individual requirements and specifications. Please do not hesitate to inform us of your wishes and we will let you know what sort of individual solution we can provide for you.

When it comes to precision and accuracy

Volumetric containers by VITLAB stand out due to the greater precision and accuracy ensured by the measuring scale. On request, VITLAB can also print a customized scale on your product. The quality colours guarantee that the scale remains readable and does not wear off.


Small gifts keep the friendship alive

Plastic products have a high level of usability and are ideally suited as advertising articles or for promotional campaigns for your products. You can have your company name and logo or other motifs printed on them permanently and thus advertise with your good name.

A unique position thanks to an unmistakable design

VITLAB provides advertising materials of the very highest quality with round, conical or flat printing, using screen or pad printing systems and with particularly durable and luminescent colours according to the Pantone and HKS colour table. Various marking techniques, such as laser printing and heat embossing, provide you with an unmistakable design.



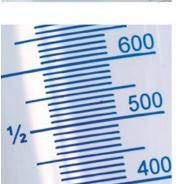





Would you like to have more information?

Please do not hesitate to contact us!

Please do not hesitate to contact us for advice on the selection, design and colour of your plastic products.

A personal consultant ensures that you receive competent advice from the first meeting to the delivery of the product.

Our contact information:


Tel: +49 6026 977 99-0


Fax: +49 6026 977 99 -30


E-mail: info@vitlab.com

www.vitlab-promotional.com

	Page
Fluoroplastic PFA	122
VITLAB® UV-protect	124
Plastics and their properties	125
 Classification and type description of plastics 	125
Chemical stability of plastics	129
Physical properties of plastics	138
Cleaning and maintenance of plastics	139
• Sterilisation of laboratory equipment made of plastics	140
Disposal and recycling	142
Suitability for foodstuffs	142
• CE-IVD	143
Accuracy of volumetric instruments	144

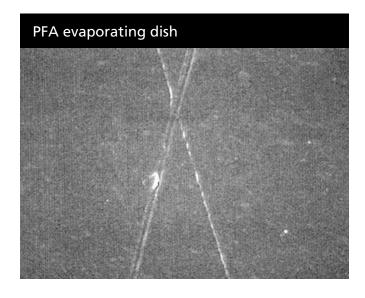

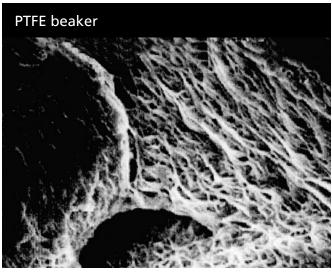

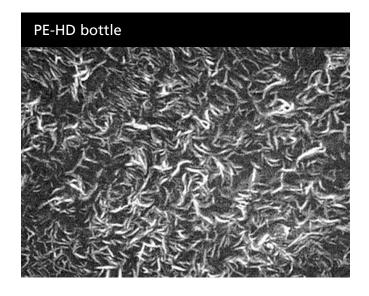

Pictures of the surface structure of PTFE, PFA, and PE-HD with a scanning electron microscope (8000 times magnification).

Fluoroplastic PFA

Today, trace analysis operates with concentrations in the range of ng/g (ppb) and pg/g (ppt). Thus, all modern processes require a corresponding trace analytical laboratory hygiene. However, the analytical accuracy of the measurement depends not only on the accuracy of the analytical instrument, but also directly on the judicious selection of container materials and the preparation of the samples themselves. Under these conditions, the fluoroplastic PFA truly shows its worth.

Labware made from polyolefins, such as polypropylene (PP) or polyethylene (PE), has found broad application in modern laboratories. Since catalysts (e.g., Ziegler-Natta or Philipps) are used during the manufacturing process, the constituent elements (frequently Al, Cr, Mg, Si, Ti or Zn) can still be detected in the ultra-trace range, and thus might affect the analysis results. In direct comparison to these manufacturing processes, PFA is especially suitable for trace element analysis because it is manufactured without the use of additives, and therefore cannot become contaminated by the additive components.

In addition to this advantage, PFA has other remarkable properties. PFA can be used for a broad range of applications by virtue of its resistance to almost all organic and inorganic compounds. PFA is one of the most resistant materials after platinum, and is practically chemically inert. Moreover, PFA stands apart with extraordinarily high thermal stability, making it possible to operate in the temperature range of -200 to +260 °C.

For all PFA products, VITLAB uses only high-purity PFA, which is particularly well-suited for trace analysis. For less critical uses, e.g. if the main requirement is a high chemical resistance, VITLAB offers bottles of "economy grade" PFA, which are partly made of recycled PFA. These are favourably priced and also environmentally friendly.

VITLAB's PFA vessels have unique, extremely smooth, liquid-repellent surfaces that are made possible through modern manufacturing processes and acknowledged expertise (see pictures on "Surface Structure"). This is especially significant in the illustrative comparison. The pictures taken in a scanning electron microscope show uneven and irregular surfaces for PE-HD and PTFE, and deep pores and concavities can be identified in the PTFE surface. In contrast to the uneven surfaces, PFA evaporating dishes had to be marked (X) and exhibit a completely smooth, even and uniform surface structure.

Due to this characteristic, all PFA labware is particularly easy to clean and presents hardly any interactions with samples as compared to conventional container materials. With these advantages combined with the minimal water absorption by PFA (< 0.03%), even samples at very low concentration can be stored for long periods of time in PFA containers without changes in their concentration (see the application example on Mercury Standards).

2,25 2,00 1,75 1,50 1,25 1,00 start 2 4 6 8 10 12 14 20 24 28 32 36 40 44 48 52

Concentration 2 ppb (ng/g) each

2,50


→ HG- not stabilized→ HG- stabilized (dichromate/HNO₃)

Application Example - Mercury Standards: Storage of an Hg standard in high-purity PFA containers (concentration 2 ppb (ng/g) each). Source: GIT Laboratory Trade Magazine 1/95

Cleaning after contamination

The valuable properties of PFA, primarily the near absence of memory effects, ensure the reliability of trace analytical results. In comparison with a commercially available glass flask, the simple cleaning after contamination with lead (Pb) solution at a concentration of 1000 ppb (ng/g) is a telling example (see the section on Cleaning of Volumetric Flasks). The cleaning of the glass and PFA volumetric flasks was conducted by shaking with 65% HNO₃ *Suprapur® (Pb < 0.005 ppm) at room temperature. With PFA volumetric flasks, the minimum concentration of 0.003 ppb is reached after three rinses, while substantially higher lead concentrations can be measured in glass flasks even after four rinses. The experiment also shows that PFA labware does not require the usual time-consuming boiling.

Mean of 4 Graduated Flask 500 ml each

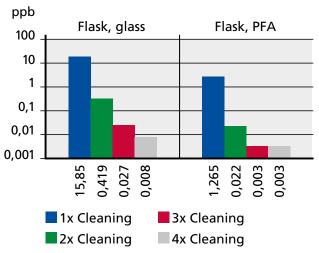

Cleaning of glass and PFA volumetric flasks after contamination. Source: Kali-Forschungsinstitut, K. Mangold

^{*} Suprapur® is a trademark of Merck KGaA.

VITLAB® UV-protect

Greater protection for light-sensitive substances.

Light shielding is an important factor for numerous laboratory applications (see Fields of Application). Accordingly, volumetric flasks made of amber glass are part of the basic equipment of a laboratory. Alternatively, coloured plastic products are used. Both solutions have disadvantages in practice: The glass flasks are fragile and most plastic products are completely opaque so that the filling level can not be read.

The UV-protect volumetric flasks are very close to the shade of amber glass, they are transparent and provide comparable protection for sensitive substances in the visible light spectrum (750 – 400 nm). In the UV range (380 – 200 nm), they shield even better.

For the user, this means easier sample preparation through improved readability of the volume mark combined with higher UV absorption to protect the sample contents. Added to this are the higher break resistance and the very good chemical resistance of polymethylpentene. As PMP also has high temperature stability, the volumetric flasks are also ideal for use in biological laboratories, where products are autoclaved regularly at 121° C before use.

This makes VITLAB® UV-protect the ideal alternative to conventional amber glass and conventional coloured plastic containers for most applications.

Fields of Application

In many fields of application, light-sensitive substances are used –, i.e. those which change or decompose through light. These include, for example, silver salts (silver chloride, silver nitrate) and iodine solutions. In laboratories, there are various sources of light (daylight, fluorescent tubes, UV lamps), from which these media must be effectively protected.

The fields of Life Sciences and Microbiology, for example, deal with living organisms that are sensitive to light, such as light-sensitive microorganisms or cell cultures. These can come from an environment in which there is little to no light – , for example, from the deep sea or the soil. In environmental analysis, sampling under light protection is also often very important –, for example, in water samples. The sample is taken on site, sealed, transported to the laboratory and analysed at a later time. In the chemical-analytical field, custom-made standard solutions should last as long as possible in order to be able to carry out exact comparative analyses. Light irradiation can change the composition of a sample and lead to incorrect results.

Classification and type description of plastics

In general, plastics can be divided into the three groups. The abbreviations of the described plastics used are those according to DIN 7728.


Thermoplastics

Plastics with a linear molecular structure, with or without side chains, can be reversibly moulded through heat treatment without changing their thermoplastic properties. Thermoplastics are frequently used in the production of plastic labware. Hence we provide here a brief description of some of the more important plastics, and explain their molecular structure, as well as their mechanical, chemical and physical properties. The most frequently used thermoplastics are polyolefins, such as polyethylene and polypropylene.

PFA (perfluoroalkoxy copolymer)




- Highly transparent, elastic thermoplastic with high-molecular, semi-crystalline structure
- Very good temperature stability
- Broad range of application, from -200 °C to +260 °C
- Virtually chemically inert, excellent chemical stability against practically all chemicals
- Very low water absorption (< 0.03%)
- Ultra-smooth, anti-adhesive surface with unique surface structure
- Typical products are, for example: Class A volumetric flasks, bottles, sample containers
 - ➡ Especially suitable for use in trace analysis, and for the storage of low-concentration solutions

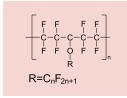

PTFE Polytetrafluoroethylene




- Non-transparent, white, elastic thermoplastic with high-molecular, semi-crystalline structure
- Very good temperature stability
- Broad range of application, from -200 °C to +260 °C
- Virtually chemically inert, excellent chemical stability against practically all chemicals
- Anti-adhesive surface
- Very good sliding properties and electrical insulating capability (very low friction coefficient)
- Typical products are, for example: Bottles, beakers, sheathing for magnetic stir bars

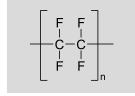

FEP Tetrafluoroethylene-perfluoropropylene copolymer



- Translucent, white, thermoplastic copolymer with high-molecular, semi-crystalline structure
- Anti-adhesive surface
- Very good temperature stability
- Broad range of application, from -100 °C to +205 °C
- Excellent chemical stability

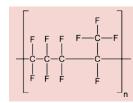

ETFE Ethylene-tetrafluoroethylene copolymer




- Translucent, white copolymer from ethylene and tetrafluoroethylene
- Very good temperature stability
- Broad range of application, from -100 °C to +150 °C
- Very good chemical stability
- Typical products are, for example: Thread adapters, Griffin beakers, threaded connectors



PMP Polymethylpentene



- Crystal-clear, stiff thermoplastic
- Similar structure to PP, with the methyl group replaced by an isobutyl group
- Good thermal stability
- Range of application from 0 to +150 °C
- Good tensile strength and dimensional stability
- Good chemical stability
- Typical products are, for example: Class A volumetric flasks, Class A measuring cylinders
 - → For the storage of light-sensitive substances, also available in highly transparent, UV-absorbing VITLAB® UV-protect

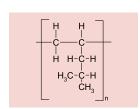


Classification and type description of plastics

PP Polypropylene




- Highly transparent, elastic thermoplastic
- Similar structure to PE, with a methyl group attached to alternate carbon atoms in an isotactic arrangement
- Good thermal stability
- Range of application from 0 °C to +125 °C
- Good tensile strength and dimensional stability
- Good chemical stability, comparable to PE
- Typical products are, for example: Class B volumetric flasks, Class B measuring cylinders, measuring pitchers, sample containers, funnels

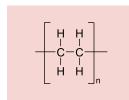


PE-HD High-density polyethylene

- Transparent, elastic thermoplastic
- Good thermal stability
- Range of application from -50 °C to +105 °C
- Compact, with increased tensile strength due to less cross-linking compared with PE-LD
- Good chemical stability
- Better chemical stability to organic solvents compared to PE-LD
- Typical products are, for example: Bottles, buckets, scoops
 - → For the storage of light-sensitive materials, also available in brown-dyed version

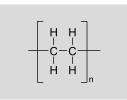

PE-LD Low-density polyethylene

- Highly transparent, elastic thermoplastic
- Moderate thermal stability
- Range of application from -50 °C to +80 °C
- · Very good flexibility
- Good chemical stability
- Typical products are, for example: Wash bottles, dropping pipettes

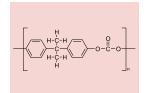


PC Polycarbonate

- Transparent, stiff thermoplastic
- Linear polymer of carbon dioxide



- Very good temperature stability
- Broad range of application, from -130 °C to +125 °C
- Good tensile strength and impact resistance
- Moderate chemical stability
- Typical products are, for example: Desiccators
- Note: Polycarbonates can lose their tensile strength through autoclaving or exposure to alkaline detergents

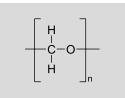

POM Polyoxymethylene

- Non-transparent, white, stiff, and high-molecular thermoplastic
- Good thermal stability



- Broad range of application, from -40 °C to +130 °C
 Use hardness and dimensional stability
- High hardness and dimensional stability
- Good sliding properties and abrasion resistance
- Good chemical resistance to aliphatic, aromatic, and halogenated hydrocarbons, and alkalis. Unstable to acids and esters
- Typical products are, for example: Slide boxes, staining jars
 - ⇒ Especially good chemical stability to organic solvents
 - → POM can replace metal in many applications

Classification and type description of plastics

PA Polyamide

- Linear polymers with regularly repeating amide bonds along the main chain
- Good thermal stability
- Range of application from -40 °C to +100 °C
- Outstanding durability and tensile strength, often used as construction material and for metal coatings
- Good chemical resistance to organic solvents
- Readily attacked by acids and oxidising agents
- Typical products are, for example: Spatula

SAN Styrene-acrylonitrile copolymer

- Crystal clear, stiff thermoplastic copolymer
- Moderate thermal stability
- Range of application from -40 °C to +70 °C
- Brittle and dimensionally stable
- Low tendency to form stress cracks
- Moderate chemical stability, SAN is slightly more chemically stable than PS
- Typical products are, for example: Graduated beakers, Class B graduated cylinders

PMMA Poly(methyl methacrylate)


- Crystal clear ("organic glass"), dimensionally stable thermoplastic
- Moderate thermal stability
- Range of application from -50 °C to +65 °C
- Very good UV radiation stability
- Poor chemical stability
- Typical products are, for example: Cuvettes


PS Polystyrene


- Crystal clear, stiff, amorphous or semi-crystalline thermoplastic
- Moderate thermal stability
- Range of application from -20 °C to +70 °C
- Hard, brittle, and dimensionally stable
- Tendency to form stress cracks
- Moderate chemical stability
- Typical products are, for example: Containers, cuvettes


PVC Poly(vinyl chloride)


- Amorphous thermoplastic, transparent with a slight blue tint
- Moderate thermal stability
- Range of application from -20 °C to +80 °C
- Good chemical stability, especially resistant to oils
- The addition of plasticisers opens up many useful applications, ranging from artificial leather to injection moulding components
- Typical products are, for example: Drawer organisers, dishes, trays

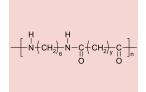

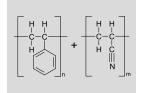

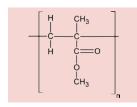

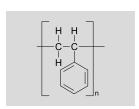

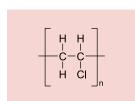

Classification and type description of plastics

Thermosets


Plastics with densely cross-linked molecules, which are very hard and brittle at normal temperatures. Heating causes irreversible hardening. These plastics are rarely used for plastic labware. The best known thermosets are the melamine resins. Melamine resin is produced by polycondensation of melamine with formaldehyde.

MF Melamine-formaldehyde resin

- Colourless thermoset, also belongs to the aminoplast group
- Good thermal stability
- Broad range of application, from -80 °C to +120 °C
- High surface hardness, abrasion resistance, and fire resistance
- Good electrical insulator, high creep resistance
- Good chemical stability
- Typical products are, for example: Trays, plates, mixing bowls
- Caution necessary when used in a microwave oven: Heating can release amounts of melamine and formaldehyde that can be harmful to health!


Elastomers

Plastics with loosely cross-linked molecules that exhibit rubber-like elasticity at room temperature. Heating causes irreversible integration (vulcanisation). The best known elastomers are natural rubber and silicone rubber.

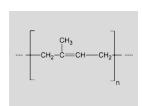


NR Natural rubber

- Elastomer obtained from latex (milk-like sap from rubber tree bark) and vulcanised with sulphur to improve the elasticity
- Composed of polymerised isoprene, with an extremely uniform structure
- Moderate thermal stability, not stable to UV
- Range of application from -40 °C to +80°C
- High tensile strength and elongation at rupture
- Poor chemical stability
- Typical products are, for example: pipette filler bulbs

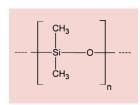

SI Silicone rubber

- Synthetic elastomer in which silicon atoms are linked together by oxygen atoms
- Includes polyorganosilanes that have groups such as hydrogen atoms, hydroxyl groups or vinyl groups for cross-linking reactions
- Very good thermal stability, and resistance to UV, ozone, and weather
- Broad range of application, from -60 °C to +180 °C
- High dimensional stability, even at high temperatures
- Very good durability, even at low temperatures
- Poor chemical stability

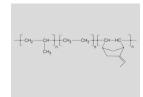

EPDM Ethylene-propylene-diene rubber

- Synthetic terpolymeric elastomer
- Manufactured with metallocene or Ziegler-Natta catalysts that utilise vanadium compounds and aluminium alkyl chlorides
- Good thermal stability
- Range of application from -40 °C to +130 °C
- High elasticity, even at low temperatures
- Stable to UV and ozone, and weather-resistant
- Very good chemical stability

Chemical resistance of plastics

With regard to chemical stability, plastics are classified as follows:



Very good chemical resistance

Within 30 days, continuous exposure to media causes no damage to the plastic. The plastic may remain resistant for years.



Good to limited chemical resistance

Within 7-30 days, continuous exposure to media causes minor damage (e.g. swelling, softening, loss of mechanical strength, discolouration), some of which is reversible.



Not suitable for continuous exposure to media. Immediate damage may occur (e.g. loss of mechanical strength, deformation, discolouration, cracking, liquification).

Chemical resistance of plastics to various substance classes

Classes of substances at 20 °C	PFA	PTFE	FEP	ETFE	PMP	ЬР	PE-HD	PE-LD	PC	POM	PA	SAN	PMMA	PS	PVC	MF	NR	S	EPDM	FKM
Alcohols, aliphatics	+	+	+	+	+	+	+	+	+	+	0	+	-	+	+	+	+	+	+	-
Ethers	+	+	+	+	-	0	0	0	-	+	+	-	-	-	-	-	-	-	-	-
Aldehydes	+	+	+	+	0	+	+	0	0	0	0	-	0	-	-	+	0	0	+	+
Esters	+	+	+	+	0	0	0	0	-	-	+	-	0	-	-	+	0	0	0	-
Hydrocarbons, aliphatic	+	+	+	+	0	+	+	0	0	+	0	-	+	-	+	+	-	-	-	0
Hydrocarbons, aromatic	+	+	+	+	-	0	+	0	-	+	0	-	-	-	-	+	-	-	-	0
Hydrocarbons, halogenated	+	+	+	+	-	0	0	0	-	+	0	-	-	-	-	+	-	-	-	0
Ketones	+	+	+	0	0	0	0	0	-	+	+	-	-	-	-	+	-	-	0	-
Alkalis	+	+	+	+	+	+	+	+	-	+	+	+	+	+	+	-	+	0	+	0
Acids, strong or concentrated	+	+	+	+	+	+	+	+	-	-	-	-	-	0	+	-	-	-	+	0
Acids, weak or diluted	+	+	+	+	+	+	+	+	0	-	-	0	-	0	+	0	0	0	+	+
Oxidising acids, oxidising agents	+	+	+	+	-	-	-	-	-	-	-	-	-	-	-	-	-	-	0	0

The carefully prepared recommendations listed here are based on technical literature and made available by the manufacturers of raw materials, to provide information and advice. However, nothing can replace suitability tests conducted by the end user under the actual conditions of use.

Medium	PFA	/FEP	РТ	FE	ET	ΓFE	PI	MP	F	PP	Н	OPE	LDPE	
Acataldahyda		50 °C	20 °C	50 °C										
Acetaldehyde	+	+	+	+	+	0	0	-	+	-	+	0	+	-
Acetic acid (glacial acetic acid), 100%	+	+	+	+	+	+	+	0	+	0	+	+	+	0
Acetic acid, 50%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Acetic anhydride	+	+	+	+	+	+	+	0	0	0	0	0	-	-
Acetone	+	+	+	+	+	0	+	+	+	+	+	+	+	0
Acetonitrile	+	+	+	+	+	+	0	-	+	0	+	0	+	0
Acetophenone	+	+	+	+	+	+	0	-	0	0	0	0	-	-
Acetyl chloride(acetic acid chloride)	+	+	+	+	+	+			+		+		+	
Acetylacetone	+	+	+	+	+	+	+		+		+		+	
Acrylic acid (2-propenic acid)	+	+	+	+	+	+	+		+		+		+	
Acrylonitrile	+	+	+	+	+	+	-	-	0	-	+	+	+	+
Adipic acid	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Allyl alcohol (2-propan-1-ol)	+	+	+	+	+	+	+	0	+	+	+	+	+	+
Aluminium chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Aluminium hydroxide	+	+	+	+	+	+	+	0	+	+	+	+	+	+
Amino acids	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Ammonium chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Ammonium fluoride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Ammonium hydroxide, 30% (ammonia)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Ammonium sulphate			+	+			+	+	+	+	+	+	+	+
n-Amyl acetate (pentyl acetate)	+	+	+	+	+	+	+	0	0	-	+	0	0	-
Amyl alcohol (pentanol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Amyl chloride (chloropentane)	+	+	+	+	+	+	-	-	-	-	-	-	-	-
Aniline	+	+	+	+	+	0	+	0	+	+	+	+	+	0
Aqua regia	+	+	+	+	+	+	-	-	-	-	-		-	-
Barium chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Benzaldehyde	+	0	+	+	+	0	+	+	+	+	+	+	+	+
Benzene	+	+	+	+	+	+	0	0	+	0	+	+	0	-
Benzoyl chloride			+	+	+	+	0	0	+	0	+	+	0	-
Benzyl alcohol	+	+	+	+	+	+	0	-	0	-	0	-	0	-
Benzyl chloride			+	+	+	+								
Benzylamine	+	+	+	+	+	+	0		0		0		0	-
Boric acid, 10%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Bromine	+	+	+	+	+	+	-	-	-	-	-	-	-	-
Bromobenzene	+	+	+	+	0	-	-	-	-	-	-	-	-	-
Bromoform	+	+	+	+	+	+	-	-	-	-	-	-	-	-
Bromonaphthalene	+	+	+	+	+	+								
Butanediol	+	+	+	+	+	+	+	+	+	+	+	+	+	+
1-Butanol (butyl alcohol)	+	+	+	+	+	+	+	0	+	+	+	+	+	+
n-Butyl acetate (acetic acid n-butyl ester)	+	+	+	+	+	+	+	0	0	0	+	+	0	0
Butyl methyl ether	+	+	+	+	+	0	+	-	+	0	0	-	0	-
Butylamine			+	+	+	+								
Butyric acid (butanoic acid)	+	+	+	+	+	+			-	-	0	-	-	-
Calcium carbonate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Calcium chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Calcium hydroxide	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Calcium hypochlorite	+	+	+	+	+	+	+	0	+	+	+	+	+	+

_	_										-				l		l	
	C	PC		P.			λN	PM		P		P∖		MF	NR	SI	EPDM	FKM
20 °C	50 °C	20 °C																
0	-	+	+	0		-	-	-	-	-	-	-	-		-	-	0	-
	-	-	-	-						-	-			0	0	0	0	-
+	0	0	-	-	-	+	0	-	-	0	0	+	0	+	-	-	-	-
	-	-	-	0	0					-	-	-	-		0	0	0	-
	-	+	+	+		-	-	-	-	-	-	-	-	+	0	-	+	-
	-	+		+		-	-	-	_	-	-	-	-		-	-	-	-
	-	+		+		-	-	-	-	-	-	-	-	+	-	-	+	-
-	-			-	-	-	-	-	-	-	-	-	-		-	-	-	+
-	-	+				-	-	-	-	-	-	-	-		-	-	+	-
-	-	-	-			-	-	-	-	-	-	-	-		-	-	-	-
-	-	-	-	+		-	-	-	-	-	-	-	-		-	-	-	-
+	+	+	+	+		+	+	+	+	+	+	+	0		+	+	+	+
0	0	+	+	0		0	-	-	-	0	0	0	-		0	-	+	+
-	-	+	0	0	-	+	+	+	+	+	+	+	0	+	0	0	+	+
0	-	+	+	+	+	0	0	0	0	0	0	+	+		+	+	+	+
+	+	+	+	+	+	+	+	+	+	+	+	+	+		+	+	+	+
0	0	+	+	+		+	+	0	0	+	+	+	0		+	+	+	+
0	0	+	+	+		+	+	0	0	+	+	+	0		-	+	+	0
-	-	0	0	0		+	0	+	+	0	-	+	0	+	+	0	+	-
+	+	+	+	+		+	+	+	+	+	+	+	+		0	0	+	-
_	-	+	+	0		-	-	+	+	-	-	_	-		0	-	0	-
+	+	+	+	+		+	+			0	0	0	0		0	-	0	0
	-	+	+	+		-	-	-	-	-	-	-	-		-	-	-	+
0	_	0	0	0		-	-	-	_	-	_	-	_		-	-	_	-
	_	-		-		-	-	-	_	-		-	-	-	-	-	-	-
+	+	+	+	+		+	+	+	+	+	+	+	+		+	+	+	+
0	_	+	+	0		-	-	-	_	-		-	-		-	-	0	-
	_	+	0	+		-	-	-	_	-	_	_	_	+	-	-	_	0
	_	+	0	_		-	-		_				_		_	-		+
0	0	+	+	_		-	_	_	_			0	0		_	0	0	+
		+		+											_	_		+
	-	+				-	-	-	-	-	-	-	-		-	0	0	+
+	+	+	+	+		+	+	+	+	+	+	+	+		+	+	+	+
		-		_		-	-	-		-		_			-	-	-	0
	-			+		-	-	-	_	-		-	-		-	-	-	+
	_	-		-		-		-	-	-		-	_		-	-		
						-				-								
		+	+	+		-									0	-	+	-
	0	+	+	+	0	+	0	0	_	0		0	0		+	0	0	+
	-	+	0	+		-	-	-		-		-	-			-	0	-
	-	+	+	•		-			_				_		_	_	-	-
		+	+			_									_	0		_
0	-			0	0	-									-	-	-	0
+	+	+	+	+	+	+	+	+	+	+	+	+	+		+	+	+	+
+	+	+	+			+	+	+	+	+	+	0			+	+	+	+
		+	+	+		+	0	+	+	+	0	+	+		+	0	+	+
		+	+	+		+	+	0		+	+	0	-		-	0	+	+

Medium	PFA	/FEP	PT	FE	ET	ΓFE	PI	MP	F	PP 9	Н	DPE	LC	PE
Carbon disulphido		50 °C	20 °C	50 °C										
Carbon disulphide		+	+	+	+	0		-	-	-	-	-	-	-
Carbon tetrachloride	+	+	+	+	+	+	-	-	-	-	0	-	0	-
Chloroacetaldehyde, 45%			+	+	+	+								
Chloroacetic acid (monochloroacetic acid)	+	+	+	+	+	+	+	0	+	0	+	+	+	+
Chloroacetone			+	+	+	+								
Chlorobenzene	+	+	+	+	+	0	-	-	-	-	-	-	-	-
Chlorobutane	+	+	+	+	+	+	0	-	0	-	0	-	0	-
Chloroform	+	0	+	+	+	0	0	-	-	-	0	-	0	-
Chloronaphthalene			+	+	+	+								
Chlorosulphonic acid	+	+	+	+	0	-								
Chromic acid, 10%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Chromic acid, 50%	+	+	+	+	+	+	0	0	0	0	+	0	+	0
Chromic-sulphuric acid	+	+	+	+	+	+	0	-	-	-	-	-	-	-
Copper sulphate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Cresol	+	+	+	+	+	0	-	-	0	0	0	-	-	-
Cumene (isopropylbenzene)	+	+	+	+	+	+	-	-	0	-	+	0	0	-
Cyclohexane	+	+	+	+	+	0	-	-	0	-	0	-	0	-
Cyclohexanone	+	+	+	+	+	+	0	0	0	-	0	-	-	-
Cyclopentane	+	+	+	+	+	+	0	-	0	-	0	-	-	-
Decane	+	+	+	+	+	+	0		0		0	-		
1-Decanol	+	+	+	+	+	+	+		+		+			
Di(ethylene glycol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Dibenzyl ether	+	+	+	+	+	+	0		+		+			
Dibromoethane	+		+	+	0									
Dibutyl phthalate	+	+	+	+	+	+	+	0	+	0	0	-	0	-
Dichloroacetic acid	+	+	+	+	+	0	+	+	0	-	0	0	0	-
Dichlorobenzene	+	+	+	+	+	0	-	-	0	-	0	-	0	-
Dichloroethane	+	+	+	+	+	+	0	-	0	-	0	-	0	-
Dichloromethane (methylene chloride)	+	+	+	+	0	0	0	-	0	-	0	-	0	-
Diesel oil (heating oil)	+	+	+	+	+	+	0	-	+	0	+	0	0	-
Diethanolamine			+	+					0		0			
Diethyl ether	+	+	+	+	+	+		-	0	-	0		-	-
Diethylamine	+	+	+	+	+	0	0	0	0	-	0	-	-	-
1,2 Diethylbenzene	+	+	+	+	+	0	-	-	-	-	0	-	-	-
Dimethyl sulphoxide (DMSO)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Dimethylaniline	+	+	+	+	+	+								
Dimethylformamide (DMF)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
1,4-Dioxane	+	+	+	+	+	0	0	0	+	0	+	+	+	0
Diphenyl ether			+	+										
1,2 Ethanediol (ethylene glycol, glycol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Ethanol (ethyl alcohol)	+	+	+	+	+	+	+	0	+	+	+	+	+	+
Ethanolamine	+	+	+	+	+	+			+					
Ethyl acetate (acetic acid ethyl ester)	+	+	+	+	+	+	0	-	+	0	+	+	+	+
Ethyl methyl ketone (MEK)	+	+	+	+	0	0	-	-	+	0	0	-	0	-
Ethylbenzene	+	+	+	+	0	0	-	-	-	-	-	-	-	-
Ethylene oxide	+	+	+	+	+	+	0	-	0	-	0	0	0	0
Fluoroacetic acid			+											

								ı										
	C	PC		P.		SA			MA	Р		P∖		MF	NR	SI	EPDM	FKM
20 °C	50 °C	20 °C	50 °C	20 °C	50 °C	20 °C	50 °C	20 °C	50 °C	20 °C	50 °C	20 °C	50 °C	20 °C				
_	-	+	+	0		-	-	-	_	-		-	-		-	-	-	+
	-	0	0	-	-	-	-	0	-	-	-	-	-	+	-	-	-	+
				-	-	-	-	-	-	-	-							
0	-	-	-	-	-	-	-	0	-	0	-	+	0		-	-	0	0
						-	-	-	-	-	-				0	-	+	-
_	-			-	-	-	-	-	-	-	-	-	-		-	-	-	0
		-	-			-	-	-	-	-	-				-	-	-	0
-	-	-	-	0	-	-	-	-	-	-	-	-	-	+	-	-	-	0
		-	-												-	-	-	+
		-	-	-	-										-	-	-	-
+	0	0	0	-	-	-	-	0	-	-	-	+	0		-	0	-	+
0	-	-	-	-	-	0	0	-	-	-	-	+	-		-	-	-	+
-	-	-	-	-	-	0	0	-	-	0	0	+	0		-	-	-	+
+	+	+	+	+		+	0	+	+	+	+	+	0		0	+	+	+
-	-			-	-					-	-	-	-		-	-	-	+
_	-	+	-			-	-	-	-	-	-	-	-		-	-	-	+
-	-	+	+	+						-	-	-	-	+	-	-	-	+
-	-			+						-	-	-	-		-	-	-	-
_	-									-	_	-	-		-	-	-	+
0		+										0			_	0	_	+
0		+				0				0		+			0	0	+	+
0	0	+	0	0		+	+	_		0		-			+	0	+	+
		+					<u> </u>	_	_						-	-	0	-
		<u>'</u>																
	_	+	+			_		-							-	0	0	0
0				_	-					0		0			_	-	-	-
-	_	_		+		_		_		-		-			_	_		+
		_		0											_	_		0
				0	-							_			_	_	_	0
		+	+	+	+			0	_			0			_		_	+
		Т						-				0			_		0	т
	-	+	+	+	+	-		-	_	-		-	-		-	-	-	_
		т		т						0	0	_			0	0	0	-
	-									-			-		-	-	-	
	-					_		_	-	-	-	-	-		_		-	+
															-	+ 0	0	0
	-					-	-	-	-	-	-	0					0	-
	- 0	+	+ 0	+		-	-	-	-	-		0	-		0	0	0	-
0			U	+		-	-	-	-	-	-	-	-		-	-		
		0		0		-	<u>-</u>	-	-	-	-				-	-	-	0
+	+	+	+	0	0	+	+	+	+	+	+	+	+	+	-	+	+	0
	0	+	+	+		0	-	-	-	0	-	+	0	+	0	0	+	0
															-	-	+	
-	-			+						-		-	-	+	-	-	0	-
	-	-	-	+		-	-			-	-	-	-		-	-		-
-	-					-	-	-	-	-	-	-	-		-	-	-	0
0	-	+	+	0		-	-	-	-	-	-	0	-		-	-	-	-
	-	-	-			-	-	-	-	-	-	-	-		-	-	-	-

Medium	PFA	/FEP	PT	ΓFE	ET	ΓFE	PI	MP	F	PP	Н	OPE	LC	PE
	20 °C	50 °C												
Formaldehyde, 40%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Formamide	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Formic acid, 98-100%	+	+	+	+	+	+	+	0	+	+	+	+	+	+
Gasoline (petroleum spirits)	+	+	+	+	+	+	0	0	0	0	+	+	0	-
Glycerine	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Glycolic acid, 70%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Heating oil (Diesel oil)	+	+	+	+	+	+	0	-	+	0	+	0	0	-
Heptane	+	+	+	+	+	+	0	0	0	0	0	0	0	-
Hexane	+	+	+	+	+	+	0	-	+	0	+	0	0	-
Hexanoic acid			+	+										
Hexanol	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydrochloric acid, 10%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydrochloric acid, 20%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydrochloric acid, 37%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydrofluoric acid, 40%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydrofluoric acid, 70%	+	+	+	0	+	+	+	0	+	0	+	0	+	-
Hydrogen bromide	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydrogen peroxide, 35%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Hydroiodic acid	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Isoamyl alcohol (3-methyl-1-butanol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Isobutanol (isobutyl alcohol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Isooctane	+	+	+	+	+	+								
Isopropanol (2-propanol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Isopropyl ether	+	+	+	+	+	0	-	-	-	-	-	-	-	-
Lactic acid (2-hydroxypropionic acid)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Lugol's solution (iodine/potassium iodide solution)	+	+	+	+	+	+	+	0	+	+	-	-	-	-
Mercury	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Mercury chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Methanol	+	+	+	+	+	+	+	+	+	+	+	+	+	0
Methoxybenzene	+	+	+	+	+	+								
Methyl butyl ether	+	+	+	+	+	0	+	0	+	+	0	-	-	-
Methyl ethyl ketone (MEK)	+	+	+	+	0	0	-	-	+	0	0	-	0	-
Methyl formate (formic acid methyl ester)	+	+	+	+	+	+								
Methyl propyl ketone	+	+	+	+	+	+	0	0	+	0	+	+	+	0
Methylene chloride (dichloromethane)	+	+	+	+	+	+	-	-	0	-	0	-	0	-
Mineral oil (motor oil)	+	+	+	+	+	+	+	+	+	+	+	+	+	0
Monochloroacetic acid	+	+	+	+	+	+	+	0	+	0	+	+	+	+
Nitric acid, 10%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Nitric acid, 30%	+	+	+	+	+	+	0	-	0	-	0	-	0	0
Nitric acid, 70%	+	+	+	+	+	+	-	-	-	-	-	-	-	-
Nitrobenzene	+	+	+	+	+	+	-	-	-	-	0	-	-	-
Oleic acid	+	+	+	+	+	+								
Oxalic acid	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Ozone	+	+	+	+	+	+	+	+	0	-	0	-	0	-
n-Pentane	+	+	+	+	+	+								
Peracetic acid	+	+	+	+	+	+								
Perchloric acid	+	0	+	+	+	+	0	-	+	-	+	-	+	-
Perchloroethylene	+	+	+	+	+	+	-	-	-	-	-	-	-	-

															ı			
	C		M		A		AN	PMI		Р			/C	MF	NR	SI	EPDM	FKM
20 °C	50 °C	20 °C																
+	0	+	+	+	0	+	+	-	-	-	-	0	-		0	0	+	0
		-	-	+											+		0	0
+	0	-	-	-	-	0	0			+	0	-	-	+	0	-	0	-
0	_	+	+	+		-	-	+		-	-	0	-	+	-	-	-	+
+	+	0	0	+	+	+	+	+	+	+	+	+	+	+	0	+	+	0
				-	-										+	+	+	0
	_	+	+	+		-	-	0	_	-	-	0	-		-	-	-	+
+	0			+				0	-	-	-	-	-		-	-	-	+
-	-	+	+	+		+	+	0	0	-	-	0	-		-	-	-	+
								+										
								+							0	0	-	+
-	-	-	-	-	-	0	-	0	-	+	+	+		-	0	0	+	+
0	0	-	-	-	-	0	-	0	-	+	+	0		-	0	-	+	+
-	-	-	-	-	-	0	-	0	-	0	0	0	-	-	0	-	+	0
-	-	-	-	-	-	+	0	-	-	+	+	0	-	-	-	-	0	0
-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-
+	+	-	-	-	-					0	-				0	-	0	+
+	+	+	-	-	-	+	+	-	-	+	+	+	0	0	-	0	0	
																	+	+
		+	+												0	0	0	0
+	+	+	+			0	-	0	-	0	0	+	0		+	+	+	+
0				+		0	-			0	-				-	-	-	+
+	+	+	+	+		+	-	0	-	0	0	+	0		+	0	+	+
_	-			-	-					-	-	-	-		-	-	-	-
+	+	+	-	0	-	+	+	0	-	+	+	0	0		0	0	0	+
0	-	0	0			0	-	-	-	0	-	-	-		+	-	+	+
+	+	+	+	+		+	+	+	+	+	+	+	+		+	+	+	+
+	+	0	0	-	-	+	+	+	+	+	0	-	-		+	+	+	+
+	0	+	+	0		0	-	-	_	0	-	+	0		0	+	+	-
_	-	0				-	-	-	-	-	-				-	-	-	-
	-	0				-	-	-	-	-	-	-	-		-	-	-	-
_	-	-	-			-	-	-	-	-	-	-	-		-	-	0	-
	-	+				-	-	-	_	-	-				-	0	0	
	_	+	+			-	_	_	_	-		-			-	-	0	-
	-	-	-	0	-	-	-	-	_	-	-	-	-		-	-	-	0
+		+	+			+		+	+	+		+	+		-	0	-	+
0	_	-		_		-		0		0		+	0		-	-		0
+	0	-	_	_		+	0	+	0	-		+	0	-	-	-	0	0
+	0	-				0		0	0	-		0		_	_	-	-	0
	-	-	_	-		-		-	-	-		-		-	-	-	-	-
-	_	0		-		-		-	_	-		-			-	-	-	-
				+		-	-	-	-	-	-				-	-	-	0
+	+	+	+	-		+	+	+	+	+	+	+	+		0	0	+	+
						0	0	+	0	0	0	+	0		-	+	+	+
				+														+
	_	-		•												-		
	_	-		_		-		_				0			_	-	0	+
	_	+	0	_		0	0	0	_			-			_	_	-	0

Medium	PFA	/FEP	PT	FE	ET	TE	ld.	MP	F	PP	Н	OPE	LC	PE
	20 °C	50 °C												
Petroleum	+	+	+	+	+	+	0	0	0	-	0	_	0	-
Petroleum ether	+	+	+	+	+	+							0	
Phenol	+	+	+	+	+	+	0	0	+	+	+	+	+	0
Phenylethanol	+	+	+	+	+	+			0		0			
Phenylhydrazine	+	+	+	+	+	+			0		0			
Phosphoric acid, 85%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Piperidine	+	+	+	+	+	+			+		+			
Potassium chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Potassium dichromate			+	+										
Potassium hydroxide	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Potassium permanganate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Propanol	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Propionic acid	+	+	+	+	+	0	+	0	+	0	+	0	0	-
Propylene glycol (propanediol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Pyridine	+	+	+	+	-	-	+	0	0	0	+	0	+	0
Salicylaldehyde	+	+	+	+	+	-	+	+	+	+	+	+	+	+
Salicylic acid	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Silver acetate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Silver nitrate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sodium acetate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sodium chloride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sodium dichromate	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sodium fluoride	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sodium hydroxide	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sulphuric acid, 60%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Sulphuric acid, 98%	+	+	+	+	+	+	+	+	-	-	0	-	0	-
Tartaric acid	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Tetrachloroethylene	+		+	+	0									
Tetrahydrofuran (THF)	0	0	+	+	+	0	0	-	-	-	0	-	0	-
Tetramethylammonium hydroxide	+	+	+	+	+	+								
Toluene	+	+	+	+	+	+	0	-	0	-	0	0	0	-
Tri(ethylene glycol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Tri(propylene glycol)	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Trichloroacetic acid	+	+	+	+	+	0	+	+	0	-	0	0	0	-
Trichlorobenzene	+	+	+	+	+	0	0	0	-	-	-	-	-	-
Trichloroethane	+	+	+	+	+	+	-	-	-	-	0	-	-	-
Trichloroethylene	+	+	+	+	+	+	-	-	-	-	0	-	-	-
Trichlorotrifluoroethane	+	+	+	+	0	-								
Triethanolamine	+	+	+	+										
Trifluoroacetic acid (TFA)	+	-	+	0										
Trifluoroethane	+	0	+	+										
Turpentine	+	+	+	+	+	+	0	0	-	-	0	-	0	-
Urea	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Xylene	+	+	+	+	+	+	0	-	-	-	0	-	0	-
Zinc chloride, 10%	+	+	+	+	+	+	+	+	+	+	+	+	+	+
Zinc sulphate, 10%	+	+	+	+	+	+	+	+	+	+	+	+	+	+

Р	C	PC	M	F	PΑ	SA	AΝ	PM	MA	Р	S	P۱	/C	MF	NR	SI	EPDM	FKM
20 °C	50 °C	20 °C																
0	0	+	+	+				+		-	-	+	-		-	0	-	+
		+	+	+				+		-	-	0	-		-	-	-	+
-	-	-	-	-	-	-	-	-	-	-	-	-	-		-	-	-	0
				0														
															0	-	-	0
+	+	+	-	-	-	+	+	-	-	+	0	+	0	-	-	-	0	+
															-	-	-	-
+	+	+	+	+		0	0	+	+	0	0	+	0		+	+	+	+
				-	-										0	0	+	0
-	-	+	+	+		0	0	+	+	0	0	0	0	-	0	-	+	-
+	+	0	0	-	-	+	0	+	+	+	+	+	+		-	-	+	+
0		+	+	+	+	+	+	0		0		+	+		+	0	+	+
-	-	-	-	0	0					0	-	0	-		-	-	0	+
+	0	+	+	-	-	-	-	0	0	+	+	0	-		+	+	+	+
-	-	+	0	+		-	-	-	-	-	-	0	-		-	-	-	-
0	0					-	-			-	-	-	-					
		-	-	+		+	+			+	+	0	-		+	+	+	+
+	+	0	0			0	0	0	0	0	0	0	0		+	+	+	+
+	+	0	0	+		+	+	+	+	0	0	0	0		+	+	+	+
+	+	+	0	+		+	+	-	-	+	+	0	0		+	0	+	-
+	+	+	+	+		+	+	+	+	+	+	+	+		+	+	+	+
+	-	+	+	+		+	0	+	0	+	0	+	+		+	0	+	+
+	+	+	+	+		+	+	+	+	+	+	+	+		0	0	+	+
-	-	+	+	+	0	+	+			+	+	+	+	-	0	0	+	0
0	0	-	-	-	-	+	0	-	-	-	-	0	-		-	-	-	+
-	-	-	-	-	-	-	-	-	-	-	-	-	-		-	-	-	+
+	+	+	+	0	0	+	+	0	0	+	+	+	+	-	+	+	0	+
				-	-										-	-	-	0
-	-	0	0	+		-	-	-	-	-	-	-	-	+	-	-	-	-
-	-	-	-														+	-
	-	+	+	+		-	-	-	-	-	-	-	-		-	-	-	0
+	0	+	0			+	+	0	0	+	+	0	-		0	+	+	+
+	0	+	0			+	+	0	0	+	+	0	-		+	+		
0	-			-	-					0	-	0	-		0		0	-
-	-									-	-	-	-					
-	-	0	-	0		-	-	-	-	-	-	-	-		-	-	-	+
	-	-	-	0		-	-	-	-	-	-	-	-		-	-	-	0
				0		-	-			-	-						0	
						-	-			-	-				0	-	0	-
		-	-			-	-			-								-
						-	-			-	-				-	-	-	+
_	-	+	+	+		0	0	+	+	-	-	+	+		-	-	-	+
	-	+	+	+		+	+	+	+	+	+	0	_	+	+	+	+	+
	-	+	+	+		-	-	-	-	-	_	-	_		-	-	-	0
+	+	+	0	-	-	+	+	-	-	+	+	+	0		+	+	+	+
+	+	0	-	-	-	+	+	0	0	+	+	+	0		0	+	+	+

Physical properties of plastics

Plastics	Max. temperature for use °C	Brittleness temperature °C	Microwaveability*	Density g/cm³
PFA	260	-200	yes	2,17
PTFE	260	-200	yes	2,17
FEP	205	-100	yes	2,15
ETFE	150	-100	yes	1,70
PMP	150	0	yes	0,83
PP	125	0	yes	0,90
PE-HD	105	-50	yes	0,95
PE-LD	80	-50	yes	0,92
PC	125	-130	yes	1,20
POM	130	-40	no	1,42
PA	90	0	-	1,13
SAN	70	-40	no	1,03
PMMA	65 - 95	-50	no	1,18
PS	70	-20	no	1,05
PVC	80	-20	no	1,35
MF	120	-80	yes**	1,50
NR	80	-40	no	1,20
SI	180	-60	no	1,10
EPDM	130	-40	-	
FKM	220	-30	-	-

^{*} Mind the chemical and temperature suitability!

^{**} Caution necessary when used in a microwave oven: Heating can release amounts of melamine and formaldehyde that can be harmful to health!

Cleaning and maintenance of plastics

All polyolefins, such as PE-LD, PE-HD, PP and PMP, as well as the fluoroplastics PTFE, PFA, FEP and ETFE have water-repellent surfaces that are very durable and easy to clean. For cleaning, depending on the level of contamination, commercially available neutral or alkaline detergents can be used. Polycarbonate (PC) labware should not be cleaned with alkaline cleaning agents (> pH 7). Please note that no scrubbing agents or scouring pads should be used for labware made from plastics.

Cleaning with dishwashers

Labware made from the above-named plastics (except for PE-LD, due to the temperature limit) can be cleaned and dried in a laboratory dishwasher together with other apparatus. Machine cleaning with laboratory dishwashers is gentler to labware than cleaning in an immersion bath. The labware is exposed to the cleaning solution for relatively short rinsing periods when sprayed by the jet or injector nozzles. Due to their light weight, we recommend securing the apparatus to be washed with washing nets to prevent them from tumbling in the water jet. Labware is better protected against scratching when the wire baskets in the washing machine are plastic coated.

Cleaning in an ultrasonic bath

Plastic labware may be cleaned in an ultrasonic bath. However, direct contact with the acoustic membrane should be avoided.

Cleaning in trace analysis

To avoid contamination with cations and anions in trace analysis, plastic labware should be allowed to stand with a 1N HCl or HNO_3 solution for a maximum of 6 hours at room temperature, and be rinsed afterwards with purified distilled water. For trace analysis conducted in the concentration range of ng/g (ppb) or pg/g (ppt), containers made of the fluoroplastic PFA are particularly suitable, because they have a smooth surface, are easy to clean without carry-over (memory effects) and interaction with the container material.

Sterilisation of laboratory equipment made of plastics

Autoclaving

Recommended autoclaving protocol

20 minutes at 121 °C (2 bar), according to DIN EN 285

Autoclaving (steam sterilisation) is defined as the destruction or irreversible inactivation of all reproducible microorganisms under exposure to "saturated steam at a minimum of 120 °C." (DIN 58946-1, 1987). DIN EN 285 specifies a minimum exposure time (t_e) of 20 minutes (killing time and safety margin) at a sterilisation temperature of 121 °C. For the correct sterilisation procedure, including biological safety (DIN EN 285), please contact your hygiene specialist.

Prior to autoclaving plastic labware, ensure that no soiling or residual contamination remains on the equipment. Otherwise, the residual contamination will bake on solidly during the autoclaving process. Even substances that have no effect on the plastic at room temperature can still lead to destruction of the plastic during the autoclaving process. Additionally, microorganisms might not be killed effectively if they are protected by the residual contamination.

Notes on autoclaving



- Containers with screw tops or stoppers must be open during autoclaving to allow for pressure equalisation
 - → Autoclaving of a closed container will lead to the deformation or destruction of the container
- Plastic labware should be **stood upright on a level surface** during autoclaving to avoid shape deformation.
 - → Plastic labware should not be laid on its side during autoclaving
- No mechanical stresses should be present during autoclaving
 - → For example, do not stack items
- · Do not autoclave any container that contains residual contamination or even rinsing agent
- · Not all plastics are resistant to steam sterilisation! For example, polycarbonate loses its tensile strength
 - → Mind the temperature limits for the plastics
 - → Autoclavable products are identified with a "121 °C" symbol in this catalogue

The surfaces of some plastics can be attacked by chemicals present during autoclaving, which can cause persistent clouding. Some transparent plastics can absorb minute quantities of steam, which can lead to reversible clouding. This clouding disappears upon drying, which can be accelerated through the use of a drying oven.

Note!

For gas sterilisation, dry heat, and prior to heating in a microwave oven, all closures and stoppers must be removed as well.

Heating plastics in microwave ovens

Many plastics are suitable for use in microwave ovens. More accurate information can be obtained from the Table "Physical Properties of Plastics" on page 138. In this connection, it is important to be mindful of the chemical and temperature stability of the various plastics, and to ascertain whether the particular article and its contents are compatible with the given temperature. When aggressive acids, alkalis, or solvents are to be heated, the use of fluoroplastics is recommended. It is very important to provide for adequate ventilation (e.g., fume hoods).

Prior to the use of plastic labware in a microwave oven, the closures and stoppers must be removed from the apparatus.

Sterilisation* of plastics

Plastics	Autoclave 121 °C, t _e 20 min according to DIN	Heated air 160 °C (dry)	Gas (Ethylene oxide)	Chemical (Formalin, ethanol)	β-/γ-radiation 25 kGy
PFA	yes	yes	yes	yes	no
PTFE	yes	yes	yes	yes	no
FEP	yes	yes	yes	yes	no
ETFE	yes	no	yes	yes	no
PMP	yes	no	yes	yes	yes
PP	yes	no	yes	yes	yes (limited)
PE-HD	no	no	yes	yes	yes
PE-LD	no	no	yes	yes	yes
PC	yes ¹⁾	no	yes	yes	yes
POM	yes ¹⁾	no	yes	yes	yes (limited)
PA	no	no	yes	yes	yes
SAN	no	no	yes	yes	no
PMMA	no	no	no	yes	yes
PS	no	no	no	yes	yes
PVC	no	no	yes	yes	no
MF	no	no	yes	no	no
NR	no	no	yes	yes	no
SI	yes	-	yes	yes	no
EPDM	yes	-	yes	yes	-
FKM	yes	-	yes	yes	-

^{*} Before sterilisation, labware must be carefully cleaned and rinsed with distilled water. Always remove covers from containers!

¹⁾ Frequent autoclaving may reduce mechanical stability!

Suitability of plastics for foodstuffs



The marked products comply with the lawful regulations of the German Consumer Goods Ordinance and/or Directives (EC) No. 1935/2004, (EC) No. 975/2009 and (EU) No. 10/2011 as amended.

In the testing for compliance with the threshold values for the global migration (or respectively, the specific migration threshold values), no determinations exceeded the allowed values. In addition, sensory testing found no olfactory and flavour-related impairments. The testing was implemented according to the 82/711/EEC and 85/572/EEC Guidelines by an independent, accredited institute.

All source materials used in the manufacturing of the products are listed in the German Consumer Goods Ordinance as at 20.12.2006, or respectively, Directive (EU) 10/2011, in accordance with the present attestation. Therefore, they represent permissible source materials in accordance with food law and may be used in the production of food commodities in accordance with the specified restrictions concerning migration threshold values and permissible residual content in the end product.

Marked PP products are suitable for contact with all foodstuff categories providing that a contact period of 24 hours and a contact temperature of 40 °C are not exceeded. Marked SAN products are suitable for contact with all aqueous, alcoholic and fatty foods, providing that a contact period of 24 hours and a contact temperature of 40 °C are not exceeded.

Disposal and Recycling of Plastics

If the disposal of a piece of plastic labware is unavoidable, regional laws and regulations must be observed. Recycling centres can be found in many cities, which are designed for the disposal of recyclable materials. To simplify the task of sorting in these recycling centres, the majority of labware from VITLAB can be easily identified and presorted by using the engraved or imprinted recycling code. Prior to disposal, plastic labware must be cleaned, and sterilised if necessary, according to the currently valid regulations.

To simplify the separation of plastics for recycling, so that these can later be reused as raw material for manufacture, plastic identification labels (number 01-07) have been introduced. This identification scheme was published in 1988 under the title "SPI resin identification coding system" by the Society of the Plastics Industry (SPI). For coding, the commonly used shorthand for plastics according the DIN 7728 is also used.













SPI number 07 stands for "other". This is used to indicate other plastics such as PMP, PFA, PTFE etc. VITLAB doesn't use the "0"; rather, it identifies the specific raw material with the abbreviation according to DIN 7728 to simplify identification of the plastic for the end user.

CE mark / CE-IVD Guidelines

IVD Guidelines of the EU

On 7 December, 1998, the EU "Guidelines for In Vitro Diagnostic Devices" (IVD Guidelines) were published in the Official Journal of the European Communities, and thus came into force. The Guidelines were transposed into German National Law on 1 January, 2002, as a corresponding modification of the German Medical Devices Act (MPG). Consequently, in vitro diagnostic devices are considered medical devices.

Definition: Medical devices*

Medical devices are all instruments, apparatus, devices, materials, or other objects including software that are intended by the manufacturer for use in humans:

- for the purpose of detection, prevention, monitoring, treatment, alleviation or compensation of diseases, injuries or disabilities;
- for the purpose of investigation, replacement or modification of the anatomy or of a physiological process;
- for the purpose of control of conception. Pharmacologically or immunologically active agents are excluded, as these are regulated by the German Pharmaceuticals Law.

Definition: In vitro Diagnostic Devices (IVD)*

"In vitro diagnostic devices" are medical devices that are used for in vitro investigations of samples derived from the human body, including donated blood and tissue. Included are reagents, calibration substances or devices, control substances or devices, equipment, instruments, apparatus, systems, or also sample containers, if they are specifically intended by the manufacturer for use in medical tests. "In vitro diagnostic devices" serve mainly to provide information on:

- physiological or pathological conditions;
- · congenital anomalies;
- monitoring of therapeutic measures.

CE Mark

With the CE mark on a product, the manufacturer affirms that the product complies with the requirements for products of that type established by the EU Guidelines and, as necessary, has undergone the required testing. The manufacturer applies this mark to the product and additionally produces a Conformity Declaration that certifies the conformity of the product with the cited guidelines and standards.

The medical products supplied by VITLAB are all included in the class of in vitro diagnostic (IVD) devices.

This includes, for example:

- Pipette tips
- Urine bottles

^{*} See the definitions according to MPG § 3 (Definition of Terms)

Accuracy

What do "tolerance, accuracy, coefficient of variation, and precision" mean in volumetric measurements?

Graphic illustration of precision and accuracy

The dart board simulates the volume range around the centred nominal value, the white dots simulate the different measured values of a specified volume.

Good accuracy: All hits are near the centre, i.e., the nominal value.


Good precision: All hits are close together. **Result:** The manufacturing process is well controlled by an accompanying quality assurance program. Minimal systematic deviations and a narrow variance in products. The permissible limits are not exceeded. There are no rejects.



Good accuracy: On average, the hits are evenly distributed around the centre. **Poor precision:** No substantial errors,

but hits widely scattered.

Result: All deviations are "equally probable". Instruments exceeding the permissible tolerance should be rejected.

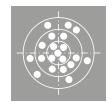


Poor accuracy: Although all hits are close together, the centre (nominal value) is still missed.

Good precision: All hits are close together. **Result:** Improperly controlled production, with systematic deviation. Instruments exceeding the permissible tolerance should be rejected.



Poor accuracy: The hits are far removed from the centre.

Poor precision: The hits are widely scattered. **Result:** These volumetric instruments are of inferior quality.

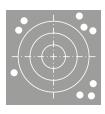

Calculation formulae

The accuracy of glass volumetric instruments is commonly defined by "Tolerance Limits", whereas for liquid handling instruments the statistical terms "Accuracy [%]" and "Coefficient of Variation [%]" have been established.

Tolerance


The term "tolerance" (tol.) in the corresponding standards defines

the maximum permissible deviation from the nominal value.



Accuracy

Accuracy (A) indicates the closeness of measured mean volume to the

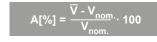

nominal value, i.e., systematic measurement deviation. Accuracy is defined as the difference between the measured mean volume (\overline{V}) and the nominal value $(V_{nom.})$, related to the nominal value in percent.

Coefficient of Variation

The coefficient of variation (CV) indicates the closeness of values

 $CV[\%] = \frac{s \cdot 100}{\overline{V}}$

from repeated measurements, i.e., random measurement deviation. The coefficient of variation is defined as standard deviation in percent, related to the mean volume.

Partial volumes

(analogous to CV_{part.} %) Generally, A and CV are based on



the nominal volume ($V_{nom.}$). These data in percent must be converted to partial volumes ($V_{part.}$). By contrast, there is no conversion for partial volumes if A and CV are stated in volume units (e.g. ml).

Tolerance from A and CV

To a good approximation, the tolerance, e.g. for the nominal



volume $(V_{nom.})$, can be calculated from the accuracy and coefficient of variation.

Precision

If the variance in the individual measurement results about the mean volume \overline{V} is given in units of volume, this relates to precision.

Cat. No. Index

33331	23	53510	73	65004	42	71598	115	80343	94
33332	23	53610	73	65014	42	71698	115	80345	94
33333	23	56695	79	65091	44	71798	115	80346	94
36491	112	56795	79	65104	42	71898	115	80347	94
39194	58	56895	79	65114	42	71998	115	80348	94
39294	58	56995	79	<u>65191</u>	44	72098	115	80375	110
39394	58	57095	79	65204	42	72198	115	80408	107
39494	58	60503	76	65214	42	72298	115	80409	107
39594	58	60603	76	65291	44	72398	115	80410	107
39694	58	60695	76	65304	42	72498	115	80411	107
39794	58	60703	76	65414	42	72898	83	80412	107
39894	58	60795	76	65965	90	72998	83	80413	107
39994	58	60803	76	65975	90	73498	83	80437	84
40093	59	60895	76	65980	90	73598	83	80438	84
40193	59	60903	76	67104	39	73698	83	80439	84
40293	59	60995	76	67195	40	73898	83	80440	84
40393	 59	61003	76	67204	39	75991	87	80441	84
40493	 59	61095	76	67295	40	76299	116	80442	84
40593	59	61103	76	67304	39	77094	85	80443	84
40693	59	61195	76	67395	40	78294	86	80445	84
40793	59	61203	76	67404	39	78394	86	80452	78
40894	61	61295	76	67495	40	78593	114	80454	78
41094	61	61403	76	67504	39	 78794	114	80455	78
41194	61	61495	76	67595	40	79194	50	80456	78
41294	61	61503	76	67604	39	79790	115	80550	90
41394	61	61595	76	67695	40	79890	115	80551	90
41494	61	61603	76	67704	39	79990	115	80553	90
41594	61	61695	76	67795	40	80130	96	80554	90
41694	61	61703	76	67895	60	80131	96	80555	91
41794	62	61795	76	67995	60	80134	96	80556	91
41894	62	61803	76	68099	60	80139	33	80557	91
41994	62	61895	76	68199	60	80140	33	80560	96
42099	62	64091	45	68299	60	80162	84	80562	96
42294	62	64191	45	68399	60	80164	84	80593	60
42393	62	64291	45	68594	93	80165	84	80594	60
42594	112	64391	45	68894	93	80213	116	80595	60
42694	112	64491	45	69194	93	80215	51	80596	60
42794	112	64591	45	69294	93	80217	51	80602	111
42894	112	64604	42	69394	78	80218	51	80603	111
42994	112	64614	42	69493	78	80219	51	80604	111
43094	112	64704	42	70494	62	80221	51	80730	93
43194	112	64714	42	70594	62	80222	51	80731	93
43610	111	64804	42	70694	62	80229	91	80732	93
44091	56	64814	42	70794	61	80230	90	80733	93
44191	56	64891	44	70894	61	80231	90	80734	93
44291	56	64904	42	70994	61	80271	88	80735	93
44391	56	64914	42	71094	61	80340	94	80736	93
44491	56	64991	44	71194	61	80342	94	80737	93

80828	60	92889	107	100389	108	106599	32	112797	75
81056	87	93389	105	100394	109	106699	32	112897	75
81213	116	93489	105	100489	108	106799	32	113197	78
81219	51	93788	69	100494	109	106899	32	113297	78
81640	110	93789	105	100589	108	107097	38	113397	78
81642	110	93793	68	100594	109	107197	38	113497	78
81644	110	93794	104	100689	108	107297	38	122097	82
81646	110	93988	69	100694	109	107397	38	122197	82
81660	110	93989	105	100789	108	107497	38	122297	82
81662	110	93993	68	100889	108	107597	38	122597	104
81664	110	93994	104	100989	108	107797	88	122697	104
81666	110	94188	69	101589	108	107897	88	122797	104
83300	70	94189	105	101594	109	107997	88	122897	104
83301	70	94193	68	101689	108	108297	102	130294	93
83302	70	94194	104	101694	109	108397	102	130297	92
83303	70	94587	71	101789	108	108497	102	130394	93
83304	70	94588	69	101794	109	109297	102	130397	92
83306	71	94589	105	101889	108	109397	102	130494	93
83307	71	94687	71	101894	109	109497	102	130497	92
83308	71	94688	69	101989	108	109597	102	130594	93
83310	105	94689	105	102089	108	109697	102	130597	92
83311	105	94987	71	102189	108	109797	102	131097	114
83312	105	94988	69	102697	103	110194	101	132193	71
83313	105	94989	105	102897	103	110204	75	132293	71
83314	105	94993	68	103297	87	110205	74	132393	71
83315	105	94994	104	103397	87	110294	101	132493	71
83316	105	95087	71	103797	94	110304	75	132593	71
83317	105	95088	69	103897	94	110305	74	132603	69
83330	101	95089	105	104097	95	110394	101	132605	69
83340	79	95093	68	104099	50	110404	75	132606	69
83350	79	95094	104	104197	95	110405	74	132608	69
90694	109	95187	71	104199	50	110494	101	132693	71
90794	109	95188	69	104297	95	110594	101	132703	69
90894	109	95189	105	104397	95	110604	75	132705	69
90994	109	95193	68	104497	95	110605	74	132706	69
91094	109	95194	104	104597	95	110704	75	132708	69
91194	109	95286	73	104997	95	110804	75	132793	71
91294	109	96093	57	105097	95	110904	75	132803	69
91394	109	96293	57	105197	95	110905	74	132805	69
91494	109	96393	57	105297	95	111004	75	132806	69
91789	107	96593	57	105597	95	111005	74	132808	69
91989	19, 107	96694	57	105599	32	111194	101	133181	68
92089	19, 107	96794	57	105697	95	112197	75	133281	68
92189	19, 107	97791	56	105699	32	112297	75	133381	68
92194	109	97891	56	105799	33	112397	75	134293	70
92489	19, 107	99099	98	105899	33	112497	75	134393	70
92689	107	99199	98	106399	32	112597	75	134493	70
92789	107	99299	98	106499	32	112697	75	134593	70

Cat. No. Index

135181	68	159670	85	307797	81	395940	58	569941	79
135281	68	159697	89	307897	81	395950	58	570941	79
135381	68	160110	47	307997	81	396940	58	605081	77
136693	98	160119	47	308097	81	396943	58	606081	77
138093	106	160210	47	308197	81	396944	58	606941	77
138193	106	160219	47	308297	81	396946	58	607081	77
138293	106	160510	47	308397	81	396950	58	607941	77
138393	106	160519	47	308497	81	396952	58	608081	77
138493	106	161010	47	308597	81	396955	58	608941	77
138593	106	161019	47	308897	81	396956	58	609081	77
138693	106	162510	47	308997	81	397940	58	609941	77
138793	106	163094	46	309097	81	397950	58	610081	77
138893	106	163194	46	309297	81	406944	59	610941	77
139393	106	163294	46	309397	81	406950	59	611081	77
139493	106	163394	46	310197	82	407944	59	611941	77
139593	106	163594	46	310297	82	407950	59	612081	77
139693	106	164094	46	310397	82	440081	54	612941	77
139793	106	164194	46	310497	82	440941	54	614081	77
139893	106	164294	46	311097	80	441081	54	614941	77
146294	27	164394	46	311197	80	441941	54	615081	77
146494	27	164494	46	311297	80	442081	54	615941	77
148893	72	164594	46	311397	80	442941	54	616081	77
148894	26	165094	112	311497	80	443081	54	616941	77
148993	72	165194	112	311597	80	443941	54	617081	77
148994	26	165294	112	311697	80	444081	54	617941	77
149093	72	165394	112	311797	80	444941	54	618081	77
149193	72	165494	112	316097	81	445081	54	618941	77
149194	27	165594	112	316197	81	445941	54	640941	44
149293	72	300497	80	316297	81	446081	54	641941	44
149393	72	300597	80	316397	81	446941	54	642941	44
149693	72	300897	80	316497	81	447081	54	643941	44
149794	26	301097	80	318293	83	447941	54	644941	44
149893	72	301197	80	318393	83	480941	55	645941	44
149993	72	301597	80	318597	82	481941	55	646081	43
149994	26	301697	80	324594	50	481942	55	646941	43
150194	27	301797	80	324694	50	481943	55	647081	43
150294	27	301897	80	324794	50	481944	55	647941	43
150394	27	301997	80	325095	84	481945	55	648081	43
151594	111	<u>3</u> 02097	80	325195	84	482941	55	648941	43
155400	26	302197	80	325295	84	482942	55	649081	43
155494	26	302297	80	325395	84	482943	55	649941	43
155600	26	<u>3</u> 02397	80	325495	84	482944	55	650081	43
155694	26	302497	80	326094	90	482945	55	650941	43
155900	27	302597	80	326194	90	483941	55	651081	43
155994	27	303097	80	326294	90	484941	55	651941	43
159497	89	303197	80	326496	89	566941	79	652081	43
159597	89	303297	80	326596	89	567941	79	652941	43
159665	85	307697	81	326696	89	568941	79	653081	43

653941	43	834331	67	1352869	66	1625506	12	1671683	17
670950	39	834332	67	1352879	66	1625507	12	1672000	24
671891	41	834341	67	1352889	66	1625508	12	1672002	24
671895	40	834342	67	1352899	66	1627515	14	1672010	24
671941	41	834351	67	1352909	66	1627525	14	1672012	24
671950	39	834352	67	1352939	66	1627535	14	1672015	24
672891	41	834413	67	1352949	66	1627545	14	1672016	24
672895	40	952861	73	1352959	66	1630500	49	1678132	16
672941	41	978081	56	1352969	66	1631500	48	1678134	16
672950	39	978941	56	1352979	66	1641000	22	1678136	16
673891	41	1037979	94	1353819	66	1641002	22	1678210	17
673895	40	1038971	94	1353839	66	1641004	22	1678212	17
673941	41	1327111	69	1353849	66	1641006	22	1678214	17
673950	39	1328111	69	1353859	66	1641008	22	1678216	17
674891	41	1331819	65	1353869	66	1641010	22	1678218	17
674895	40	1331839	65	1431829	65	1641012	22	1678220	17
674941	41	1331849	65	1432829	65	1670066	16, 31	3125970	82
674950	39	1331859	65	1432989	65	1670067	16, 31	3126970	82
675891	41	1331869	65	1433829	65	1670068	16, 31	3190940	97
675895	40	1332819	65	1433959	65	1670072	16	3190943	97
675941	41	1332839	65	1451829	66	1670080	16, 31	3190948	97
675950	39	1332849	65	1452829	66	1670100	16, 31	3191940	97
676891	41	1332859	65	1452989	66	1670105	16	3191943	97
676895	40	1332869	65	1453829	66	<u>1670110</u>	16, 31	3191948	97
676941	41	1332879	65	1453959	66	1670115	16, 31	3192940	97
676950	39	1332889	65	1610501	18	1670120	16, 31	3192943	97
677891	41	1332899	65	1610502	18	1670125	16, 31	3192948	97
677895	40	1332909	65	1610503	18	1670145	19	3193940	97
677941	41	1332939	65	1610504	18	1670150	16	3193943	97
760941	45	1332949	65	1610506	18	1670155	16, 31	3193948	97
834102	67	1332959	65	1611503	18	1670180	16, 31	3194940	97
834111	67	1332969	65	1611506	18	1670647	49	3194943	97
834112	67	1332979	65	1611508	18	1670648	49	3194948	97
834121	67	1333819	65	1620506	30	1670650	49	3195940	97
834122	67	1333839	65	1620507	30	1670660	49	3195943	97
834131	67	1333849	65	1621503	13	1670702	17	3195948	97
834132	67	1333859	65	1621504	13	<u>1671085</u>	31	3197940	97
834141	67	1333869	65	1621505	13	1671090	15	3197943	97
834142	67	1351819	66	1621506	13	1671095	31	3197948	97
834151	67	1351839	66	1621507	13	<u>1671116</u>	15	3198940	97
834152	67	1351849	66	1621508	13	<u>1671500</u>	15, 31	3198943	97
834213	67	1351859	66	1622502	13	1671505	19	3198948	97
834302	67	1351869	66	1622504	13	1671506	19	3961111	58
834311	67	1352819	66	1622505	13	<u>1671510</u>	15, 31	4811111	55
834312	67	1352839	66	1625503	12	<u>1671515</u>	15	4821111	<u>55</u>
834321	67	1352849	66	1625504	12	1671520	<u>15</u>		
834322	67	1352859	66	1625505	12	1671682	17		

Alpha - Index

Α		E	
Accessories for VITLAB pipeo® & maneus®		Erlenmeyer flasks	7
Accessories for VITLAB® microlitre pipettes		Evaporating dishes, PFA	8
Adapter for VITLAB® continuous E/RS		_	
Adapter for VITLAB® dispensers		F	
Adapter for VITLAB® piccolo		Filling tubes for VITLAB® dispensers	
Analytical funnels, PP		Filtering rack	
Autosampler-vials, PFA	95	Flexible discharge tubes for VITLAB® dispensers	
В		Forceps	
Beakers	56 74 79	Funnels	
Bottle-top burette VITLAB® continuous		Funnels, PP	6
Bottle-top dispenser VITLAB® genius²		G	
Bottle-top dispenser VITLAB® piccolo		Gas wash bottles, PFA	Q
Bottle-top dispenser VITLAB® simplex² / simplex _{fix}		Graduated beakers, coloured, PP	
Bottle-top dispenser VITLAB® TA ²		Graduated beakers, PPGraduated beakers, PP	
Bottles		Graduated beakers, PPGraduated beakers, SAN	
Bottles for VITLAB® continuous E/RS		Graduated beakers, stackable, PP	
Bottles for VITLAB® piccolo, PE-HD		Graduated cylinders, PMP, Class A, tall shape	
Bottles, brown glass, for VITLAB® genius² / simplex²		Graduated cylinders, PP, Class B, short shape	
Bowl, PP, with lid		Graduated cylinders, PP, Class B, silort shape	
Bowls, PP	112	Graduated cylinders, FF, Class B, tall shape	
Brown glass bottles for VITLAB® genius² / simplex²	15	Graduated cylinders, SAN, Class B, short shape	
Büchner funnels, PP	84	Graduated cylinders, 3AN, Class B, tall shape	
Bucket with spout, PP	57	Graduated pipettes. PP, Class B	
Buckets, PE-HD		Griffin beakers, ETFE	
Bulb pipettes, PP, Class B		Griffin beakers, PFA	
Burette clamps, PP		Griffin beakers, PMP	
Burette stopcocks, PMP/PTFE		Griffin beakers, PP	
Burette VITLAB® continuous		Griffin beakers, PTFE	
Burettes, borosilicate glass 3.3		Gilliii beakers, i ii E	/
Burettes, VITLAB® Dr. Schilling	32	Н	
C		Hydrometer cylinder, PP	4
Cans	92 - 95		
Caps with dropper inserts, PE-HD		I	
Catchment trays, PP		Imhoff cone, SAN	8
Chemical waste disposal system, PE/PP		1	
Collectors, PP or SAN		L	
Coloured wash-bottles, PE-LD/PP		Lab bottles VITgrip™, PP, GL 45	
Connectors	114	Laboratory trays, PP	
Connectors, universal PP	114	Large funnels, PE-HD	
Coplin staining chamber, PP	98	Large funnels, PP	
_		Lids for trays, PS	11
D		M	
Desiccator plates, PP and porcelain			0
Desiccators with stopcock, PC		Magnetic stir bar, triangular, PTFE Magnetic stir bars	
Desiccators with stopcock, PP/PC		Magnetic stirring bar retrievers, flexible, PTFE	
Desiccators, PP/PC		Magnetic stirring bar retrievers, rexibie, FTFE	
Discharge tubes, flexible for VITLAB® dispensers		Magnetic stirring bar retrievers, PTFE	
Dispenser VITLAB® genius²		Magnetic stirring bars, barbell, PTFE	
Dispenser VITLAB gerius		Magnetic stirring bars, cross shape, PTFE	
·		Magnetic stirring bars, cross shape, FTFE	
1)ishansar VIII AR simpley / simpley		agricus suring sais, octagonai, i il L	
Dispenser VITLAB® simplex² / simplex² _{fix}		Magnetic stirring bars, oval PTFF	, a
Dispenser VITLAB® TA ²	14	Magnetic stirring bars, oval, PTFE	
Dispenser VITLAB® TA²	14 17	Magnetic stirring bars, polygonal, PTFE	8
Dispenser VITLAB® TA²	14 17 47	Magnetic stirring bars, polygonal, PTFE maneus [®]	8ı
Dispenser VITLAB® TA²	14 17 47	Magnetic stirring bars, polygonal, PTFE maneus [®] Measuring cup, PP	8 4 7
Dispenser VITLAB® TA² Dispensing cartridge for VITLAB® TA² Disposable pipettes, PS Double spatula. Dropping bottles, PE-LD	14 17 47 60	Magnetic stirring bars, polygonal, PTFE	8 4 7 5
Dispenser VITLAB® TA²	14 17 47 60 71	Magnetic stirring bars, polygonal, PTFE maneus® Measuring cup, PP Measuring scoops, coloured, PP Measuring scoops, PP	
Dispenser VITLAB® TA² Dispensing cartridge for VITLAB® TA² Disposable pipettes, PS Double spatula Dropping bottles, PE-LD. Dropping bottles, PE-LD/PE-HD	14 17 60 71 71	Magnetic stirring bars, polygonal, PTFE	
Dispenser VITLAB® TA² Dispensing cartridge for VITLAB® TA² Disposable pipettes, PS Double spatula Dropping bottles, PE-LD Dropping bottles, PE-LD/PE-HD Dropping pipettes		Magnetic stirring bars, polygonal, PTFE maneus® Measuring cup, PP Measuring scoops, coloured, PP Measuring scoops, PP Microlitre pipettes	

N		Stopper, PP, Standard joint	10)9
Narrow-mouth bottles	102 - 107	Storage bottles, PE-HD, with tap		
Narrow-mouth bottles, PE-HD	107	Storage bottles, PE-HD, without tap	11	0
Narrow-mouth bottles, PE-LD	105, 106	T		
Narrow-mouth bottles, PFA		T		_
Narrow-mouth bottles, PP	104	Tap for storage bottles, PP		
Non-return valves	114	Telescopic filling tube for VITLAB® continuous E/RS		
		Telescopic filling tubes for VITLAB® dispensers		
0		Test tube racks, coloured, PP		
O-rings for desiccators	91	Threaded bottles for VITLAB® continuous E/RS		
_		Tip-Box		
P		Transport container, PE-HD		
Pasteur pipettes	72	Trays, MF		
Pestles, MF		Tubing connectors	11	4
pipeo [®]	48	U		
Pipette baskets	51			
Pipette fillers		Urbanti funnels, PMP	8	54
Pipette jars	51	V		
Pipette stands	50	-		٠,1
Pipette tips, PP	26 - 27	Valve for desiccators		
Pipette washer	51	VENT-CAP wash bottle caps		
Pipettes	46, 47, 72	Ventilation plug for micro filter for VITLAB® genius² / simplex²		
Pipettes, Piston-operated pipettes	22	VITGrip TM lab bottles, PP, GL 45		
Piston-operated pipettes	22	VITLAB maneus®		
Plastic stand for VITLAB® dispensers	15	VITLAB maneus® pipette helper		
Powder funnels, PP	61	VITLAB pipeo®		
PTFE-tape	114	VITLAB pipeo® pipette helper		
		VITLAB® continuous E/RS		
R		VITLAB® Dr. Schilling burettes		
Reagent bottles, PP	108, 109	VITLAB® genius²		
Reagent tube racks	96, 97	VITLAB® micropipette Starter-Sets		
Round-bottom flask stand	88	VITLAB® micropipettes		
Round-bottom flasks, PFA	88	VITLAB® piccolo		
		VITLAB® simplex² / simplex _{fix}		
S		VITLAB® TA ²		
Safety wash bottles, VITsafe™		VITLAB® UV-protect volumetric flasks, PMP, Class A		
Sample containers	92 - 95	VITsafe™ safety wash bottles, narrow-mouth		
Sample containers, PE-LD	93	VITsafe™ safety wash bottles, wide-mouth		
Sample containers, PFA	92	Volumetric flasks, PFA, Class A		
Sample containers, PP		Volumetric flasks, PMP, Class A		
Sample tubes, PFA	94	Volumetric flasks, PMP, Class B		10
Sample vials, PFA	95	Volumetric flasks, PP, Class B	4	11
Scoops	58, 59	187		
Scoops, PE-HD	59	W	_	
Screw caps, GL 45, PP	79	Wash bottle caps, PP		
Screw caps, PFA	103	Wash bottle caps, VENT-CAP		
Screw caps, PP		Wash bottles with imprint, PE-LD/PP		
Sealing ring for the valve block for VITLAB® dispensers	17	Wash-bottles		
Sealing tape, PTFE	114	Wash-bottles, PE-LD		
Sedimentation cone, SAN		Wash-bottles, PE-LD/PP		
Sedimentation rack	87	Wash-bottles, PP		
Spatula		Wash-bottles, VITsafe™		
Spatula spoon		Watch glasses, PP		
Spray bottles		Watch glasses, PTFE	7	8'
Spray bottles, PP		Water-jet pump, PP	88	35
Staining chamber for slides, POM		Weighing jars, PP	9)4
Staining jar for slides, POM		Wide-mouth bottles	102 - 10)7
Staining rack for slides, POM		Wide-mouth bottles, PE-HD	10)7
Standard joint funnels, PP		Wide-mouth bottles, PE-LD	105 - 10)7
Standard joint stoppers, PP		Wide-mouth bottles, PE-LD, with eye closure	10)7
Starter-Sets VITLAB® micropipette		Wide-mouth bottles, PFA		
Stir bars, magnetic		Wide-mouth bottles, PP	10)4
Stirring rod, PP		Wide-mouth bottles, PTFE	10)4
January 100, 11				

General Terms and Conditions of VITLAB GmbH

General

- These General Terms and Conditions (GT&C) are intended for use in commercial transactions
- These GT&C shall apply for all, including future, contracts with the customer. Other terms and conditions shall not become part of the contract, even if VITLAB does not expressly object to them. Amendments to and changes of the contract must be in written form. The waiver of the requirement for written form shall only be possible in writing. This shall not apply to individual contractual agreements. The language of the contract shall be German and/or English. In the event of a discrepancy between the German language version of these GT&C and a version in any language, the German language version shall prevail.
- VITLAB offers are subject to change and non-binding. VITLAB reserves the right to make technical improvements to VITLAB products.
- VITLAB may electronically store and process data necessary for the purpose of processing the
- A set-off by the customer shall not be permitted unless the counterclaims are undisputed or legally established, or pecuniary counterclaims arising from the right to refuse payment pursuant to Section 320 Bürgerliches Gesetzbuch (BGB) (German Civil Code). Orders with a goods value of below \leqslant 250 shall be subject to a minimum order surcharge of
- € 50. Delivery shall be undertaken generally in packaging units (PU) according to the currently valid price list. For deliveries within five (5) working days or for order values up to € 500, VITLAB reserves the right to waive an order confirmation.
- For commercial transactions with customers having no general place of jurisdiction in Germany and between businesses, public law persons or special funds under public law the place of jurisdiction shall be the court responsible in Aschaffenburg, Germany. VITLAB shall also be entitled to appeal to the court responsible for the head office of the customer. VITLAB shall, furthermore, as plaintiff have the court responsible for the nead office of the customer. VII.AB shall, furthermore, as plantint made the right to invoke the Arbitration Court at the Chamber of Commerce and Industry in Frankfurt am Main, Germany. The Arbitration Court shall, in this case, make the final judgment in accordance with the Rules of Arbitration of the Chamber of Commerce and Industry in Frankfurt am Main without recourse to the ordinary courts of law. The instigation of legal dunning proceedings by VITLAB shall not signify the exertion of its right of choice; it shall be admissible in all cases.

 German law shall apply exclusively under the exclusion of the conflict of laws principles of Private than the Chamber of the Conflict of Chamber of Conflict (CGC).
- International Law and the UN Convention on Contracts for International the Sale of Goods (CISG).

- **Delivery**The place of performance shall be the factory of VITLAB in Großostheim, Germany. The risk shall transfer to the customer when the goods for delivery are packed and ready for pick-up (EXW (Incoterms® 2010 ex works)). This shall also apply to partial deliveries or where VITLAB has performed additional services, such as shipping; costs for transport, packaging or insurance; exportation and installation. This shall also apply in case of delivery to a consignment warehouse of the customer.
- Insofar as VITLAB has agreed to orders on call, the customer must take delivery prior to the agreed date, otherwise within six (6) months.

Delivery Period, Force Majeure, Delay

- Delivery times shall be ex works. Delivery periods or delivery dates are subject to the customer providing relevant information and documents such as drawings, permits or approvals, opening letters of credit as agreed, making agreed advanced payments and complying with any obligation on time. In particular, the customer is obliged to immediately provide all information required for export, import or shipment (e.g. final recipient, final destination and intended use), documents, approvals and certificates which may be required for the fulfillment of VITLAB's obligations. If delays occur due to required approval, examination or information procedures, delivery periods and delivery dates shall be extended accordingly if VITLAB is not solely responsible for the delay. Delivery shall be subject to VITLAB receiving its own supplies punctually and in good order.
- If delays of delivery periods are due to Force Majeure, e.g. natural disasters, mobilization, war, riot or similar events beyond control of VITLAB, e.g. strike or lock-out, delivery periods shall be extended by the time during which the aforementioned event or its effects persist.
- VITLAB shall immediately inform the customer about delivery barriers or delays and their expected
- VITLAB shall be considered to be in default of delivery only if the customer has issued VITLAB with a reminder, has set a reasonable extension period which has elapsed. The customer shall be
- obliged to immediately inform VITLAB in writing of any likely consequences of delay. In the case of delay damages, VITLAB's liability for compensation shall be limited to 10% of the value of the delayed delivery/service. The limitation shall not apply in cases of willful intent, gross negligence and/or injury to life, limb or health.

- Prices, Terms of Payment
 Prices shall be EXW (Incoterms® 2010 ex works) Großostheim and exclusive of statutory VAT, if applicable. Costs of packaging, transportation, freight and insurance shall be borne by the customer. Prices shall also be exclusive of the cost of returning and recycling/disposing of old
- Invoices shall be payable to VITLAB account in EUROs (€) without deductions and free of charges and expenses. Payment shall be made immediately or by the date stated. The determinant factor shall be the receipt of payment. Cheques and bills of exchange shall only be accepted on account of performance and at the cost of the customer.
- of performance and at the cost of the customer. In the case of customers, with whom VITLAB is working for the first time or with whom VITLAB does not work regularly, after delay in payment or in the case of reasonable doubt as to the creditworthiness of the customer, VITLAB shall reserve the right to make individual deliveries dependent on a pre-payment or a security deposit to the value of the invoice amount. Should the period between conclusion of contract and agreed delivery exceed four (4) months, so may VITLAB, at its discretion, demand a reasonable additional charge equivalent to the increase in its costs up until delivery. 4.5 In the case of an agreed return of goods that are free of defects, the customer shall be charged a checking and processing fee of 20 % of the invoice amount (minimum € 50).
- Should the customer be in arrears with payment, VITLAB debt claims against him shall be due immediately, and VITLAB shall not be obliged to make any further deliveries based on current
- If payment is delayed, VITLAB shall charge notwithstanding further damage compensation claims – interest on arrears at the statutory rate.

 VITLAB may offset amounts payable to the customer (e.g. from credit notes) against VITLAB
- claims against the customer.

- Retention of Title, Assignment of Future Claims
 The goods delivered shall remain property of VITLAB until the complete and unconditional payment.
 Should VITLAB still have further claims against the customer arising from the business relationship,
- Should VILAB still have further claims against the customer arising from the business relationship, VITLAB shall then retain its property rights until payment of such claims has been effected. The customer may neither use goods subject to retention of title nor combine them with other objects to which a third party may have rights. Should, however, goods subject to retention of title become, through their combination with other objects, part of a new (complete) item, VITLAB shall be a proportional co-owner of this new item directly, even if this latter component is regarded as the main component. VITLAB'S proportion of co-ownership shall be determined by the ratio of the invoice value of the goods to the value of the new item at the time of combination.
- The customer may resell goods subject to retention of title in the course of his normal business as long as he has not assigned, pledged or otherwise encumbered his claims from the resale.
- The customer shall assign to VITLAB in advance as collateral any claims against his customers from the resale of the goods subject to retention of title (see Clause 5.3) and/or newly formed items (see Clause 5.2) to the value of VITLAB's invoice for the goods subject to retention of title. As long as the customer is not in default of payment for the goods subject to retention of title, he may collect the assigned claims in the ordinary course of business. He may, however, only use the proportional proceeds for the payment to VITLAB for the goods subject to retention of title.
- At the customer's request, VITLAB shall release collateral at its discretion if and to the extent that the nominal value of the collateral exceeds 120% of the nominal value of our outstanding debt claims against the customer

- 5.6 The customer shall be required to inform VITLAB immediately of any attachments, seizures or any other third-party dispositions relating to the goods that are reserved or co-owned by VITLAB.
 5.7 In the event of delay in payment, failure to pay bills of exchange or cheques, or failure or recall of a payment made via SEPA Direct Debit Scheme, suspension of payments or insolvency of the of a payment made via SEPA Direct Debit Scheme, suspension or payments or insolvency of the customer or of the end buyer, the rights of the customer under Clause 5.3 shall no longer be valid. The customer must then immediately inform the buyer of VITLAB's extended retention of title; he may use the proportional proceeds relating to the assignment only to pay for the delivered goods. VITLAB shall be entitled to collect the assigned receivables itself. In the event of customer's culpable breach of contractual obligations, in particular for the cases covered in Clause 5.7, VITLAB shall be entitled to withdraw from the contract and/or, without withdraw from the contract and/or without with the contract and with the contract and without the contract and with the contract and with the without the contract and with the without the contract and without the contract and without the without the contract and without the contract and without the without the without the contract and without the witho
- withdrawing from the contract, demand the return of any goods subject to retention of title still in possession of the customer and to collect the assigned receivables itself. In order to ascertain the rights of VITLAB, VITLAB shall be entitled to have all of customer's documents/books concerning the reserved rights examined by a person who is subject to the professional duty of confidentiality

Warranty, Limitation of Liability

- VITLAB warrants that its delivered goods (including any agreed installation) are free of defects at time of risk transfer. The required quality, durability and use of VITLAB's delivered goods are based solely on the agreed written specification, product description and/or operating manuals. Any information beyond this, in particular in preliminary discussions, advertising and/or referencing industrial standards shall only become part of the contract if they are expressly referenced in writing. Should the customer require the delivered goods for purposes other than those agreed, he must
- Should the customer require the delivered goods for purposes of oner than those agreed, he must take responsibility himself for examining their special suitability for this also in terms of product safety and ensure their compliance with all relevant technical, legal or regulatory provisions before the intended use. VITLAB shall not be liable for any usability that was not expressly confirmed by VITLAB in writing. In the case of material or design requirements of the customer, VITLAB shall accept no liability for the suitability or permissibility of the desired materials or designs. Compliance with safety-related and occupational health regulations depends on the location and operating conditions of which VITLAB has no prior knowledge. Action for ensuring compliance shall therefore be the responsibility of the customer or his layer.
- compliance shall therefore be the responsibility of the customer or his buyer.

 VITLAB shall not be liable for the consequences of improper handling, use, maintenance and operation of the delivered goods; the consequences of normal wear and tear, in particular of wearing parts, such as pistons, seals, valves; the breakage of glass, plastic or ceramic parts; for the consequences of chemical, electrochemical or electrical influences; or non-observance of the operating instructions.

 If a notice of defect is justified, VITLAB shall initially only be required to provide cure. Any such
- rain another of beginning with a small minday only be required to provide Circ. Any sold free cure shall be, at the discretion of VITLAB, either rectification of the defect or delivery of goods free of defects. Further warranty claims shall only apply in the event of rejection, impossibility or failure of the cure. If cure is possible only at disproportionate expense pursuant to Section 439 Para. 4 BGB (German Civil Code), VITLAB will bear the cost necessary for the cure according to Section 439 Para. 2, 3 BGB (German Civil Code) up to a limit of 150 % of the value of the delivered goods free from defects.
- The customer must, immediately upon receipt of the goods, inspect them carefully, also in terms of product safety, and notify obvious defects immediately in writing; any hidden defects must be immediately notified upon discovery. The customer must notify the carrier immediately of any transport damage. Failure to observe the testing and notification obligation shall void any customer claims for defects.
- Customer claims for defects.

 VITLAB's liability for slight negligence shall be limited to claims for injury to life, limb or health, to claims under the Produkthaftungsgesetz (German Product Liability Act) or to claims of culpable breach of fundamental contractual obligations through which the purpose of the contract is endangered. Otherwise, our liability for slightly negligent breach of fundamental contractual obligations is limited to the typically occurring damages which VITLAB could have foreseen when the contract was concluded.
- Should the customer use the delivered goods in conjunction with environmentally harmful, toxic, radioactive or otherwise hazardous materials, he must notify VITLAB about any such materials and decontaminate the delivered goods prior to returning them to VITLAB. If applicable, VITLAB may charge any necessary costs for decontamination/cleaning and disposal to the customer's account.

Limitation Period

The warranty period shall be one year and starts from the date of delivery of the goods to the customer. The same shall apply for claims for damages, irrespective of their legal basis. The limitation periods of Section 438 Para. 1 Nos. 1 and 2 and Section 634a Para. 1 No. 2 of the BGB (German Civil Code) shall remain unaffected. The warranty period for any rights of recourse pursuant to Section 445a BGB (German Civil Code) shall be one year and starts from the date of delivery of the goods to the customer. The suspension of expiry pursuant to Section 445b BGB (German Civil Code) shall end no later than three years after VITLAB has delivered the goods to the customer. The restrictions of the limitation periods shall not apply to claims based on fraudulent concealment of a defect, for claims under the Produkthaftungsgestez (German Perdukt highlight Act has the restrictions of the limitation periods shall not apply to claims based on fraudulent concealment of a defect, for claims under the Produkthaftungsgestez (German Perdukt highligh Act highligh Act has the produkt highlight and those should be produkt highlight and those should be producted the produkt highlight and those should be producted by the product highlight and those productions are producted by the production of the product highlight and those productions are producted by the production of t Product Liability Act) or for damages resulting from injury to life, limb or health and other damages based on intent or gross negligence. The limitation period in respect of replaced or repaired goods shall only commence anew if VITLAB admitted the defectiveness of the replaced or repaired goods.

Software Use

- If software is included in the scope of a delivery, the customer shall be granted a non-exclusive right to use the software and its associated documentation. It is provided for use on the designated delivery item. The use of the software on more than one system shall be prohibited.
- designated delivery item. The use of the software on more than one system shall be prohibited. The customer shall only be entitled to copy, transfer or translate the software or to convert it from object code to source code to the extent permitted by law (Sections 69a et seq. Urheberrechtsgesetz German Copyright Act). The customer undertakes to refrain from removing manufacturer information, in particular copyright notices, or from changing these without VITLAB's prior express consent or the prior express consent of the software supplier. All other rights to the software and the documentation including copies thereof shall remain with VITLAB and/or the software supplier. The issue of sublicenses is not permitted.

Spare Parts, Maintenance/Repair and Calibration

- For spare parts and maintenance, repair and calibration services, the current repair service and
- For spare parts and maintenance, repair and calibration services, the current repair service and spare parts price list shall apply. Insofar as there is an obligation on the part of VITLAB to maintain/supply spare parts, then this obligation shall be limited to a period of five (5) years from the date of delivery. If spare parts are not manufactured by VITLAB, or are no longer available on the market, for example electronic components, or if the raw material for their production is no longer available, the obligation of VITLAB to deliver spare parts shall be considered. VITLAB to deliver spare parts shall lapse
- For calibration and maintenance, expendable items from VITLAB production are normally used. Maintenance and calibration services can only be provided if the customer has declared the devices sent to be safe to work on from a health hazard perspective.
- For repair/service values of up to € 50, VITLAB reserves the right not to provide a separate cost

- Legal Reservation, Industrial Property Rights, Confidentiality
 VITLAB reserves ownership and all industrial property rights and copyrights to all moulds, tools or other devices, samples, pictures, and business and technical documents produced or provided by VITLAB. This also applies where the customer has wholly or in part taken on the costs hereof. The customer may use these only in the manner agreed with VITLAB. Without VITLAB's written consent, the customer may not himself manufacture contractual objects delivered or have the
- same manufactured by third parties.

 10.2 Insofar as VITLAB delivers goods according to the designs or other requirements specified by the customer (models, patterns etc.), the customer shall be liable to VITLAB by default for ensuring that, through the manufacture and delivery of these goods, the industrial property rights or other rights of third parties are not infringed. If the customer is at fault he shall reimburse VITLAB all damage resulting from any such infringement of rights.

 10.3 Any information acquired from this business relationship and not deemed to be public knowledge
- must not be disclosed by the customer to third parties

Status as of: October 2018

Edition 2019

VITLAB GmbH

Linus-Pauling-Str. 1 63762 Grossostheim

Germany

tel: +49 6026 9 77 99-0 fax: +49 6026 9 77 99-30 e-mail: info@vitlab.com www.vitlab.com USt.-IdNr. / VATREG NO

DE 1116669 59

Banking detailsSparkasse Aschaffenburg
Volksbank Main-Tauber eG

Volksbank Main-Taub Commerzbank AG Deutsche Bank AG WEEE-Reg.-Nr. DE 30031601

IBAN
DE91 7955 0000 0000 0003 15
DE03 6739 0000 0034 2765 01
DE63 7908 0052 0309 9404 00
DE49 5087 0005 0010 5619 00

SWIFT-BIC

BYLA DE M1 ASA

GENO DE 61 WTH

DRES DE FF 790

DEUT DE FF 508